

Quantification of Urea in Ethanol by Raman Spectroscopy

Urea is widely employed as a nitrogen-release fertilizer with more than 90% of urea production destined for agricultural applications[1]. Urea is also known to form complexes with fatty acids[2], which have been employed for separation of complex mixtures and purification

processes[3]. In this application note, we present the quantification of the concentration of urea in ethanol by Raman Spectroscopy and show how this method can be employed for determining the percentage of urea in a solid inclusion compound with stearic acid[4].

Find more information in the video:

EXPERIMENTAL

Chemicals: Urea (Aldrich, >99%); Stearic acid (SA) (Aldrich, >95%); Ethanol (Soria)

Instrument: [i-Raman® Plus 785S](#)

Stock solutions of urea (0.0420 gr/gr ethanol) and SA (0.04126 gr/gr ethanol) were prepared. Standard solutions were prepared by mixing these stock solutions in different proportions in order to obtain urea concentrations between 0 and 0.042 gr/gr ethanol and keeping the total mass concentration (urea + SA) approximately

constant. Raman spectra of 0.5 mL of each solution in aluminum containers were measured employing a 785 nm laser (power: 30%, ~ 90 mW) and 5000 msec of acquisition time (20 repetitions). No effects of heating or evaporation were observed. Spectra of the standard solutions that were background corrected employing the [BWSpec®](#) software are shown in **Figure 1**.

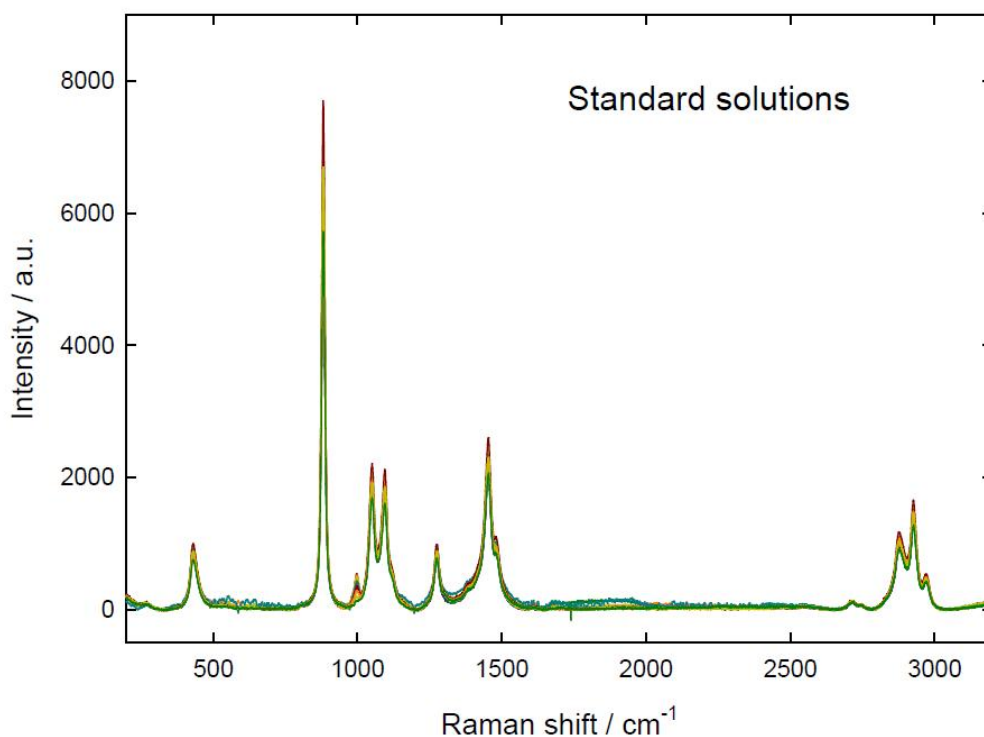


Figure 1. Dark-subtracted baseline-corrected Raman spectra of the standard solutions of urea and SA in ethanol.

RESULTS

Spectra were normalized by employing the intensity of the ethanol band at $1049\text{--}1050\text{ cm}^{-1}$. Normalized spectra (**Figure 2**) clearly show that the only appreciable change due to the increment in relative urea concentration occurs at the band at $996\text{--}997\text{ cm}^{-1}$, which is attributed to urea. This band corresponds to the symmetric C-N stretching[5] that is experimentally and theoretically reported at about 1010 cm^{-1} for the solid urea,[6] but shifts to lower wavenumbers in solution.[5,7]

For quantification purposes, spectra were deconvoluted, fitting the experimental results in the region $950\text{--}1200\text{ cm}^{-1}$ by 4 Lorentzian functions. These curve fitting results are shown in **Figure 3** for some of the standard solutions.

The ratio of the intensities of the fitted peaks assigned to urea at 996 cm^{-1} (peak 1, a1) and ethanol at 1049 cm^{-1} (peak 2, a2) was employed

as the analytical parameters. The dependence of this ratio on the urea concentration of the samples is presented in **Figure 4**. The calibration curve plotted in this figure shows a good linear behavior which indicates this parameter could be employed for urea quantification.

For the determination of the urea content of real samples containing both urea and SA, the solid sample was dissolved in ethanol (0.04299 gr/gr ethanol) and Raman spectra were recorded in the same conditions. From the values of the ratio $a1/a2$ for this sample, obtained by fitting of the peaks at 996 and 1049 cm^{-1} (**Figure 5**), the urea concentration of the solution was determined to be $0.03274\text{ gr urea /gr ethanol}$. Thus, the urea content of the sample resulted to be 76% w/w. This value is consistent with other reported values for the inclusion compounds formed by stearic acid and urea (ca. 75%).[2,4]

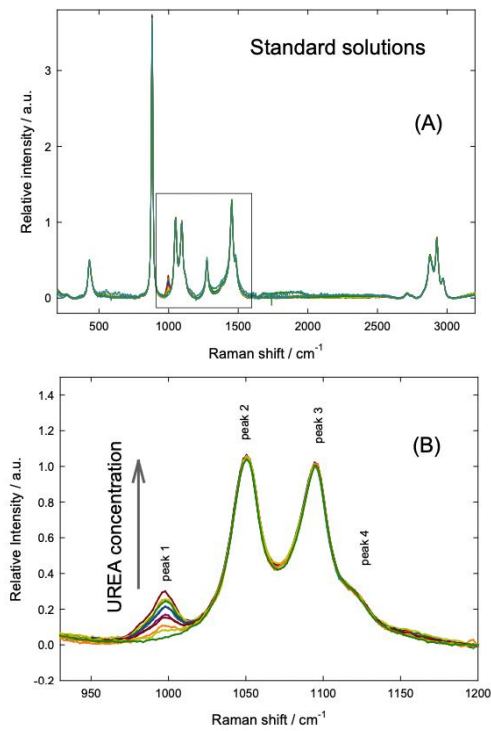


Fig.2: Normalized spectra of the standard solutions of urea + SA in ethanol. (A) Complete spectra (B) Analyzed region.

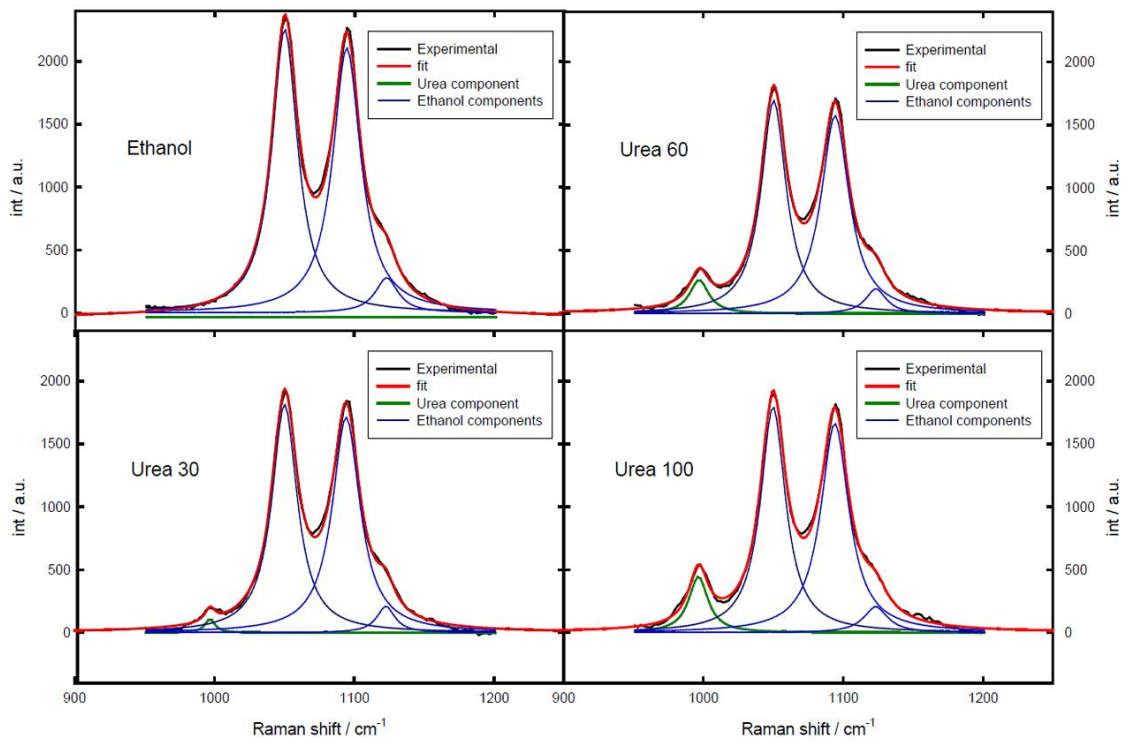


Figure 3. Fitting of the Raman spectra in the region 950-1200 cm^{-1} . Urea 30, Urea 60 and Urea 100 mean solutions 0.0123, 0.0248 and 0.0413 gr urea/gr ethanol respectively.

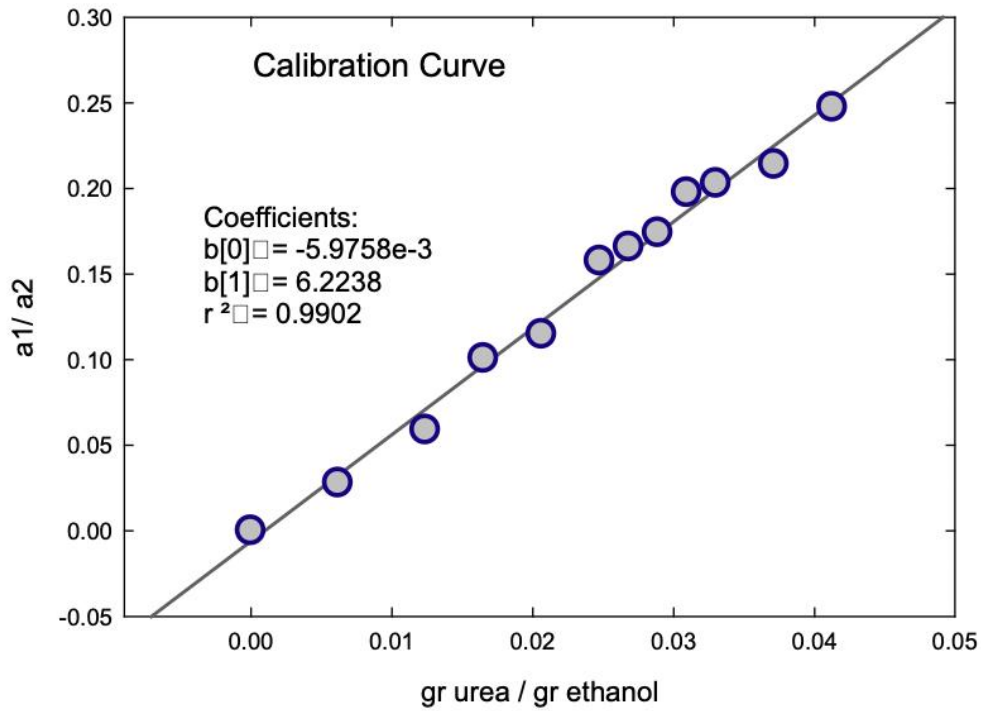


Figure 4. Calibration curve for the urea quantification in ethanol. Ratio of the intensities of the fitted bands of urea (a1) and ethanol (a2) as a function of the urea content of the standard solutions.

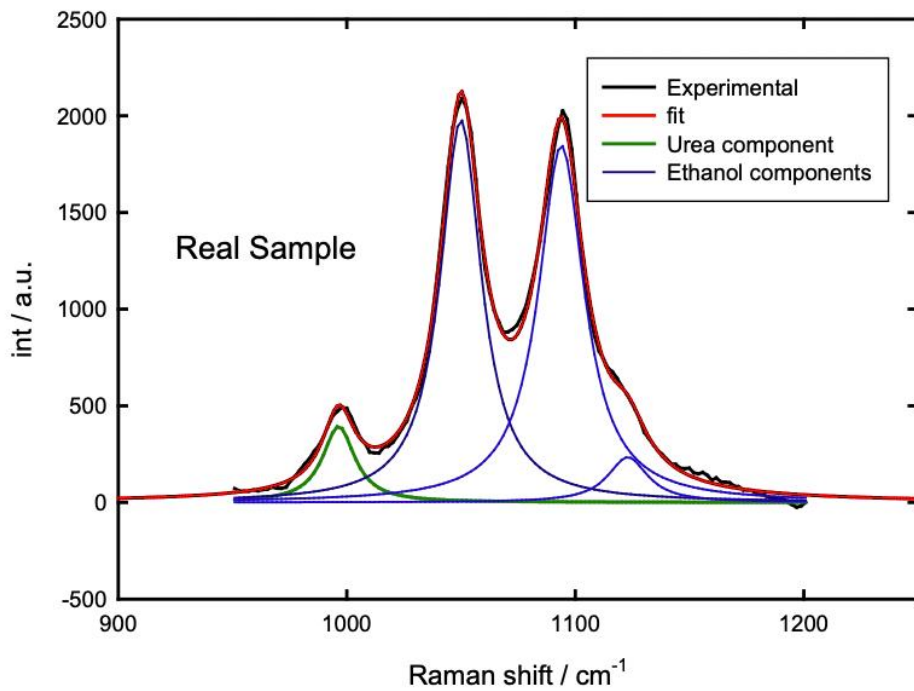


Figure 5. Fitting of the sample spectrum.

CONCLUSION

We presented a simple method for quantification of urea concentration in ethanolic solutions by Raman spectroscopy. The calibration curve presents good linearity in the concentration range analyzed (up to 0.042 gr urea / gr ethanol). The presence of stearic acid in

the samples does not modify appreciably the Raman spectra (at least up to 0.042 gr / gr ethanol), so this method allows the quantification of urea in solid binary samples containing both urea and stearic acid.

FURTHER INFORMATION

Related application notes

[Choosing the Most Suitable Laser Wavelength](#)

[Quantification of methanol in contaminated spirits with Raman](#)

ACKNOWLEDGEMENTS

We would like to thank Dr. Waldemar A. Marmisollé of the Soft Matter Laboratory INIFTA – Universidad Nacional de La Plata (UNLP),

Buenos Aires, Argentina – CONICET for sharing these research results.

REFERENCES

1. J. H. Meessen, H. Petersen, in Ullmann's Encycl. Ind. Chem., Wiley-VCH Verlag GmbH & Co. KGaA, 2000.
2. H. Schlenk, R. T. Holman, J. Am. Chem. Soc. 1950, 72, 5001–5004.
3. D. G. Hayes, Y. C. Bengtsson, J. M. Van Alstine, F. Setterwall, J. Am. Oil Chem. Soc. 1998, 75, 1403–1409.
4. A. Strocchi, G. Bonaga, Chem. Phys. Lipids 1975, 15, 87–94.
5. R. Keuleers, H. O. Desseyn, B. Rousseau, C. Van Alsenoy, J. Phys. Chem. A 1999, 103, 4621.
6. B. Rousseau, C. Van Alsenoy, R. Keuleers, H. O. Desseyn, J. Phys. Chem. A 1998, 102, 6540–6548.
7. D. Gangopadhyay, S. K. Singh, P. Sharma, H. Mishra, V. K. Unnikrishnan, B. Singh, R. K. Singh, Spectrochim. Acta Part A Mol. Biomol. Spectrosc. 2016, 154, 200–206.

CONTACT

メトロームジャパン株式会社
143-0006 東京都大田区平
和島6-1-1
null 東京流通センター アネ
ックス9階

metrohm.jp@metrohm.jp

CONFIGURATION



i-Raman Plus 785S

i-Raman® Plus 785S は、弊社の革新的かつインテリジェントなスペクトロメーターテクノロジーを備えた i-Raman ホータフル型ラマンスペクトロメーターの受賞歴のあるシリーズの構成要素です。このホータフル型ラマンスペクトロメーターは、高い量子効率、TE 冷却、ならびに高いタイナミックレンジを備えた CCD アレイ検出器を用い、積分時間でさえ最長30分、低ノイズの傑出した性能を提供します。こうして、弱いラマン信号も測定することかできます。

i-Raman Plus 785S は、 65 cm^{-1} から 3350 cm^{-1} までの測定を可能にするコンフィクレーションを有する幅広いスペクトル範囲と高分解能のユニークな組み合わせを特徴としています。システムの小さな底面、軽量構造、そして低消費電力により、とこでもラマン分析を研究レベルで実施することかできます。i-Raman Plus には、サンプル採取を簡単にする光ファイバーフローセルが装備されており、キュベットホルター、ヒテオマイクروسコーフ、フローホルター付き XYZ スライドテーブル、ならびに弊社独自の多変量解析ソフトウェア BWIQ® および同定ソフトウェア BWID® と共に使用することかできます。i-Raman Plus により、品質分析および定量分析のための高精度のラマンソリューションを常に使用することかできます。



9.5 mm

BCR100A ラマンキュベットホルターにより、ラマンフローセルをホルターに固定することで液体および粉末の容易なラマンスペクトル測定かできるようになります。この付属品は、いまたかつてない再現性のために三点精密ねし付きの内蔵ミラーを使用しており、これによりラマン信号を従来のキュベットホルターの3倍まで増強します。これは、フローシャフトがキュベットに直接接触しないように設計されており、またバックラウンド蛍光を減量するためのライトトラップも含まれています。BCR100A は、直径 9.5 mm または 12 mm のフローセルのための設計モデルで入手可能で、外径 12.5 mm x 12.5 mm の各標準キュベット (波長 1 cm) とともに液体や粉末のサンプル採取に用いられます。