



Application Note AN-T-195

Determination of nitrogen content

Kjeldahl determination according to USP general chapter <461>

SUMMARY

The Kjeldahl method is used to determine the nitrogen content in organic and inorganic samples. The Kjeldahl analysis consists of three steps: digestion, distillation, and titration. During the catalytic digestion step, organic nitrogen (except nitro- and azo-groups and nitrogen in rings) is converted into ammonium. Sodium hydroxide is added just before the distillation step in order to convert the ammonium into ammonia. Through steam distillation, the latter is

transferred into the receiver vessel containing an absorbing agent (e.g., boric acid). Finally, the separated ammonia is titrated with sulfuric acid.

Protein content in samples can also be determined from the nitrogen content obtained by Kjeldahl setup. USP<461> describes the titration method to determine nitrogen content in organic products using Kjeldahl nitrogen setup. This Application Note illustrates nitrogen determination in heparin sodium.

The analysis is demonstrated on heparin sodium. An appropriate amount of heparin sodium is weighed accurately and transferred into a clean 2-neck round bottom flask. Sodium sulfate, copper sulfate, and sulfuric acid are added for the digestion process. The content is gently heated below the boiling point until the frothing ceases. It is then heated again at a higher temperature until the content boils and the solution becomes brown. The contents are cooled and carbon

dioxide-free water is added carefully while thoroughly mixing. Sodium hydroxide solution is added through the side-neck of the round bottom flask. Granulated zinc is added and the flask is connected immediately to the Kjeldahl nitrogen distillation setup. The outlet of the setup is put into a solution of boric acid. The distillation is carried out until approximately 80% of the total volume is transferred into the boric acid solution.

EXPERIMENTAL

The analysis is performed automatically on a Titrando system consisting of a 905 Titrando. The Unitrode is used for the indication of the titration curve.

The prepared sample is titrated potentiometrically against standardized sulfuric acid until after the first equivalence point.



Figure 1. Example of a Titrando system consisting of a 905 Titrando and a 900 Touch Control. Alternatively the 905 Titrando can also be connected to a PC and controlled by tiamo.

RESULTS

Sharp titration curves are obtained where the equivalence point is reliably determined by the Touch Control or *tiamo*TM.

The determined nitrogen content of heparin sodium

is 1.581% (SD(rel) = 1.48%, n = 5), which is within the nitrogen content specified by USP (1.3% to 2.5%) for heparin sodium.

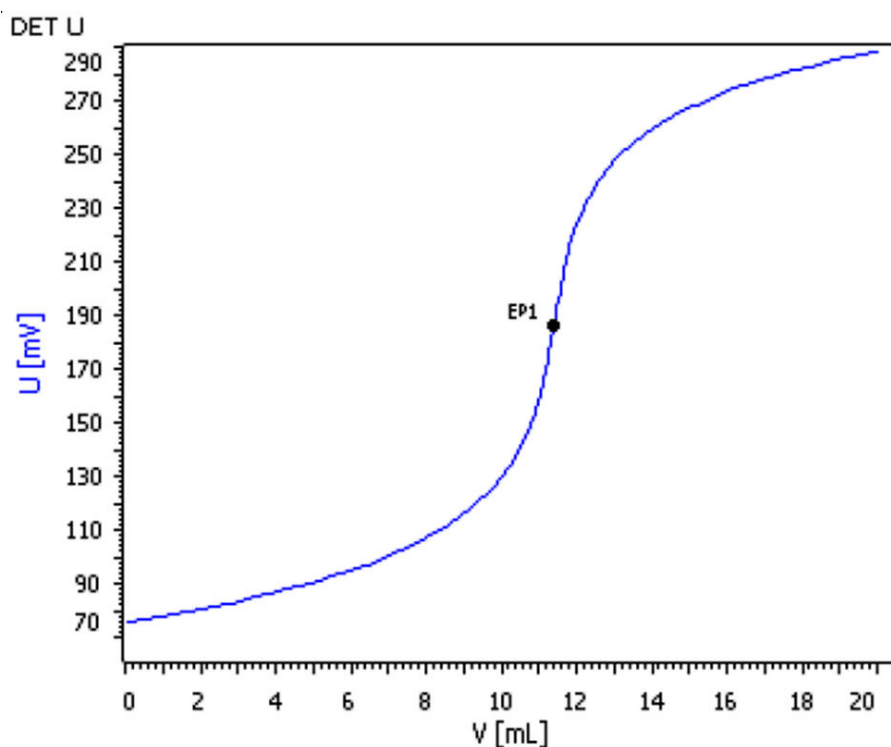


Figure 2. Titration curve of the Kjeldahl determination of heparin sodium according to USP general chapter <461>.

CONCLUSION

This method shows the possibility to determine the nitrogen content in various kinds of samples automatically, accurately, and reliably by titration according to the **USP general chapter <461>**.

Aside from heparin sodium, the following compounds can also be analyzed with this method:

- Antithrombin III Human
- Beta Glucan
- Cellulose sodium phosphate
- Chlorophyllin copper complex sodium
- Colloidal oatmeal
- Copovidone
- Crospovidone
- Dalteparin sodium
- Dextran 1
- Dextrin
- Dihydroxyaluminum aminoacetate,
- Enoxaparin sodium etc.
- Guar gum
- Mecamylamine Hydrochloride tablets
- Melphalan
- Polyvinyl acetate dispersion
- Povidone iodine
- Povidone
- Psyllium hemicellulose
- Pullulan
- Racemic calcium pantothenate
- Ralbumin human
- Scaffold bovine dermis
- Spirulina tablet
- Taurine
- Thioguanine
- Trehalose
- Wheat bran
- Zein

CONTACT

Metrohm Nordic ApS
Lyskær 3C. 1. th.
2730 Herlev

mail@metrohm.dk

CONFIGURATION



905 Titrando

High-end titrator for potentiometric titration with one measuring interface for use with Dosino dosing systems.

- up to four dosing systems of the 800 Dosino type
- dynamic (DET), monotonic (MET), and endpoint titration (SET)
- Measurement with ion-selective electrodes (MEAS CONC)
- Dosing functions with monitoring, Liquid Handling
- four MSB connectors for additional stirrers or dosing systems
- "iTrode" intelligent electrodes
- USB connector
- For use with OMNIS Software, *tiamo* software, or Touch Control unit
- Compliance with GMP/GLP and FDA regulations such as 21 CFR Part 11, if required



iUnitrode with Pt1000

Intelligent, combined pH electrode with integrated memory chip for storing sensor data and Pt1000 temperature sensor. This electrode is particularly suitable:

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte , storage in Idrolyte .

iTrodes can be connected to Titrando, Ti-Touch or 913/914 meters.