

Raman spectroscopy analyses of road construction materials

Various solid materials used in road construction were analyzed using a handheld Raman spectrometer. The investigated materials are commonly used pigments and resins such as CaCO_3 , TiO_2 , and DEGALAN®. The

measured spectra differed remarkably from each other. To evaluate the main differences in the chemical structures, the peaks of the different spectra were assigned to the functional groups causing them.

INTRODUCTION

Different materials such as paints, (white) pigments, and resins are commonly used in the construction of roads. Together with glass beads (used for night-visible road marking) and various other materials, they help us to get safely from A to B.

In this study, different roadmaking materials were

analyzed using the Mira M-1 handheld Raman analyzer. The gathered spectra were compared with each other to see the main differences in the functional groups. The analysis proved that Mira M-1 is suited for the differentiation of such materials.

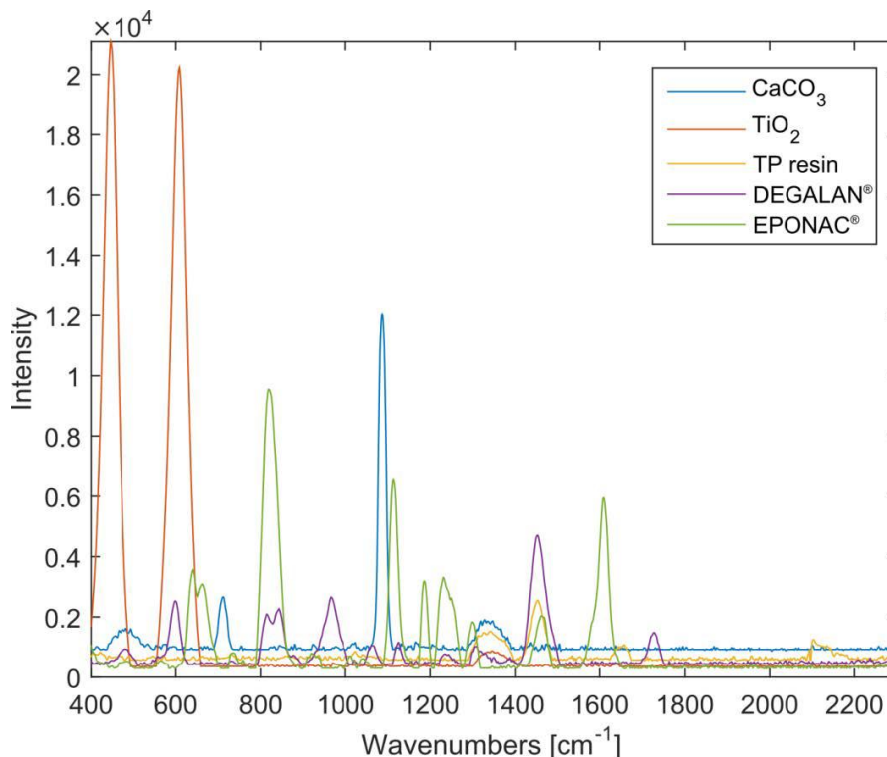


Figure 1. Overlay of spectra of different roadmaking materials

EXPERIMENTAL

All spectra were measured using the Mira M-1 Raman spectrometer in auto-acquisition mode, i.e., integration times were determined automatically. A laser wavelength of 785 nm and the Orbital-Raster-Scan (ORS) technique were applied. The measurements were done in small sample vials with the vial holder adapter.

The following samples were investigated:

1. Chalk (CaCO_3)
2. Titanium dioxide (TiO_2)
3. EPONAC® resin
4. TP resin
5. Pigment yellow
6. Pigment blue
7. Pigment red
8. DEGALAN® resin

RESULTS AND DISCUSSION

The measurement of calcium carbonate furnished a clear spectrum with two main peaks at 712 cm^{-1}

(symmetric O–C–O bending vibration) and 1087 cm^{-1} (symmetric stretching vibration).

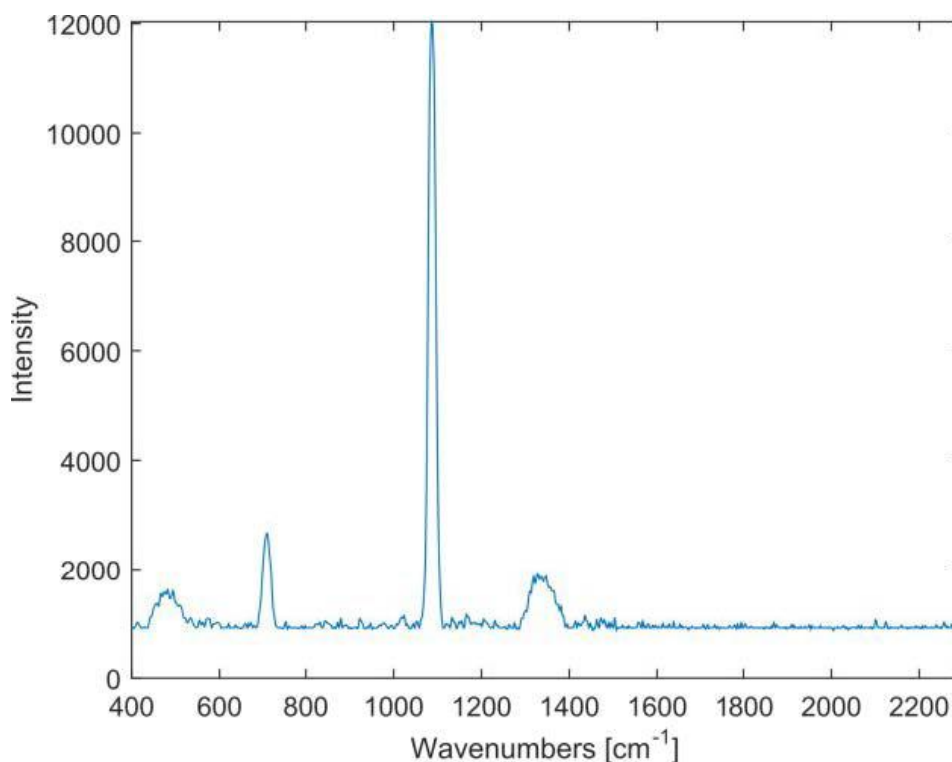


Figure 2. Spectrum of calcium carbonate.

Looking at the titanium dioxide spectrum, one can see two main peaks which give information on the present crystal modification (rutile or anatase). Two peaks represent rutile, while three peaks are typical for anatase because of its crystal symmetry.

Both peaks are symmetric stretching vibrations and belong to O–Ti–O (446 cm^{-1}) and Ti–O (609 cm^{-1}). In the case of anatase, the peak at 446 cm^{-1} is split in two. This makes it easy to differentiate between rutile and anatase.

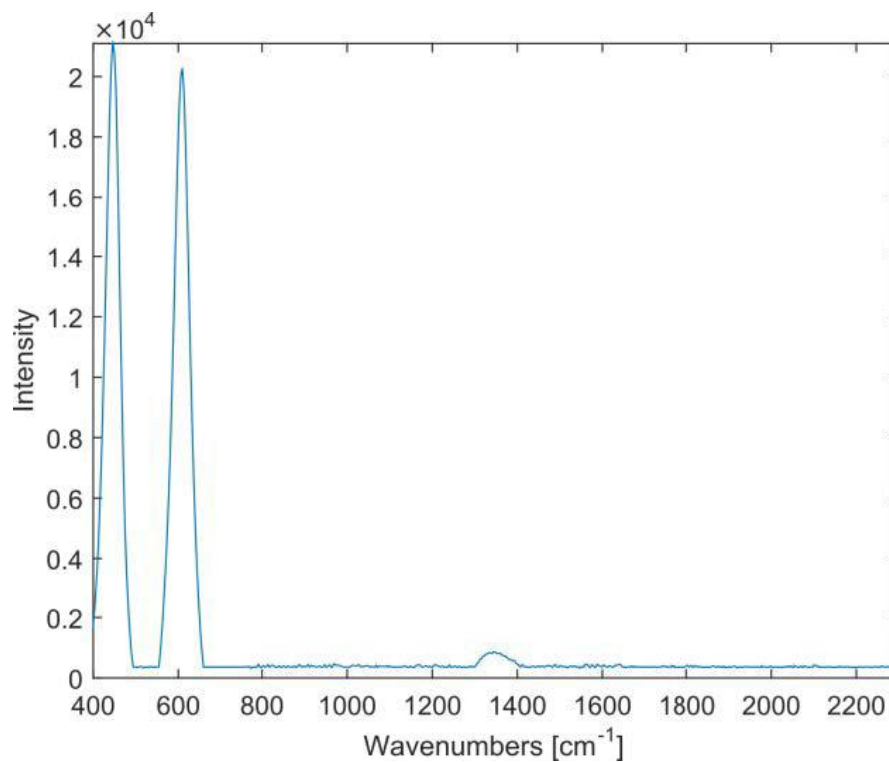


Figure 3. Spectrum of TiO₂ (rutile).

A glance at the spectra of typical resins such as EPONAC® or DEGALAN® reveals nicely separated

peaks which can be assigned to their corresponding functional groups (see below).

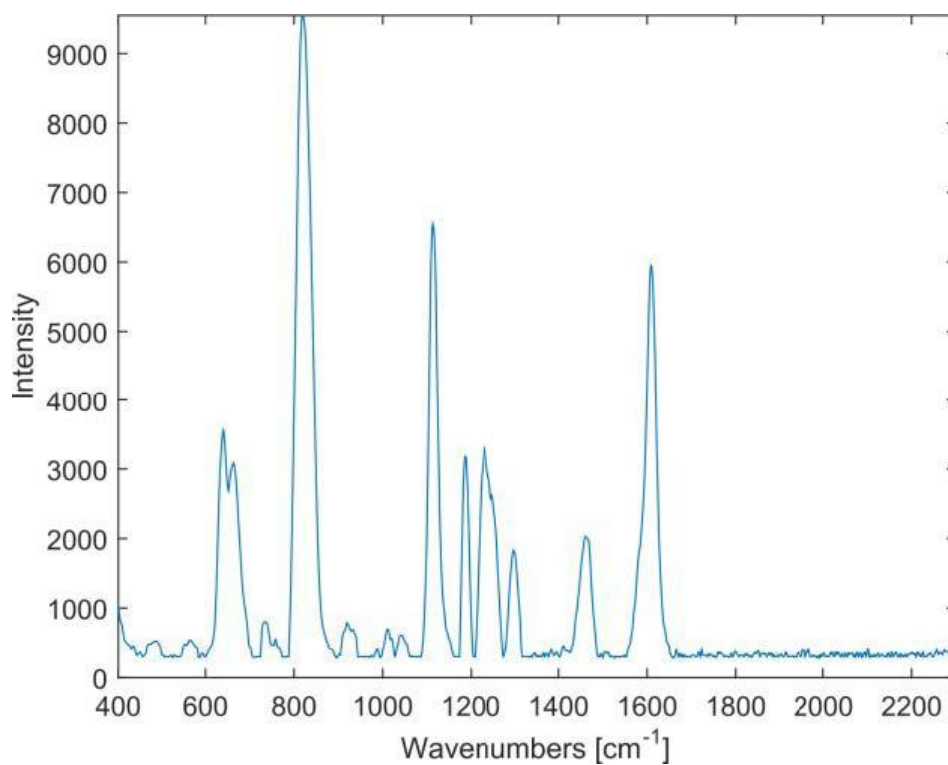


Figure 4. Spectrum of the EPONAC® resin.

Table 1. Peaks observed in the EPONAC® spectrum.

Peak [cm ⁻¹]	Description
640	Cyclic vibration (para-substituted benzene)
819	C–H bending vibration (para-substituted benzene)
1000	Various C–C stretching vibrations
1113	C–OH stretching vibration
1189	C–(CH ₃) ₂ stretching vibration
1231	C–O
1248	C–H (benzene)
1298	–CH ₂ – twisting vibration
1461	–CH ₂ – bending vibration
1609	C=C

The peaks are very characteristic for para-substituted benzene and confirm that EPONAC® is a copolymer of bisphenol A (BPA) and another component.

When comparing the EPONAC® spectrum with the DEGALAN® spectrum (see **Figure 5**), it is obvious that

the peak of the benzene at 1600 cm^{-1} is missing. The peak around 1700 cm^{-1} , together with the peaks slightly below and above 1200 cm^{-1} , is characteristic for carbonyl groups. In addition to that, the C–C peaks are more distinct for EPONAC® than for DEGALAN®.

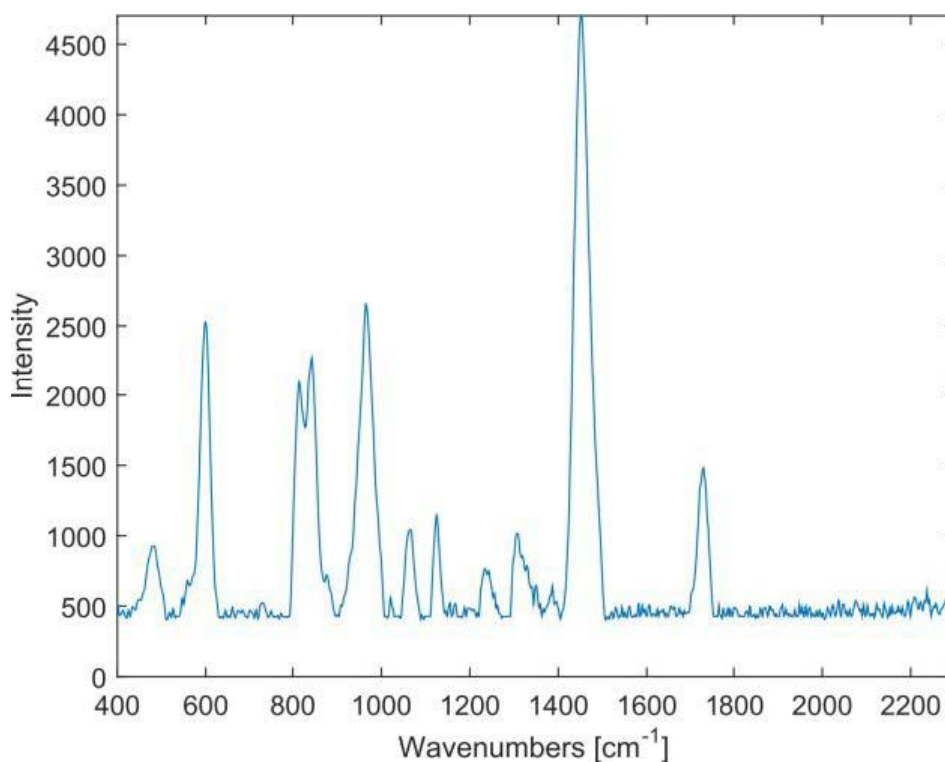


Figure 5. Spectrum of the DEGALAN® resin

Table 2. Peaks observed in the DEGALAN® spectrum.

Peak [cm ⁻¹]	Description
599	–COO bending vibration
814	Propionate bending vibration
843	C–CH ₃ stretching vibration
965	C–C stretching vibration
1065	C–COO stretching vibration
1125	C–O
1234	C–O
1308	–CH ₂ – twisting vibration
1452	–CH ₂ – bending vibration
1728	C=O carbonyl group, ester

Using these differences, it is easy to discriminate between the resins and also between resins and

pigments (see Figure 6).

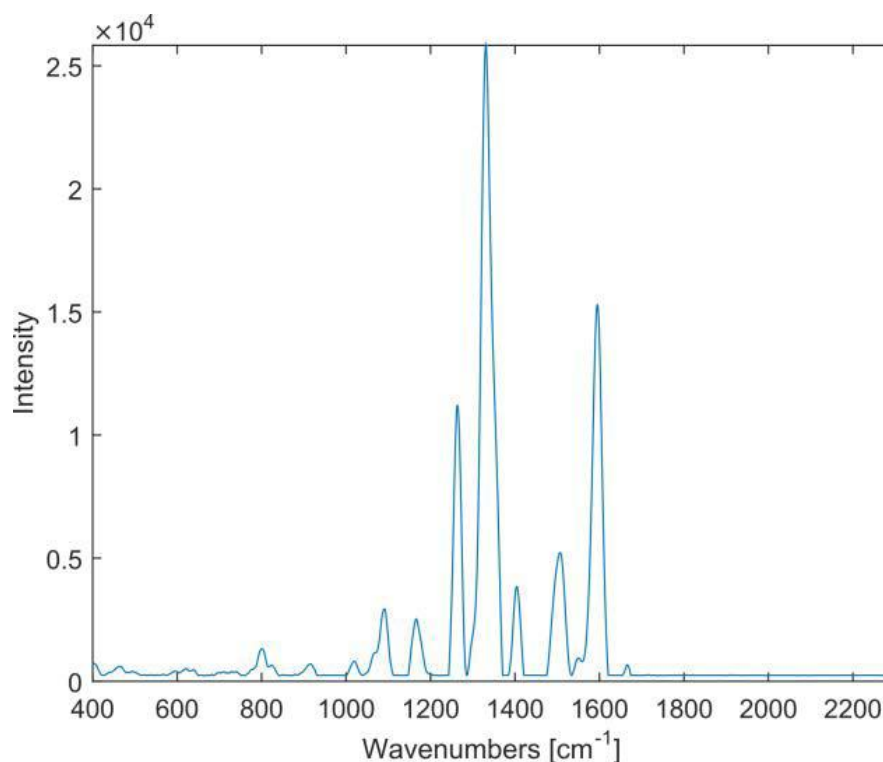


Figure 6. Spectrum of the yellow pigment «Hansa®-Brilliantgelb»

CONCLUSIONS

Because of the great differences in their spectra, handheld Raman spectroscopy is ideally suited for the analysis of materials used in road construction. The investigation of the spectra showed that there are

significant differences in the functional groups of the materials, thus allowing the identification with handheld Raman systems such as Mira M-1.

CONTACT

Metrohm France
13, avenue du Québec - CS
90038
91978 VILLEBON
COURTABOEUF CEDEX

info@metrohm.fr

CONFIGURATION



MIRA P Advanced

Le Metrohm Instant Raman Analyzer (MIRA) P est un spectromètre Raman portable performant qui s'utilise pour les déterminations rapides et non destructives et le contrôle des matériaux les plus divers, comme les principes actifs pharmaceutiques et les excipients. De très petite taille, le MIRA P est pourtant très robuste et dispose d'une structure de spectrographe haute efficacité, équipée de notre technologie « Orbital Raster Scan » (ORS) inédite. MIRA P satisfait aux prescriptions FDA 21 CFR partie 11.

Le Advanced Package comprend une lentille avec laquelle les matériaux peuvent être analysés directement ou dans leur conditionnement (classe de laser 3b), ainsi qu'un support de flacon pour analyser les échantillons dans des flacons en verre (classe de laser 1).