



Application Note AN-NIR-110

Quality control of sugarcane juice

Multiparameter determination within one minute using NIRS

Sugarcane (*Saccharum* spp.) is a very important crop for the global economy. It is often used as a raw material for the production of sugar, alcohol, yeast, and more. Brix (° Brix), Pol (%), juice purity (%), reducing sugars (%), and total recoverable sugars (Kg t^{-1}) are some of the many quality control (QC) parameters that must be analyzed in sugarcane juice.

Numerous methods based on several analytical

techniques are available for sugarcane juice quality control. These methods can be quite time-consuming since sample treatment is a prerequisite. A quicker alternative to these other methods is near-infrared spectroscopy (NIRS). NIRS allows the simultaneous determination of several QC constituents, without chemicals or sample preparation, in less than one minute.

Sugarcane juice was analyzed by NIRS, and a total of 100 spectra were collected to create a prediction model for quantification of several QC parameters. All samples were measured with a Metrohm NIRS DS2500 Liquid Analyzer (400–2500 nm) in transmission mode with a

DS2500 Holder Flow Cell (**Figure 1**). A flow cell with 1 mm pathlength was used in this study. This flow cell was filled via peristaltic pump. The Vision Air Complete software package from Metrohm was used for all data acquisition and prediction model development.

Table 1. Hardware and software equipment overview.

Equipment	Article number
DS2500 Liquid Analyzer	2.929.0010
DS2500 Holder Flow cell	6.7493.000
NIRS quartz cuvette flow 1 mm	6.7401.310
Vision Air 2.0 Complete	6.6072.208





Figure 1. Metrohm NIRS DS2500 Liquid Analyzer and DS2500 Holder Flow Cell used for the fast determination of several QC parameters in sugarcane juice.

The obtained Vis-NIR spectra (**Figure 2**) were used to create a prediction model for quantification of Brix ($^{\circ}$ Brix), Pol (%), juice purity (%), reducing sugars (%), and total recoverable sugars (Kg t^{-1}). The quality of the prediction model was evaluated using

correlation diagrams which display a very high correlation between the Vis-NIR prediction and the reference values. The respective figures of merit (FOM) display the expected precision of a prediction during routine analysis (**Figures 3–8**).

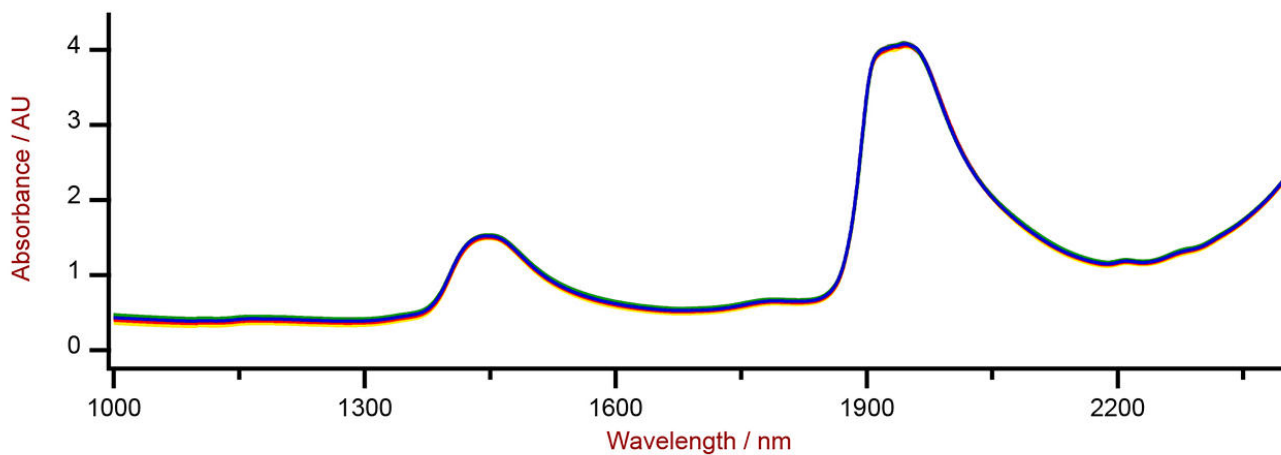


Figure 2. Selection of Vis-NIR spectra of sugarcane juice samples analyzed on a DS2500 Liquid Analyzer with a 1 mm pathlength flow cell.

RESULT BRIX

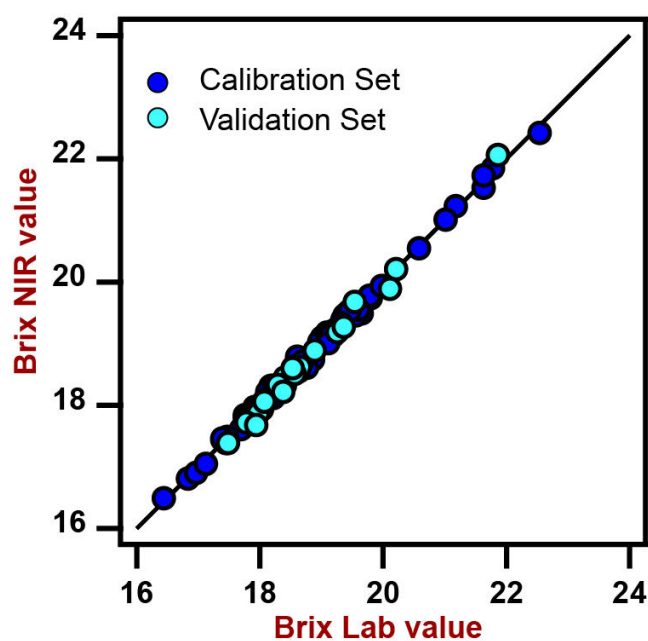


Figure 3. Correlation diagram and the respective FOMs for the prediction of Brix in sugarcane juice using a DS2500 Liquid Analyzer. Laboratory values were evaluated using a refractometer.

Figures of Merit	Value
R^2	0.9875

Standard Error of Calibration	0.1323 (° Brix)
Standard Error of Cross-Validation	0.1467 (° Brix)
Standard Error of Validation	0.138 (° Brix)

RESULT POL

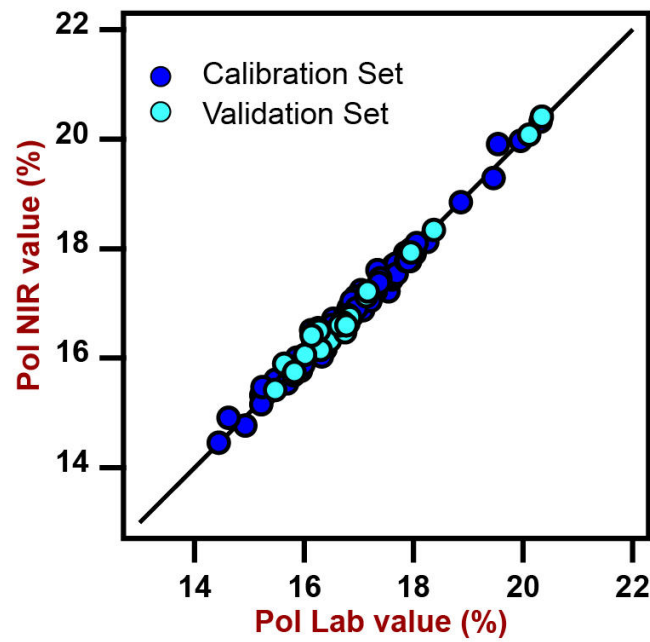


Figure 4. Correlation diagram and the respective FOMs for the prediction of Pol in sugarcane juice using a DS2500 Liquid Analyzer. Laboratory values were calculated from the sucrose reading, Brix, and a few constants.

Figures of Merit	Value
R ²	0.9833
Standard Error of Calibration	0.1506%
Standard Error of Cross-Validation	0.1851%
Standard Error of Validation	0.1388%

RESULT PURITY

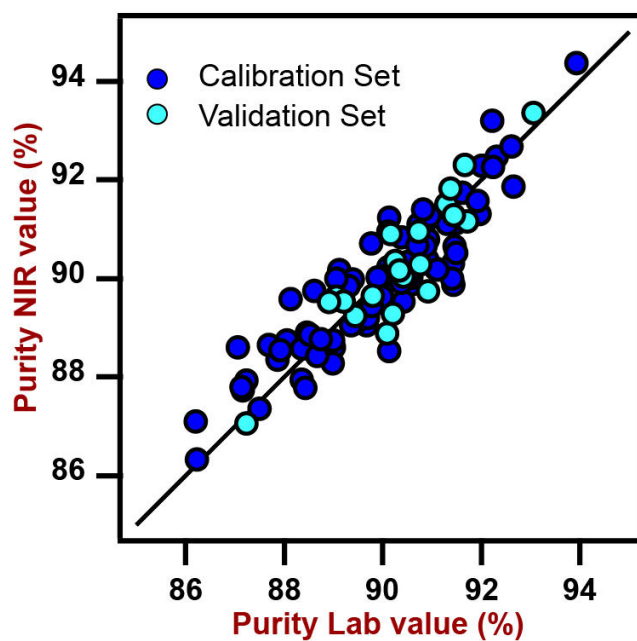


Figure 5. Correlation diagram and the respective FOMs for the prediction of sugarcane juice purity using a DS2500 Liquid Analyzer. Laboratory values were calculated using the results from Pol and Brix determinations: Purity = $100 \times (\text{Pol}/\text{Brix})$.

Figures of Merit	Value
R^2	0.8194
Standard Error of Calibration	0.7202%
Standard Error of Cross-Validation	0.7596%
Standard Error of Validation	0.564%

RESULT REDUCING SUGARS

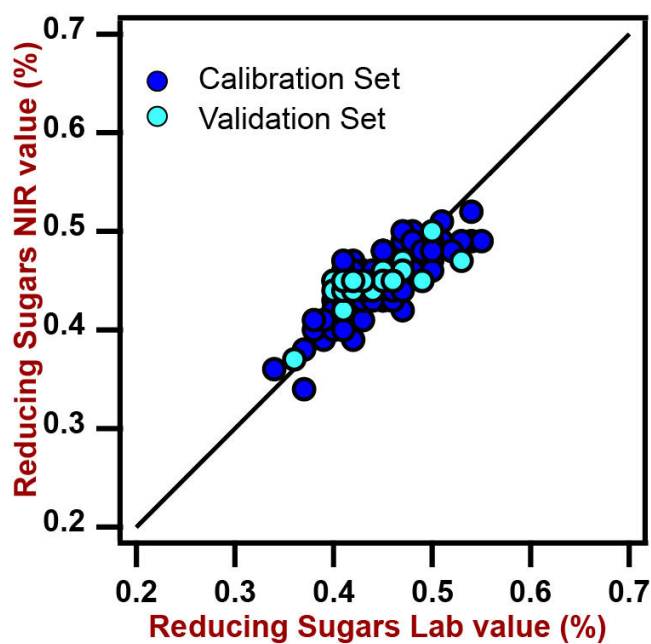


Figure 6. Correlation diagram and the respective FOMs for the prediction of reducing sugars in sugarcane juice using a DS2500 Liquid Analyzer. Laboratory values were measured with ion chromatography (IC).

Figures of Merit	Value
R^2	0.6497
Standard Error of Calibration	0.0263%
Standard Error of Cross-Validation	0.0291%
Standard Error of Validation	0.0249%

RESULT SUCROSE

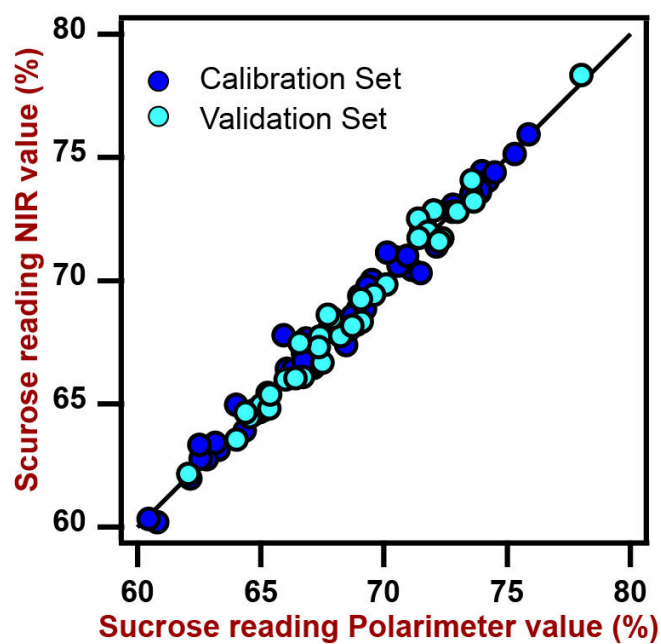


Figure 7. Correlation diagram and the respective FOMs for the prediction of sucrose reading in sugarcane juice using a DS2500 Liquid Analyzer. Laboratory values were evaluated with a polarimeter.

Figures of Merit	Value
R^2	0.9911
Standard Error of Calibration	0.5388%
Standard Error of Cross-Validation	0.6604%
Standard Error of Validation	0.497%

RESULT TOTAL RECOVERABLE SUGARS

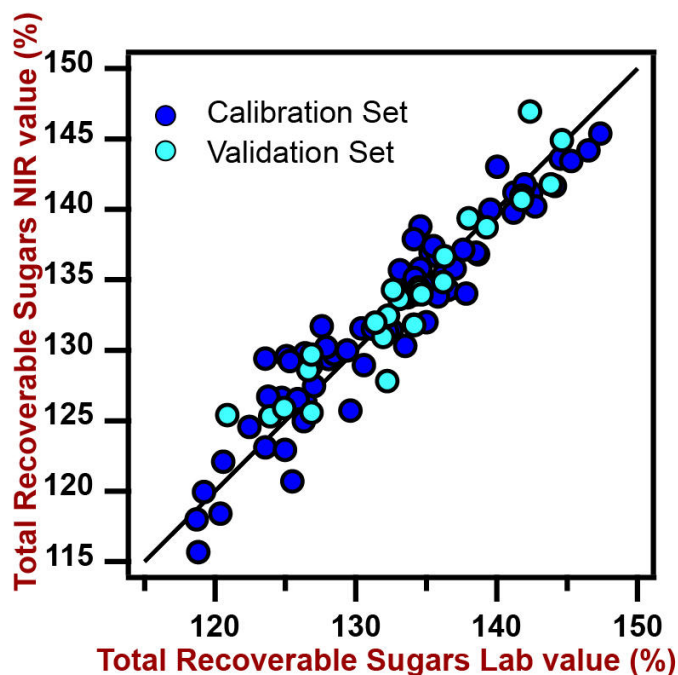


Figure 8. Correlation diagram and the respective FOMs for the prediction of total recoverable sugars in sugarcane juice using a DS2500 Liquid Analyzer. Laboratory values were evaluated using Pol and reducing sugars values: $TRS = (9.5263 \times Pol) - (9.05 \times RS)$.

Figures of Merit	Value
R^2	0.9463
Standard Error of Calibration	2.2985%
Standard Error of Cross-Validation	2.5118%
Standard Error of Validation	1.9074%

CONCLUSION

This Application Note demonstrates the feasibility to determine Brix, Pol, juice purity, reducing sugars, and total recoverable sugars in sugarcane juice with NIR spectroscopy. Vis-NIR

spectroscopy enables a fast and highly accurate alternative to other standard methods (Table 2). No sample preparation is required, and results are delivered in less than a minute.

Table 2. Time to result overview for the different quality control parameters typically measured in sugarcane juice.

Parameter	Method	Time to result
Brix	Refractometer	~1 min
Pol	Calculated from Pol and Brix, as well as the application of a few constants	~10 min sample preparation (clarification & filtration) + ~1 min polarimeter + ~1 min refractometer
Purity	Calculated from Pol and Brix	$\text{Purity} = 100 \times (\text{Pol}/\text{Brix})$
Reducing sugars (RS)	Ion Chromatography	10 min sample preparation (clarification & filtration) + ~40 min IC
Sucrose reading	Polarimeter	~10 min sample preparation (clarification & filtration) + ~1 min polarimeter
Total recoverable sugars (TRS)	Calculated from Pol and reducing sugars	$\text{TRS} = (9.5263 \times \text{Pol}) - (9.05 \times \text{RS})$

CONCLUSION

Internal reference: AW NIR CH-0073-042023

CONTACT

メトロームジャパン株式会社
 143-0006 東京都大田区平和島6-1-1
 null 東京流通センター アネックス9階

metrohm.jp@metrohm.jp

CONFIGURATION



DS2500 Liquid Analyzer

ラボおよび生産環境における品質管理用の堅牢な近赤外分光法。

NIRS DS2500 Liquid Analyzerは、生産チェーン全体に沿った液体のルーチン分析に実績のあるフレキシブルなソリューションです。頑丈な仕様により、NIRS DS2500 Liquid Analyzerは粉塵、湿気、および振動に強い為、過酷な生産環境での使用に理想的です。

DS2500 Liquid Analyzer は400~2500 nmのスペクトル範囲全体をカバーし、サンプルを80° Cまで加熱し、様々な使い捨てハイアルやクォーツキュベットとの互換性を有します。このようにして個々のサンプル要件に対応し得るDS2500 Liquid Analyzerは、1分未満で正確かつ再現性ある結果を得られるようサポートします。さらに一体型のサンプルホルダー検出、および説明不要のVision Airソフトウェアを用いることで、ユーザーが簡単かつ安全に操作できることが保証されます。

サンプル量が多い場合、Metrohmサンプルロケットと組み合わされたフローセルの導入によって、生産性を著しく向上させることができます。