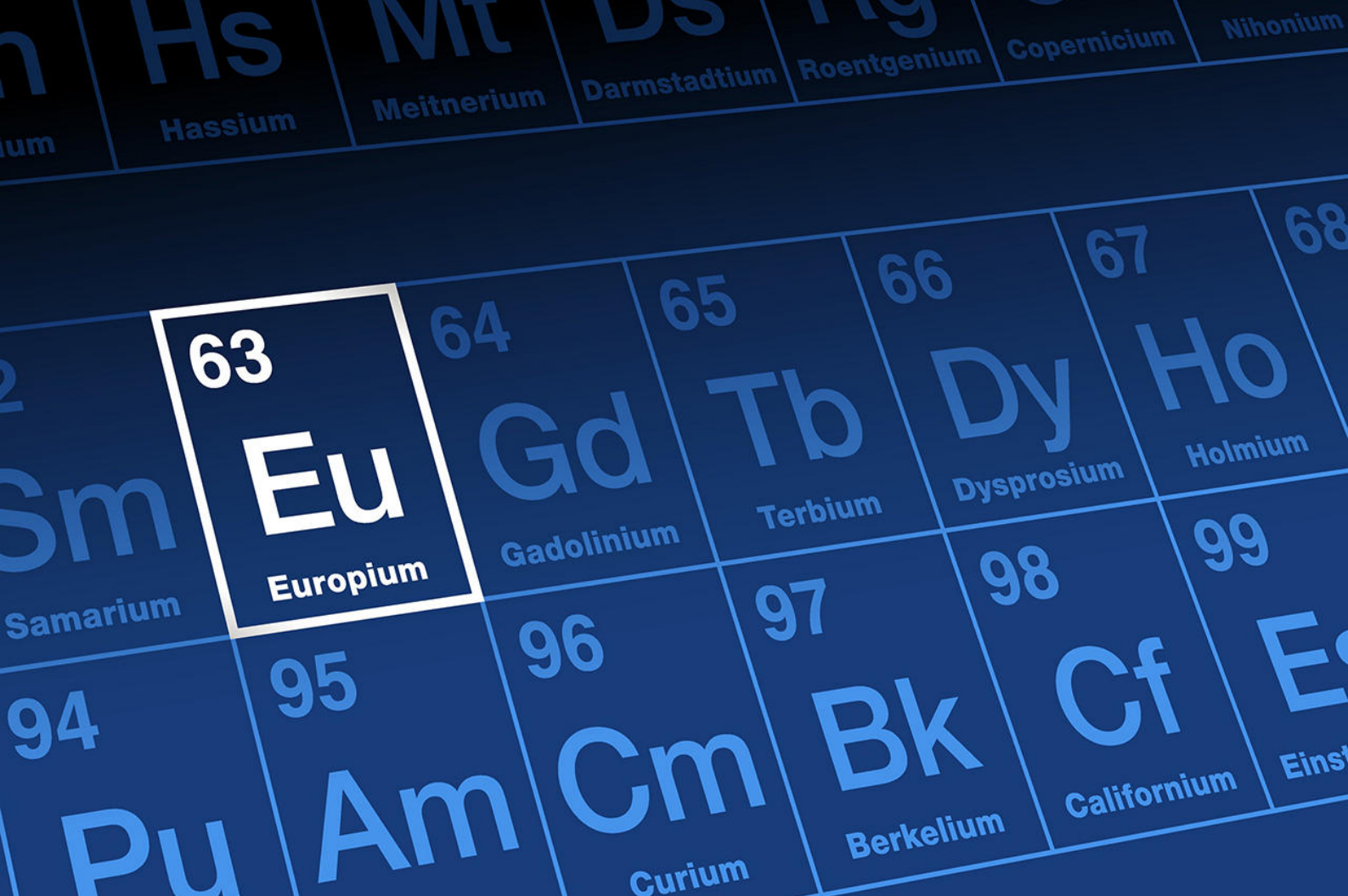
Application Note AN-T-250

Análise potenciométrica de elementos de terras raras (REEs)

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metrohm@metrohm.com.br

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OMNIS Sample Robot S – WSM (1T/2P)

OMNIS Sample Robot S – WSM, equipped with an OMNIS Workstation Module with 2 pumps for cleaning and extraction of the sensors and sample vessels, a workstation, rod stirrers, and extensive accessories for getting started directly with fully automated titration. The system provides space in two sample racks for 32 sample beakers of 120 mL each. This modular system is supplied completely installed and can thus be put into operation in a very short time. The system can also be extended upon request to include 2 additional peristaltic pumps and another workstation, thus doubling the throughput. If additional workstations are required, then the Sample Robot S can be extended to become an L-sized OMNIS Sample Robot, thus enabling samples from 7 racks to be processed in parallel on up to 4 workstations to quadruple the sample throughput.

OMNIS Professional Titrator without stirrer

Innovative, modular potentiometric OMNIS Titrator for endpoint titration and equivalence point titration (monotonic/dynamic). Thanks to 3S Liquid Adapter technology, handling chemicals is safer than ever before. The titrator can be freely configured with measuring modules and cylinder units and can have a stirrer added as needed. Including "Professional" function license for parallel titration with additional titration or dosing modules.

- Actuation via PC or local network
- Connection option for up to four additional titration or dosing modules for additional applications or auxiliary solutions
- Can be supplemented with magnetic stirrer and/or rod stirrer
- Various cylinder sizes available: 5, 10, 20 or 50 mL
- Liquid Adapter with 3S technology: Safe handling of chemicals, automatic transfer of the original reagent data from the manufacturer

Measuring modes and software options:

- Endpoint titration: "Basic" function license
- Endpoint and equivalence point titration (monotonic/dynamic): "Advanced" function license
- Endpoint and equivalence point titration (monotonic/dynamic) with 5-way parallel titration: "Professional" function license

OMNIS Dosing Module without stirrer

Dosing module for connection to an OMNIS Titrator for extending the system to include an additional buret for titration/dosing. Can be supplemented with one magnetic stirrer or rod stirrer for use as separate titration stand. Freely selectable cylinder unit with 5, 10, 20 or 50 mL.

Ion-selective electrode, Cu

Copper-selective electrode with crystal membrane.

This ISE has to be used in combination with a reference electrode and is suitable for:

- ion measurements of Cu^{2+}
- ion measurements in small sample volumes (minimum immersion depth 1 mm)
- Complexometric titrations with CuEDTA

Thanks to the robust/break-proof plastic shaft made of EP, this sensor is mechanically very resistant.

The polishing set supplied enables easy cleaning and renewing of the electrode surface.