



Application Note 410000006-B

# Low Frequency Raman Spectroscopy

Raman spectroscopy is an analytical tool that permits measurement of the molecular structure and identification of materials based on the rotational and vibrational modes of a molecule. Most commercial laboratory Raman systems cover the fingerprint spectral region ranging from 200–3400  $\text{cm}^{-1}$ .

The B&W Tek i-Raman Plus BAC102 probe can access

lower-frequency modes down to 65  $\text{cm}^{-1}$ , providing a cost-effective solution for fuller range measurements. Access to lower-frequency regions provides key information for applications in protein characterization [1], polymorph detection and identification [2], along with material phase and structure determination.

## INTRODUCTION

The low-frequency region augments the information content from the fingerprint region of the Raman spectrum and broadens possible molecular

applications, such as detection of hydrogen bonds. The result is increased detection sensitivity and differentiation of very similar materials.


## Amino acids

Raman spectroscopy is a modern way to study the structure and conformation of the building blocks of proteins—amino acids. Specifically, the vibrational information content in a Raman spectrum can help with interpretation of molecular interactions and biological processes [3].

Unlike many substances that do not exhibit peaks below  $\sim 400\text{ cm}^{-1}$ , the low-frequency portion of the

Raman spectrum is a necessary source of information for a comprehensive study of amino acids. This is obvious in the full Raman spectrum of L-asparagine, from  $65\text{--}3200\text{ cm}^{-1}$  (Figure 1).

Figure 1 shows both the fingerprint region (blue) as well as the low-frequency Raman region (red) for L-asparagine; note the three dominant bands below  $200\text{ cm}^{-1}$ .

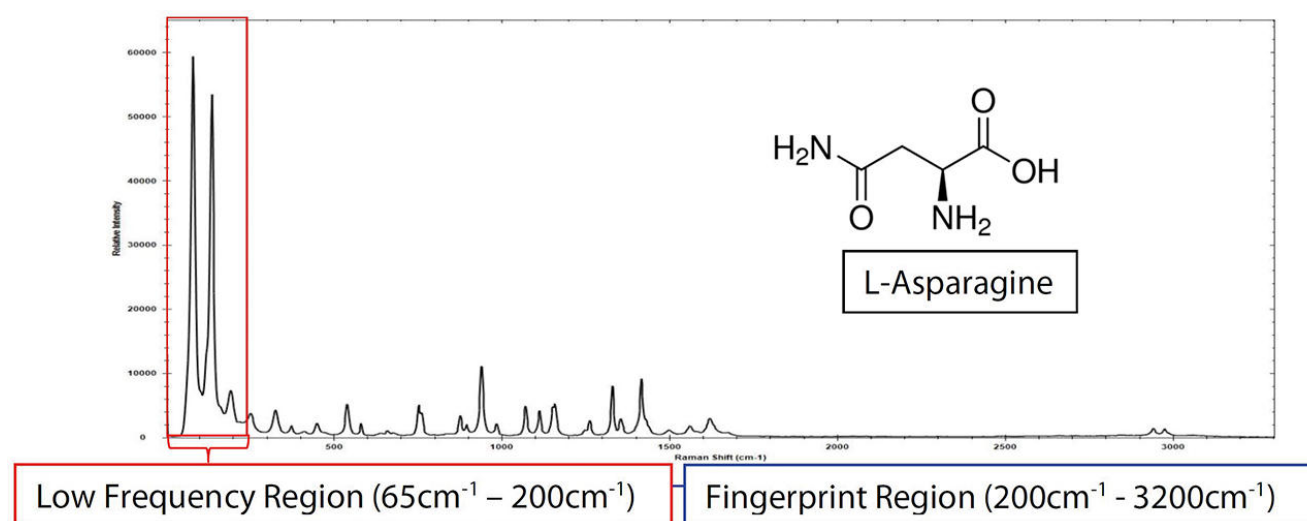

**Figure 1.** i-Raman Plus and a BAC102 E-grade probe were used to collect the low-frequency spectra of L-asparagine with a total integration time of 1.2 s.

## EXPERIMENT

i-Raman Plus 785S, B&W Tek's laboratory Raman spectrometer that utilizes patented CleanLaze®, was used in this application. This instrument features 785 nm laser excitation with a linewidth of less than 0.2 nm and maximum power output of 300 mW.

i-Raman Plus is equipped with a sensitive, TE-cooled, back-thinned CCD. A BAC102 E-grade probe, utilizing

proprietary technology, supports data collection within a full spectral range of  $65\text{--}33500\text{ cm}^{-1}$  with a spectral resolution of  $4.5\text{ cm}^{-1}$ .

Raman spectra were collected at room temperature using 300 mW laser power with integration times that ranged from 100 milliseconds to 10 seconds (Table 1).

**Table 1.** Experimental parameters.


Equipment	Acquisition settings	
i-Raman Plus 785S	Laser Power	300 mW
BAC102 probe	Integration time	1.2 s
BWSpec software	Averages	1

Determining the structural form of active pharmaceutical ingredients (APIs) is a primary concern for the pharmaceutical industry. This is especially true during drug development, manufacturing, and quality control.

APIs exhibit polymorphism—identical chemical composition but different solid-state structures. Polymorphs may affect bioavailability and therapeutic

index. The efficacy of a drug can be compromised if the wrong form is used [2]. Pseudo-polymorphs include solvents suspended in a lattice structure.

**Figure 2** is an example of the pseudo-polymorph D-glucose, demonstrating the ability of the E-grade probe to detect differences between the monohydrate and anhydrous forms at frequencies below  $200\text{ cm}^{-1}$ .

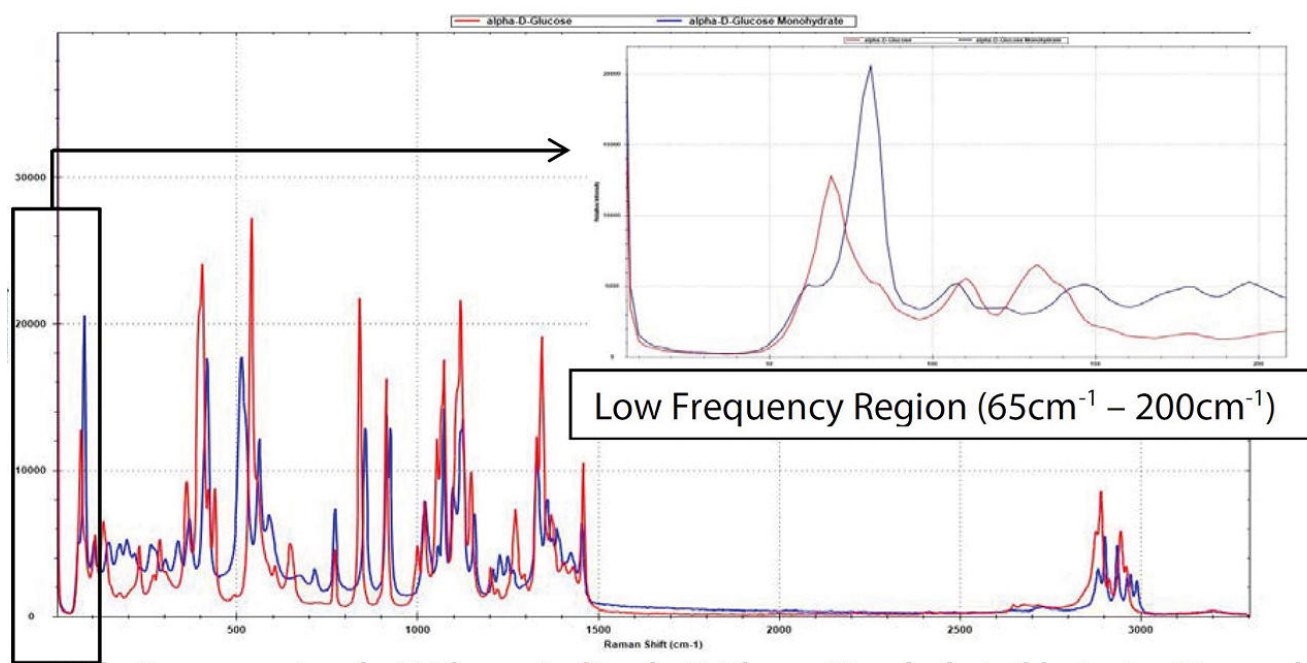


**Figure 2.** Raman spectra of  $\alpha$ -D-glucose (red) and  $\alpha$ -D-glucose monohydrate (blue), taken with 10 s integration time. Note the significant difference between the two pseudo-polymorphs within the low-frequency range (see inset).

Exceptional specificity is required to monitor phase changes like crystallization in chemical processes. The low-frequency E-grade probe can monitor such phase changes, as demonstrated for sulfur (Figure 3).

Solid  $\alpha$ -sulfur was deposited onto an aluminum tray and heated with a hot plate while the Raman spectra were collected with an E-grade probe and i-Raman Plus, using 100% laser power (~300 mW) and 0.1 s

integration time for both the solid and liquid phase.

After the sample was heated above its melting point at 115.2 °C, the low-frequency peak at 83.6  $\text{cm}^{-1}$  broadened and shifted, indicating the change from the  $\alpha$ - to the  $\lambda$ -form. Note that there are no observable changes within the fingerprint region (Figure 3).

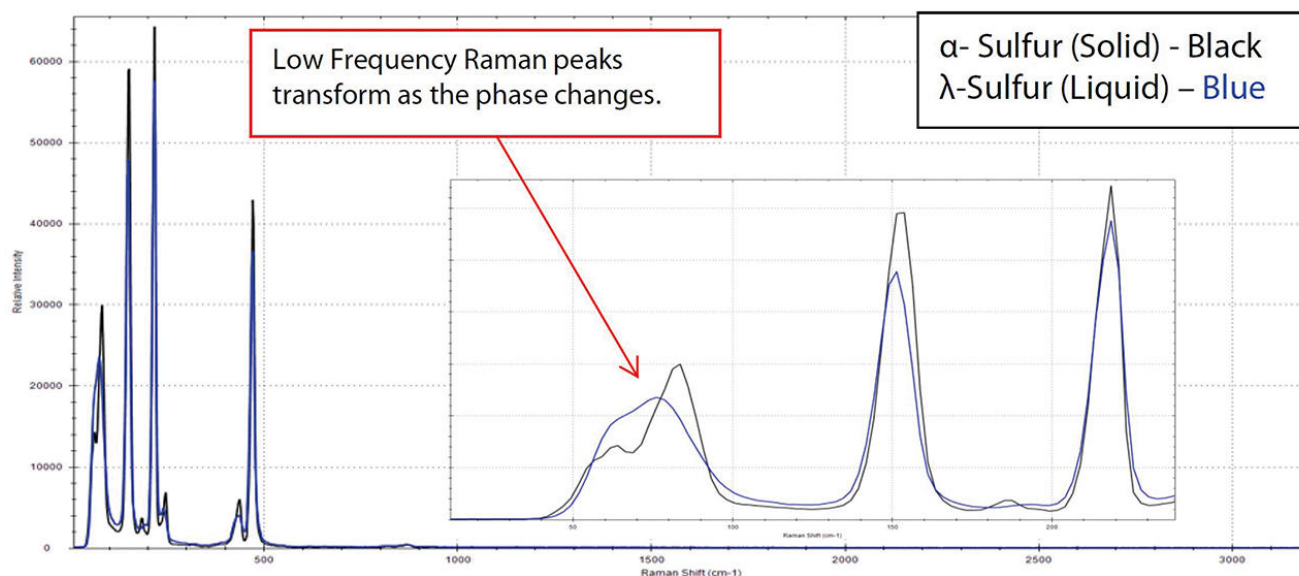

**Figure 3.** Raman spectra of sulfur transitioning from the  $\alpha$ -crystalline form to the  $\lambda$ -liquid form, taken with 0.1 s integration time. Note significant broadening in the peaks located in the low-frequency region (see inset).

## CONCLUSION

i-Raman Plus 785S Raman spectrometer, coupled with the low-frequency E-grade probe, can be a valuable tool for applications requiring low-frequency detection down to 65  $\text{cm}^{-1}$ . The ability to characterize polymorphs and solvated forms supports manufacturing and formulation processes in the

pharmaceutical and biological industries.

Along with protein, polymorph, and phase characterization, low-frequency Raman spectroscopy can also be used to study semiconductor lattices [4], carbon nanotubes [5], solar cells, and an assortment of minerals, pigments and gemstones.

## REFERENCES

1. Teixeira, A. M. R.; Freire, P. T. C.; Moreno, A. J. D.; et al. High-Pressure Raman Study of L-Alanine Crystal. *Solid State Communications* **2000**, *116* (7), 405–409.  
[https://doi.org/10.1016/S0038-1098\(00\)00342-2](https://doi.org/10.1016/S0038-1098(00)00342-2).
2. Larkin, P. J.; Dabros, M.; Sarsfield, B.; et al. Polymorph Characterization of Active Pharmaceutical Ingredients (APIs) Using Low-Frequency Raman Spectroscopy. *Appl Spectrosc* **2014**, *68* (7), 758–776.  
<https://doi.org/10.1366/13-07329>.
3. Golichenko, B. O.; Naseka, V. M.; Strelchuk, V. V.; et al. Raman Study of L-Asparagine and L-Glutamine Molecules Adsorbed on Aluminum Films in a Wide Frequency Range. *Semicond. Phys. Quantum Electron. Optoelectron.* **2017**, *20* (3), 297–304.  
<https://doi.org/10.15407/spqeo20.03.297>.
4. Smith, E.; Dent, G. *Modern Raman Spectroscopy: A Practical Approach*, 2nd ed.; John Wiley & Sons, 2019.
5. Pelletier, M. J. *Analytical Applications of Raman Spectroscopy*, 1st ed.; Blackwell Science: Oxford, 1999.

## CONTACT

Metrohm Portugal  
R. Frei Luis de Granada 14G  
1500-680 Lisboa

[vendas@metrohm.pt](mailto:vendas@metrohm.pt)