



Application Note AN-BAT-006

MacMullin数的定 加方法

子池的一个重要部件是隔膜,用于物理分正。通常隔膜由聚合物材料制成。

为了量隔膜的性能,引入了MacMullin数。MacMullin数解的子率与有隔膜解的子率的比。

本用告中,MacMullin数是通加方法算的。方法包括具有不同堆厚度的品行化学阻抗量,通增加池内隔膜的数量来。然后根据数据合算子阻,并制出与池中隔膜数量的系,斜率出了MacMullin数。

INTRODUCTION

One important component of lithium ion batteries is the separator foil – a thin, electronically insulating, porous membrane that prevents the electrodes from being in direct contact with each other while being permeable to ionic current. To maintain a good lithium ion

battery performance, the resistance of a separator foil filled with electrolyte should be low. Consequently, information about the effective ion conductivity of these systems is an important parameter for the performance of a separator material.

Usually, the focus is not on the effective ionic conductivity, but on another quantity: the so-called MacMullin number N_M [1]. The MacMullin number is the ratio of the ionic conductivity of

the pure electrolyte $\sigma_{electrolyte}$ and ionic conductivity of the separator foil filled with electrolyte $\sigma_{separator}$ (Equation 1).

$$N_M = \frac{\sigma_{electrolyte}}{\sigma_{separator}} \quad 1$$

Therefore, the MacMullin number is a quantity that describes the decrease of the effective conductivity by the presence of a separator foil. For lithium ion battery separators, often a value

between 4 and 20 is found [1–3]. Theoretically, the MacMullin number is linked to the separator foil's porosity ϵ and tortuosity τ [2, 3].

EXPERIMENTAL

Chemicals

As electrolyte, a 1 mol/L LiPF_6 (lithium hexafluorophosphate) solution in a 1:1 ratio in volume mixture of EC (ethylene carbonate) and DMC (dimethyl carbonate, purchased from Sigma-Aldrich / Merck KGaA) was used without any further purification. After receiving, the

electrolyte was stored and handled inside of an argon-filled glove box. As battery separator material, a tri-layer separator based on polyethylene/polypropylene with a thickness of 21.5 μm was used.

Sample preparation and measuring setup

Circular specimens with a diameter of 12 mm were punched out from the separator foil. To ensure an optimal wetting, the specimens were stored for at least 24 hours in the electrolyte solution.

For the electrochemical measurements, a TSC battery standard cell (Figure 1) in combination with a Microcell HC setup was used.

As current collectors, two planar stainless steel disc electrodes (diameter 8 mm) press-fitted into a PEEK sleeve were used. The cell constant is the ratio of the electrode area A (cm^2) and its thickness l (cm), and it is measured in cm^{-1} . To calculate the cell constant, the area of the stainless steel electrodes and the separator thickness were measured, resulting in a cell constant value of 0.0043 cm^{-1} .

The cell constant, together with the resistance related to the bulk ion transport inside of the porous separator network R_{ion} (Ω), is used in the calculation of the conductivity σ (S cm^{-1} , Equation 2).

$$\sigma = \frac{1}{R_{ion}} \cdot \frac{l}{A}$$

2

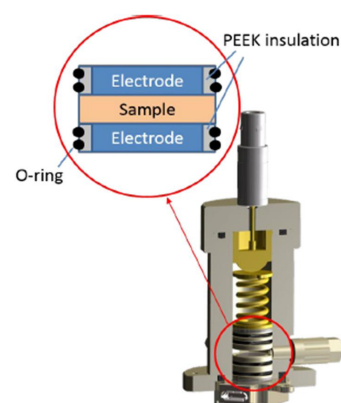


Figure 1. Schematic drawing of the measuring cell TSC battery standard. The separator foils filled with electrolyte were placed between the upper and lower electrode.

The contact pressure applied onto the electrolytesoaked separator foil stack inside of the TSC battery cell was adjusted to approximately 100 kPa using a gold-plated spring with a spring constant of 2.43 N mm^{-1} . At this contact pressure, no significant compression of the separator foil was found.

In the Microcell HC setup, the temperature is adjusted via a Peltier element, with an accessible temperature range of $-40 \text{ }^\circ\text{C}$ to $+100 \text{ }^\circ\text{C}$. To

accurately measure the temperature, a Pt100 temperature sensor is embedded in the base unit of the TSC battery cell, at a position very close to the sample. The accuracy of the temperature measurement is $0.1 \text{ }^\circ\text{C}$ with respect to the sensor position inside of the base unit. For the determination of the MacMullin number, the temperature was adjusted to $20.0 \text{ }^\circ\text{C}$.

Because the TSC battery measuring cell is

completely airtight, the experiments took place on the Microcell system outside of the glove box. For the impedance measurements, a Metrohm Autolab PGSTAT204 equipped with a FRA32M module was used. The communication with the

Measurement parameters

The impedance measurements were performed at a frequency range of 1 MHz to 100 Hz with an amplitude of 10 mV.

Impedance measurements were taken on samples with different stack thickness, achieved

temperature controller is integrated in the NOVA software enabling automated measuring routines.

by increasing the number of separators inside the cell.

After reaching the temperature set point of 20.0 ° C, a waiting time of 300 s was chosen to ensure for complete thermal equilibrium.

Table 1. Steps and relative actions per step.

Step	Action to be performed
1	Transferring the cell into the glove box and loading the cell with one specimen of the electrolyte-filled separator foil.
2	Taking the prepared cell out of the glove box and connecting it to the Microcell HC cell stand.
3	Connecting the measuring device (2-electrode configuration).
4	Setting the temperature to 20.0 ° C and waiting for 300 s.
5	Performing an impedance spectroscopy measurement.
6	After finishing step 5, repeating steps 1–5 by adding further specimens of the electrolyte-filled separator foil (one per loop) until a stack of five specimens is created.

RESULTS

The resulting impedance spectra for different stack thicknesses are shown in **Figure 2**.

At high frequencies, the impedance behavior is dominated by inductive effects caused by the connection cables. The intersection at the Z' -axis is close to the bulk resistance value for bulk ion transport inside of the porous separator network. In absence of the inductive behavior at high frequencies, they would be identical. The increase at lower frequencies is then caused by electrode polarization.

As expected, increasing the stack thickness by adding further specimens of the electrolyte-soaked separator foil leads to a shift to higher values for the intersection point due to the higher bulk resistance value.

The spectra were fitted with the equivalent circuit shown in **Figure 3**.

The inductance L_{cable} takes into account the high frequency contributions by the electrical cables and connectors, the resistor R_{ion} represents the ion migration in the porous separator network, and the constant phase element CPE_{pol} describes the electrode polarization at low frequencies.

The resulting stack thickness-dependent values for R_{ion} are given in **Table 2**. The fit error of R_{ion} was lower than 0.4% in all cases.

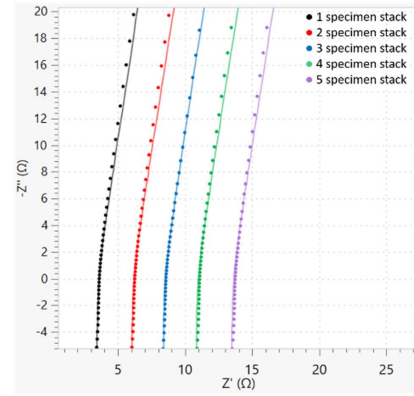


Figure 2. Impedance spectra of the soaked separator foils for different stack thicknesses. For a better visibility of the shift along the Z' -axis, the spectra were cut at 3 kHz. The points are the measured values and the lines represent the fit.



Figure 3. Equivalent circuit for fitting the impedance spectra.

Table 2. Ion migration resistance per sample number.

$N_{specimen}$	R_{ion}/Ω
1	3.38
2	5.98
3	8.34
4	10.78
5	13.43

A plot of these values shows a linear relationship, where the slope gives a value of 2.49Ω for the ionic resistance change per added specimen $\Delta R_{ion} / \Delta N_{specimen}$ (Figure 4).

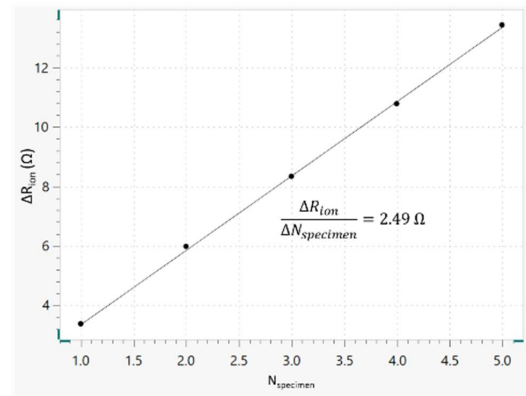


Figure 4. Plot of total ionic resistance value R_{ion} in dependence of the number of separator specimens $N_{specimen}$.

Taking into account the separator thickness d , the active area A of the stainless steel electrode, and the resistance per added specimen $\Delta R_{ion} / \Delta N_{specimen}$, the corresponding ionic

conductivity of the electrolyte-soaked separator $\sigma_{separator}$ can be calculated as shown in Equation 3.

$$\sigma_{separator} = \frac{1}{\Delta R_{ion} / \Delta N_{specimen}} \cdot \frac{d}{a} \quad 3$$

Here, a value of $\sigma_{separator} = 1.7 \text{ mS cm}^{-1}$ was measured at 20° C . The conductivity of the pure electrolyte solution at 20.0° C was determined to be $\sigma_{electrolyte} = 9.9 \text{ mS cm}^{-1}$. Therefore, the MacMullin number was calculated to be $N_M = 5.8$ which lies within the typical range for such separators [3]. It should be mentioned that Raccichini et al. found that the MacMullin number can be slightly dependent on the chosen electrolyte solution [4], although Landesfeind et al. showed

that there was no dependence of the MacMullin number on the electrolyte salt, solvent composition, or salt concentration for the samples chosen by them [3]. Therefore, further research has to be performed to shed more light on the influence of the electrolyte properties on the determined MacMullin number values for separator foil. Ideally, the MacMullin number should be determined by using the electrolyte that is used for the final application.

CONCLUSION

In this application note, the so-called stacking method is applied to determine the MacMullin number of separator foils. A Metrohm Autolab Microcell HC setup combined with a TSC battery

standard cell is used, together with a typical lithiumion battery electrolyte as well as a typical separator material.

ACKNOWLEDGEMENT

The results presented here are part of a sub-project related to the project LiMES funded by the Federal Ministry for Economic Affairs and

Energy as part of the 7. Energieforschungsprogramm which is highly appreciated.

BIBLIOGRAPHY

1. P. Arora, Z. Zhang, Chem. Rev. 2004, 104, 4419- 4462
2. M. J. Martínez, S. Shimpalee, J. W. Van Zee, J. Electrochem. Soc. 2009, 156 (1), B80-B85.
3. J. Landesfeind, J. Hattendorff, A. Ehrl, W. A. Wall, H. A. Gasteiger, J. Electrochem. Soc. 2016, 163 (7), A1373-A1387.
4. R. Raccichini, L. Furness, J. W. Dibden, J. R. Owen, N. García-Araez, J. Electrochem. Soc. 2018, 165 (11), A2741-A2749.

CONTACT

瑞士万通中国
北京市海淀区上地路1号院
1号楼7702
100085 北京

marketing@metrohm.com.cn

CONFIGURATION



Autolab PGSTAT204

PGSTAT204 合了小巧格和模化。器包括基本恒位/恒流,其从 20 V,最大流 400 mA 或 10 A,与 BOOSTER10A 合使用。此恒位可随用附加模行展,例如 FRA32M 化学阻抗(EIS)模。

PGSTAT204 是一款惠的器,可置于室的任何位置。具有模和数字入/出,可控制 Autolab 附件和外部。PGSTAT204 包括内置模分器。与高性能的 NOVA 件用,可用于大多数准化学技。




TSC Battery Standard

TSC Battery 是一封式量池,用于装 **半池** 或 **全池** 或 **有活性物和隔板的超大容量池**。此外,它用于量与板/板几何形状的金属接触的固体和凝状品。

池的配置 **具有活性**,可以 **通面孔** 插入微参比来行个和三个量。

准版本包含扁平金属。有使用元素作材料的集流体,参展版本。”

此化学池是密封的,因此与大气隔,可用于量湿度敏感的物品。

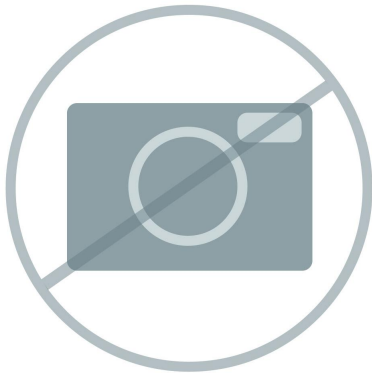

NOVA 是通 USB 接口控制所有 Autolab 器的件包。由化学家化学而,集成了超二十余年的用体和最新的 .NET 件技,NOVA 使的 Autolab 恒位/恒流有更性能和活性。

NOVA 提供了以下的独特功能:

- 功能大且活的程序器
- 重要数据一目了然
- 大的数据分析和工具
- 集成化控制外器,如万通 LQH 液体理