



Application Note AN-RS-017

## 香蕉上的痕量

Protecting consumer safety with Misa

### 概要

Thiabendazole (TBZ) is a broad-spectrum pesticide used both as a fungicide in fruits and vegetables and for controlling parasites in animal feed. While some individuals experience negative side effects from consuming foods which contain TBZ, overall toxicity to humans is very low. To ensure consumer safety, regulatory agencies establish maximum residue levels (MRL) for pesticide-treated crops based on their review of risk assessment studies. For bananas, which are either aerially sprayed or dipped in

protectant solutions of TBZ, the US FDA reports a MRL of 3  $\mu\text{g/g}$ , and the EU stipulates a MRL of 6  $\mu\text{g/g}$  by weight.

With Misa (Metrohm Instant SERS Analyzer), the rapid and sensitive detection of TBZ on bananas is demonstrated in formats easily adapted for food safety surveillance testing. SERS substrates for Misa offer unparalleled sampling flexibility—from dipping to swabbing, P-SERS can be adapted to the desired application.


This application note describes the simple extraction and detection of TBZ residue on banana peels and in a processed banana

product. Samples are applied to Silver P-SERS (Ag P-SERS) strips through both **dipping** and **swab methods**.

### REFERENCE SPECTRUM AND LIBRARY CREATION

To establish a reference spectrum for TBZ, a pure standard in methanol (10 µg/mL) was analyzed using Ag P-SERS. The unique SERS spectrum

shown in **Figure 1** can be used to create a library entry for TBZ.

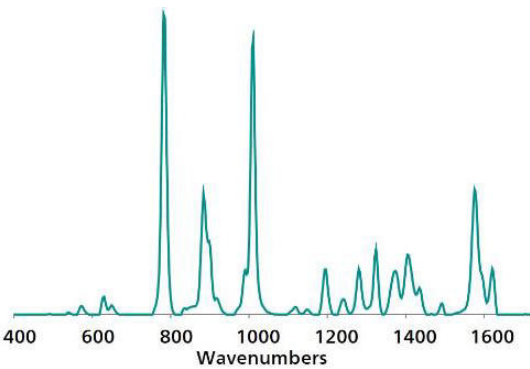

**Figure 1.** Standard reference Ag P-SERS spectrum of Thiabendazole.


**Table 1.** Experimental Parameters

Instrument		Acquisition	
Firmware	0.9.33	Laser Power	5
Software	Misa Cal V1.0.15	Int. Time	1 s
Misa Vial Attachment	6.07505.030	Averages	10
ID Kit - Ag P-SERS	6.07506.470	Raster	ON

## EXPERIMENT AND RESULTS

### Detection of TBZ in banana chips

Dried banana chips purchased from a local grocer were ground to a fine powder and spiked with TBZ in methanol to provide a range of test samples: 100–0.1  $\mu\text{g/g}$ . Samples were mixed well and dried. Afterward, 0.1 g of each sample was added to a glass vial with 1 mL of 0.1 mol/L NaOH and shaken. After settling, an Ag P-SERS strip was **dipped** into each extract, dried for 2 minutes, and immediately placed into the Misa P-SERS attachment for measurement. **Figure 2** shows the range of detection of TBZ in ground dried banana samples from 100  $\mu\text{g/g}$  to 100 ng/g.

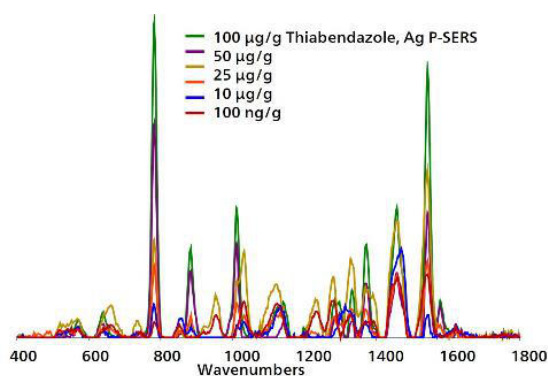


**Figure 2.** Overlaid baseline-corrected, dipped Ag P-SERS spectra for pulverized banana chips spiked with TBZ show detection down to 100 ng/g.

### Detection of TBZ on banana peel

To simulate Ag P-SERS **swab testing**, locally purchased organic bananas were washed with tap water and isopropanol/distilled water. To simulate typical application of TBZ, a 2  $\text{cm}^2$  section of peel was treated with 20  $\mu\text{L}$  of methanol containing different concentrations

(250, 100, 75, 50, 25, and 10  $\mu\text{g/mL}$ ) of TBZ. After drying, an Ag P-SERS strip wetted with methanol was used to swab the section. The strip was dried for 1 minute and placed into the Misa P-SERS attachment for measurement, with resulting spectra shown in **Figure 3**.

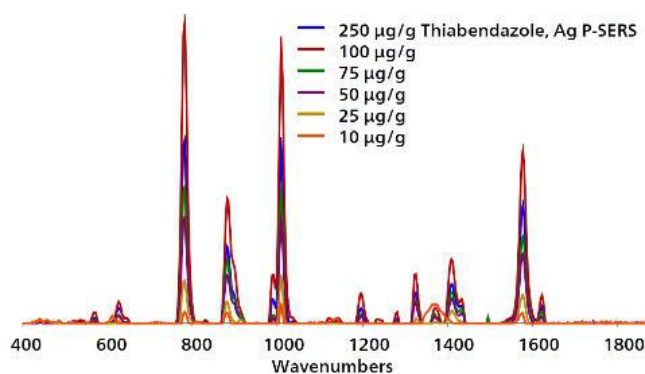


**Figure 3.** Overlaid baseline-corrected spectra acquired from Ag P-SERS swabs show detection of TBZ on banana peel to 100 ng/g.

### Direct test of a store bought banana

In a final test, a non-organically grown unwashed banana from the grocers was **swabbed** with a methanol-moistened Ag P-SERS strip. In **Figure 4**,

detection of TBZ is confirmed by peak agreement between the spectrum from the banana swab and that of the TBZ standard.

## EXPERIMENT AND RESULTS

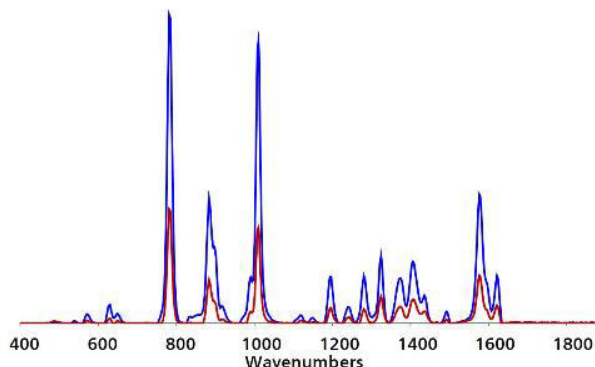

**Figure 4.** Direct testing of an untreated non-organic banana purchased from a grocery store returned an Ag P-SERS spectrum consistent with detection of TBZ.

### Detection of thiabendazole in the field

Using a disposable pipette, wet the printed end of a P-SERS strip with 2 drops of methanol. *Touch only the unprinted back side of the swab with your gloved hands.* Use the moistened, printed portion of the strip to swab the peel of a banana. Move the strip

along the peel for at least 2 cm. After drying for 2 minutes, insert the Ag P-SERS strip, with the printed portion facing down and toward the instrument, into the P-SERS attachment for measurement.

**Table 2.** Requirements for Field Test Protocol

<b>ID Kit - Ag P-SERS</b>	<b>6.07506.470</b>
includes:	Silver P-SERS
	Scoop
	Disposable pipettes
	2 mL glass vials
<b>Reagents</b>	Methanol
<b>Test settings</b>	Use <b>ID Kit OP</b> on MISA

## CONCLUSION

After Misa's ability to sensitively detect TBZ on bananas is demonstrated, a successive experiment confirms post-harvest treatment of bananas with

TBZ. Such fast and accurate test results recommend Misa as an excellent solution for streamlined, onsite identification of chemically treated foods.

## CONTACT

117702  
100085

marketing@metrohm.com.cn

## CONFIGURATION



### MISA Advanced

Metrohm Instant SERS Analyzer (MISA) 是一款高性能、便携式分析系,可快速/定非法物、食品添加和微量食品染料。MISA 的特点是配了 Metrohm 的道光栅描 (ORS) 技的光。其空需求小和并且池寿命有所延,用于或移室用。MISA 提供各 1 激光附件,可活取。

分析可通 BlueTooth 或 USB 接行。  
MISA Advanced 套件是一个完整套件,其作用是用能用 Metrohms 米粒溶液和 P-SERS 条行 SERS 分析。  
MISA Advanced 套件包含了一个 MISA 小管附件、一个 P-SERS-附件、一个 ASTM 校正准件、一个 USB 迷、一个 USB 供元和用于行 MISA 器的 MISA Cal 件。随供了一个用来安全保管器和附件的固保箱。



### ID - -P-SERS (Ag P-SERS)

ID 套件 - Ag P-SERS 包含了 Mira/Misa 用使用-P-SERS 基材行 SERS 分析所需的件。套件包含了一个一次性抹刀、一个移液管、品小瓶和 2 个-P-SERS 条。