



Application Note AN-CIC-034

通燃子色法(CIC)快速分析水中的可吸附有机素(AOX)

Measurement of AOCl, AOBr, AOI, and AOF according to DIN 38409-59 and ISO/DIS 18127

AOX (adsorbable organically bound halogens) is a complex parameter covering the sum of halogenated organic compounds adsorbable on activated carbon. Many of these organohalogens and their degradation products pose serious risks to human health and the environment [1–4]. Monitoring them is essential to ensure appropriate water quality, to trace their sources, or to investigate the efficiency of AOX removal techniques in water treatment processes. Historically, AOX was determined via microcoulometric titration after adsorption of water samples on activated carbon and

subsequent combustion (DIN EN ISO 9562 or EPA 1650) [1,2]. By definition, based on the technical setup, AOX was comprised of adsorbable organically bound chlorine (AOCl), bromine (AOBr), and iodine (AOI)—but not fluorine (AOF)—as a sum parameter and not its individual fractions. Both **DIN 38409-59** and **ISO/DIS 18127** describe a validated procedure of adsorption and analysis via **combustion ion chromatography (CIC)** to determine AOCl, AOBr, AOI, AOF, and the sum parameter **CIC – AOX_(Cl)**. This Application Note explains the CIC method used to fulfill the requirements of these

EXPERIMENTAL

This application is focused on the experimental approach of AOX analysis. More detailed information can be found in related Metrohm literature ([WP-081](#), [AN-CIC-033](#) – specifically about AOF). The complete validation dataset of DIN 38409-59 is available on the [Water Chemistry Society webpage](#).

The overall sample preparation procedure, i.e., preconcentration and adsorption of organically bound halogens, resembles that of DIN EN ISO 9562, as adsorption on activated carbon is a key point for both methods (**Figure 1**). While for AOF it is crucial that the samples are neutral to avoid adsorption of inorganic fluorine to the activated carbon, sample acidification is mandatory for the other organically bound halogens, similar to DIN EN ISO 9562. For CIC-AOX_(Cl) determination (i.e., AOCl, AOBr, and AOI), samples need to be acidified with nitric acid to pH <2 prior to preconcentration (**Table 1**).

The adsorption of the organically bound halogens is handled in a semi-automated manner using the APU sim system from Analytik Jena (**Figure 1**). Two columns filled with activated carbon (at least 50 mg in each column) are connected in series and 100 mL of sample is passed through. The organically bound halogens adsorb to the activated carbon (using dedicated disposable columns for AOF and AOX

determination, **Table 1**), while inorganic halogens are rinsed off (**Figure 1**).

After the semi-automated sample preparation is finished, the complete content of the two adsorption columns is transferred into one or two separate ceramic boats for CIC analysis. Combustion occurs at temperatures above 950 °C in the presence of argon and oxygen (**Figure 1**). For pyrohydrolytic combustion, a water stream is essential as it converts the halogens into their hydrogenous forms. Chlorine, bromine, iodine, and fluorine are volatilized in the combustion step, transported into the absorber solution (ultrapure water) with an argon/oxygen gas stream, and transferred into the liquid phase (**Figure 1**). Dosinos guarantee precise automated liquid handling, e.g., the transfer of the aqueous sample into the IC for analysis, or the water stream essential for pyrohydrolytic combustion.

The ion chromatographic separation is achieved on a Metrosep A Supp 5 - 250/4.0 column in combination with the A Supp 5 Guard/4.0. AOF (as F) elutes in under 7 minutes while AOX (i.e., Br, Cl, and I) elutes in less than 25 minutes (**Figure 2**). Automatic system calibration with [MiPT](#) (Metrohm intelligent Partial-Loop Injection Technique) is performed using inorganic anion standards for fluoride, chloride, bromide, and iodide (1 g/L

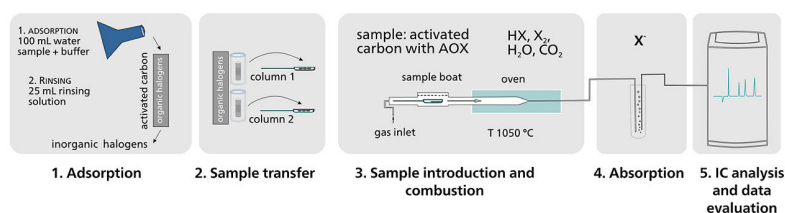


Figure 1. Schematic of the procedure for AOX (WP-081) and AOF (AN-CIC-033) analysis. The first step is adsorption performed with the APU sim (Analytik Jena) for semi-automated and standardized adsorption of up to six samples in parallel. After the second step of sample transfer into the combustion boats, the sample is automatically combusted (step 3, combustion module from Analytik Jena consisting of a combustion oven with Auto Boat Drive (ABD) and an autosampler (MMS 5000)). In the fourth step, the volatilized halogens are transported to the absorber solution via gas stream (920 Absorber Module). The last step (5) is the automatic analysis of AOBr, AOCl, and AOI, or of AOF with the IC (930 Compact IC Flex) including data evaluation. The complete CIC process is fully automated and controlled by MagIC Net software from Metrohm.

Table 1. Parameters for AOF and AOX (AOCl, AOBr, and AOI) sample preparation.

	AOF	AOCl, AOBr, AOI
pH	Neutralized	Acidified to pH <2 with nitric acid
Buffer	0.5 mL 2 mol/L sodium nitrate	0.5 mL 2 mol/L sodium nitrate, acidified with nitric acid
Sample volume		100 mL
Rinsing solution		25 mL
	0.01 mol/L sodium nitrate	0.01 mol/L sodium nitrate, acidified with nitric acid
Absorption columns	Two activated carbon tubes (disposable, from Analytik Jena)	
	402-880.616	402-880.610
Flow rate APU sim		3 mL/min

Performance checks of AOF and AOX determinations and the standard series for LOD determination (**Table 2**) are accomplished using organic reference standard solutions with varying concentrations (4-fluorobenzoic acid, 4-chlorobenzoic acid, 4-bromobenzoic acid, and 4-iodobenzoic acid), treated in the same way as the

samples.

As the procedure for the determination of AOX is rather complex, dedicated sample boats and charcoal (i.e., fluoride-free materials for AOF, **Table 1**) and blank measurements are essential to guarantee a low background and an appropriate blank correction (**Equation 1**).

RESULTS

Individual concentrations for AOCl, AOBr, AOI, and AOF are calculated according to **Equation 1**. A sum parameter for AOX (CIC-AOX_(Cl)) is calculated using **Equation 2**. However, due to the novelty of this validated approach, CIC-AOX_(Cl) has not yet replaced AOX in water or wastewater regulations.

Dedicated materials and the sensitive analysis of the halogens with suppressed conductivity detection results in low blank values. Blank values were only measurable for fluoride and chloride (**Table 2**). The requirements for DIN 38409-59 are fulfilled—in fact, the overall procedure here is even more sensitive.

RESULTS

$$c(X_{ads}) = \left(c(X^-)_{IC} \cdot \frac{V_{Abs}}{V_{Smpl}} \right) - \left(c(X^-_{BW})_{IC} \cdot \frac{V_{AbsBW}}{V_{SmplBW}} \right)$$

Equation 1.

$c(X_{ads})$	Mass concentration of individual adsorbable organically bound halogens (with X = Cl, Br, I, and F) in g/L
$c(X^-)$	Halogen concentration in the sample's absorption solution in g/L (with X = Cl, Br, I, and F) in g/L
V_{Abs}	Final volume of the absorption solution in L
V_{Smpl}	Volume of the sample that was used for adsorption; always 0.1 L
$c(X^-)_{BW}$	Halogen concentration in the absorption solution of the blank in g/L
V_{AbsBW}	Final volume of the absorption solution of the blank in L
V_{SmplBW}	Volume of the blank solution that was used for adsorption; always 0.1 L

$$c(CIC-AOX_{(Cl)}) = c(AOCl) + c(AOBr) \cdot 0.4437 + c(AOI) \cdot 0.2794$$

Equation 2.

$c(CIC-AOX_{(Cl)})$	Sum concentration of adsorbable organically bound halogens in g/L as mass concentration based on chloride
---------------------	---

During the DIN validation process, several water samples were analyzed from different laboratories

using similar setups (validation report: wasserchemische-gesellschaft.de).

Table 2. Blank, LOD (limit of detection), and DIN scope for the determination of adsorbable organically bound halogens. LODs are determined according to DIN 32645. For AOB_r and AOI, the LODs are determined using the calibration curve as no blank values were found. For AOF and AOCl, the blank method was applied (DIN 32645).

	Blank (g/L)	LOD (DIN 32645) (g/L)	Scope of DIN application (g/L)
AOF	1.1	0.38	2
AOCl	2.6	1.36	10
AOB _r	0	0.24	1
AOI	0	0.47	1

Using IC, it is now possible to not only determine the sum parameter CIC-AOX_(Cl), but also to measure the individual fractions contributing to the AOX contents

(**Figure 2, WP-081**) and to assess AOF (**AN-CIC-033**).

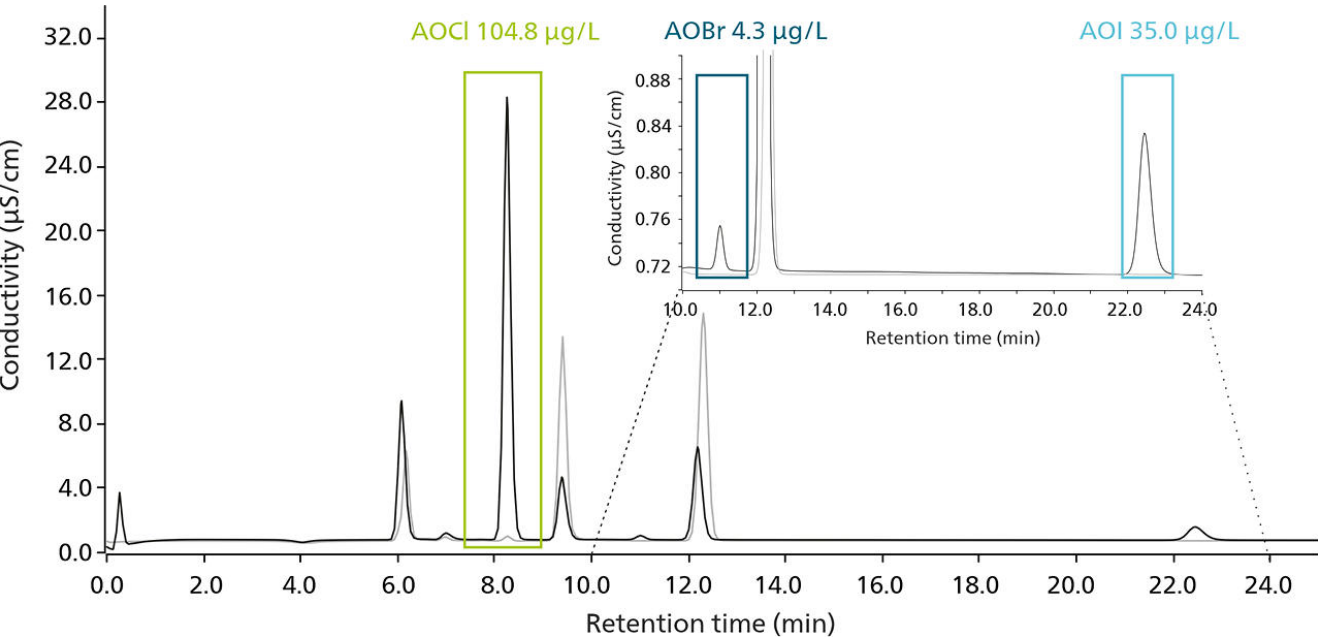


Figure 2. Chromatogram overlay of the blank and a wastewater sample for the determination of AOCl, AOB_r, and AOI measured from absorption column #1. To calculate the mass concentration of the individual AOX fractions, blank correction was performed according to Equation 1. No halogens were adsorbed on column #2, revealing the retention efficiency for AOX on column #1.

Overall, the validated procedure benefits from its easy, straightforward, and standardized handling, the precise determination of the analytes, the automatic calculation of results, and a low-maintenance, single-manufacturer setup.

A significant advantage of the standards DIN 38409-59 and ISO/DIS 18127 is that they allow the determination of adsorbable organically bound halogens as individual sum parameters, i.e., AOCl, AOBr, AOI, as well as AOF (a screening parameter for «total» PFASs, AN-CIC-033). Automation (e.g., automated eluent production, MiPT, intelligent and logical MagIC Net features) improves repeatability,

accuracy, and reliability of the results, saves valuable laboratory time for the liquid handling, standard, and eluent preparation, and allows 24/7 analysis – from which every laboratory, either research, routine, or governmental lab – can profit.

The world of organically bound halogens is so varied that these sum parameters enable insights about hot spots, transport pathways, but also particularly vulnerable regions in a very simple way, while complex targeted analysis, if at all, can resolve individual organically bound halogens for deeper investigations afterwards.

REFERENCES

1. Xu, R.; Xie, Y.; Tian, J.; et al. Adsorbable Organic Halogens in Contaminated Water Environment: A Review of Sources and Removal Technologies. *J Clean Prod* **2021**, 283.
2. Müller, G. Sense or No-Sense of the Sum Parameter for Water Soluble “Adsorbable Organic Halogens” (AOX) and “Absorbed Organic Halogens” (AOX-S18) for the Assessment of Organohalogenes in Sludges and Sediments. *Chemosphere* **2003**, 52 (2), 371–379.
3. Dann, A. B.; Hontela, A. Triclosan: Environmental Exposure, Toxicity and Mechanisms of Action. *J Appl Toxicol* **2011**, 31 (4), 285–311.
4. Xie, Y.; Chen, L.; Liu, R. AOX Contamination Status and Genotoxicity of AOX-Bearing Pharmaceutical Wastewater. *J Environ Sci* **2017**, 52, 170–177.

Internal reference: AW IC CH6-1438-042021

CONTACT

Metrohm AG
Ionenstrasse
9100 Herisau

info@metrohm.com

CONFIGURATION



Metrosep A Supp 5 - 250/4.0

Metrosep A Supp 5 - 250/4.0 , 150 mg/L 1 g/L
,Metrosep A Supp 5 - 250/4.0



Metrosep A Supp 5 Guard/4.0

Metrosep A Supp 5 Guard/4.0 IC Metrosep A Supp
5 7 ,
Metrosep A Supp 5 , PEEK(), («On Column Guard
System»), A Supp 5 Guard/4.0,



MiPT

Dosino




858 Professional Sample Processor – Pump

858 Professional Sample Processor – Pump 500 L
500 mL , 800 Dosino



930 Combustion IC PP (AJ)

930 Combustion IC PP (AJ) (), (Combustion IC),
Analytik Jena Combustion Module (2.136.0700) 920
Absorber Module 930 Compact IC Flex
Oven/SeS/PP/Deg MagIC Net, 930 Metrohm
Combustion IC-Paket Autosampler (Autosampler
MMS 5000), MagIC Net

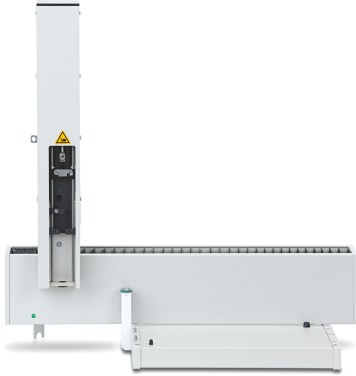

Autosampler MMS 5000 (AJ)

Analytik Jena() Autosampler MMS 5000 (AJ)
Metrohm Combustion IC Multi-Matrix, (6.7303.000)
(6.7302.000)