

# 利用拉曼光区分

This Application Note shows the ability of handheld Raman systems like the Mira M-1 to identify and differentiate between salts such as carbonates, phosphates, and sulfates. The main


focus of this work was to evaluate the influence of the cationic part and the water of crystallization on the Raman-spectroscopic identification of the mentioned salts.

## INTRODUCTION

Inorganic salts often have the same anionic part, while the cationic part differs. Many salts can be distinguished by the number of water molecules bond to the salt.

In this study, the influences on the Raman spectrum of the cationic part of the salt and the

number of water molecules bound to its anionic part were investigated. Although the differences between the salts are very small, the spectra that were recorded using a handheld Raman spectrometer differed sufficiently to differentiate between them.

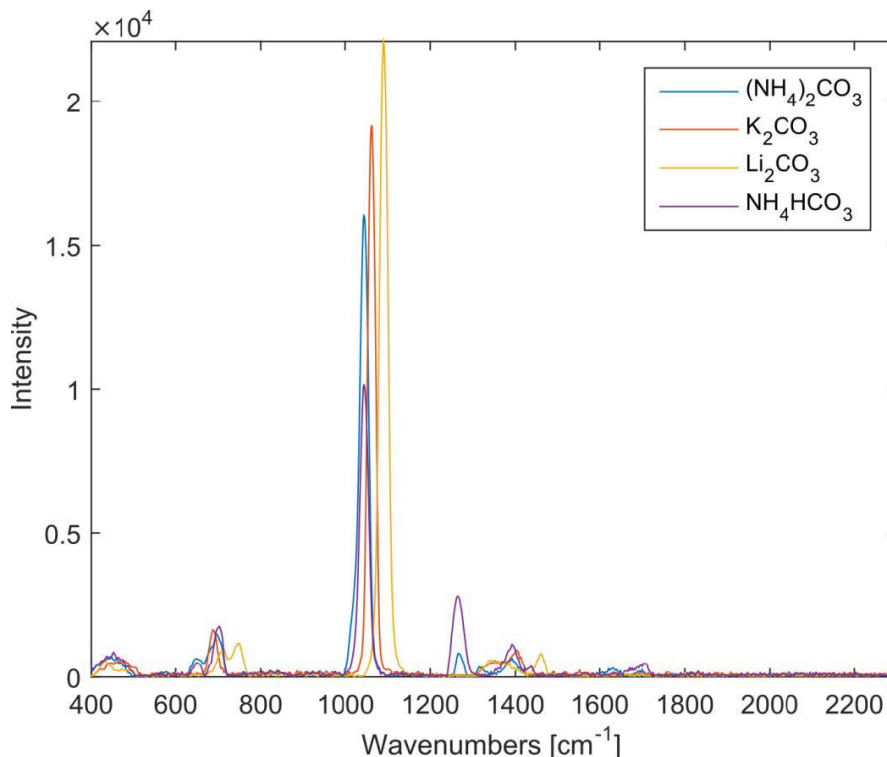

**Figure 1.** Overlaid Raman spectra of different carbonate salts.

## EXPERIMENTAL


All spectra were measured using the Mira M-1 Raman spectrometer in auto-acquisition mode, i.e., integration times were determined automatically. A laser wavelength of 785 nm and the Orbital-Raster-Scan (ORS) technique were used. The samples were analyzed in small

## RESULTS AND DISCUSSION

**Figure 2** shows how the carbonate peak in the Raman spectrum shifts depending on the counter cation. This shift is significant enough to use it to differentiate between the various carbonate salts (except for ammonium carbonate and ammonium bicarbonate, where the difference is very small).


Similar results are also observed for various phosphates; the differences in the spectra are useful for their differentiation (see **Figure 3**) – despite the fact that there are some difficulties

vials with the vial holder adapter. Three different anions that give are Raman-active were chosen: carbonate, phosphate, and sulfate. The cationic part of the salts was varied and the resulting spectral changes were investigated.

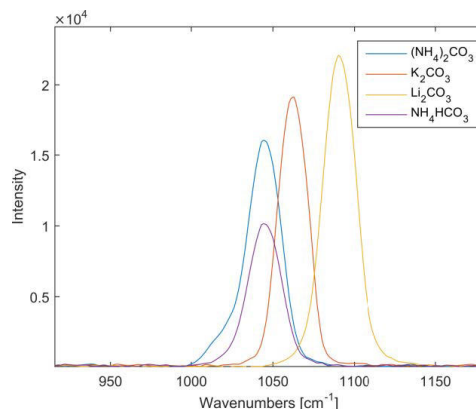

**Figure 2.** Main peaks of the various carbonates.


when trying to distinguish between diammonium phosphate and dipotassium phosphate (due to the similar ionic radius of potassium and ammonium).

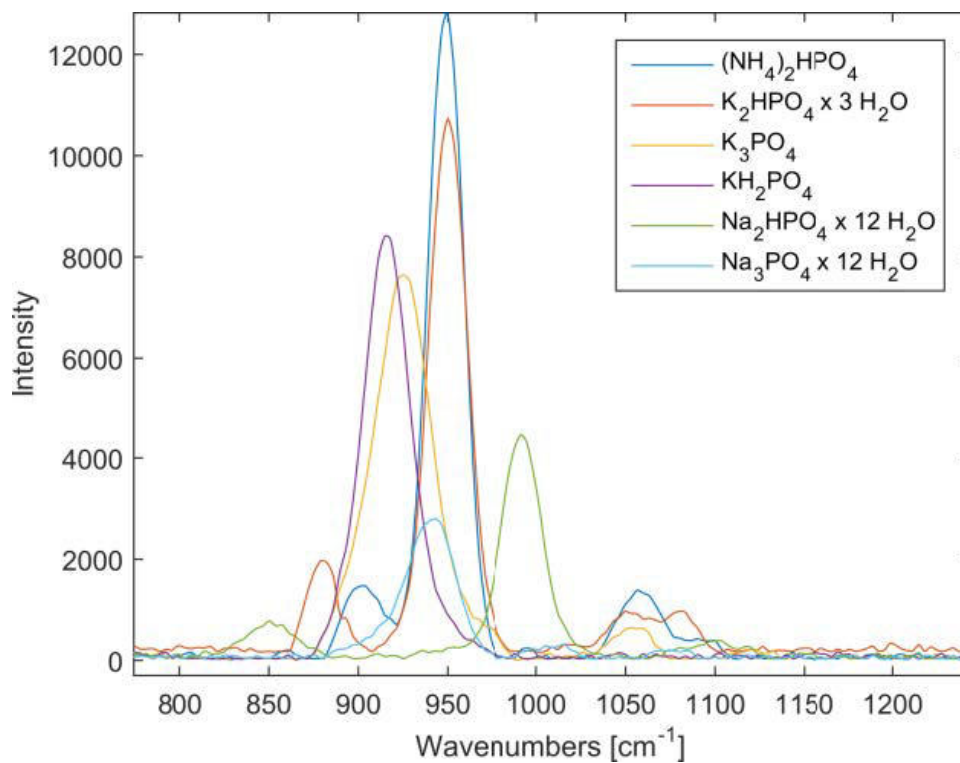

**Figure 3.** Main differences in the spectra of the phosphates.

The analyzed sulfates also differed from each other significantly, allowing their unambiguous identification with the handheld Mira M-1 (see **Figure 4**), although  $\text{CuSO}_4 \times 7 \text{H}_2\text{O}$ ,  $\text{K}_2\text{SO}_4$ ,

$\text{SnSO}_4$ , and  $\text{ZnSO}_4 \times 7 \text{H}_2\text{O}$  could not be identified unambiguously due to their similar spectra.

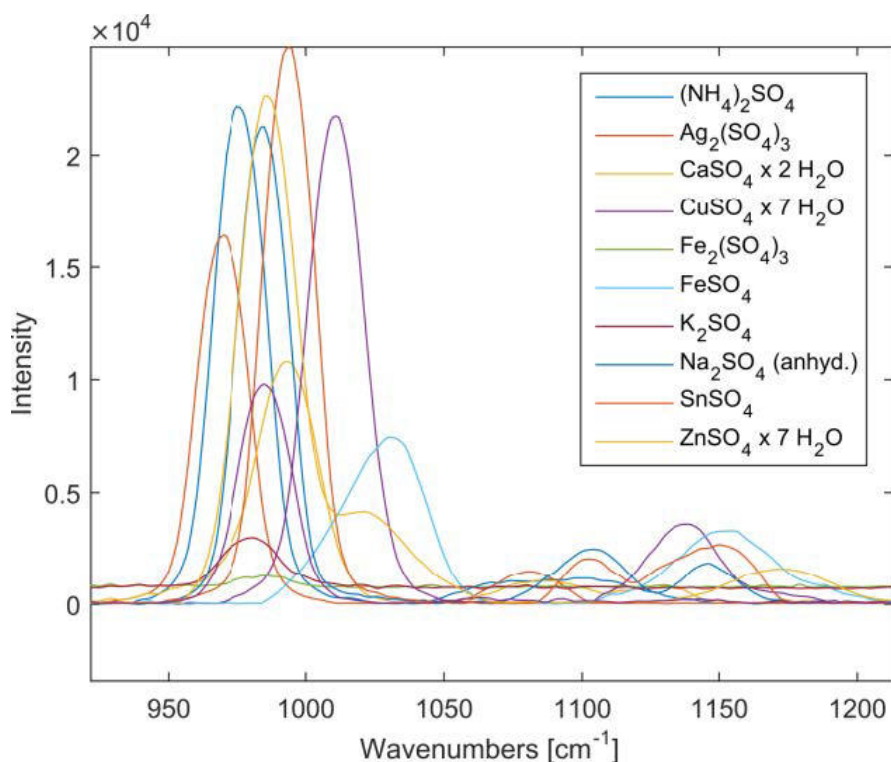

**Figure 4.** Main differences in the spectra of the sulfates.

## CONCLUSIONS

Generally speaking, handheld Raman spectrometers, such as the Mira M-1, are very useful when it comes to inorganic material identification, even when dealing with salts

having the same anionic part. Differences in the cationic parts of the salts help to unambiguously identify many salts with the Mira M-1 handheld Raman analyzer.

## CONTACT

Metrohm AG  
Ionenstrasse  
9100 Herisau

[info@metrohm.com](mailto:info@metrohm.com)

## CONFIGURATION



### MIRA P Advanced

瑞士万通快速拉曼分析 (MIRA) P 是一款性能大的手持式拉曼光,可用于各材料的快速无定和,例如物有效成分和形。MIRA P 小而固,配了瑞士万通的ORS技。MIRA P 符合 FDA 邦法 21 章第 11 款的定

。Advanced Package 包含一个附加透,可用它直接分析材料或者在材料容器中分析(3b 激光器),有一个小管支架套筒用于分析玻璃小管中的本(1 激光器)。