



## Application Note AN-RS-020

## 痕量粉中的金 O

用 Misa 保消費者安全

Auramine O (AO) is an industrial dye used for a broad range of manufactured products and as a fluorescent stain for detecting acid-fast bacteria in clinical specimens. Due to its intense yellow coloration, AO is also prized as an additive for enhancing the visual appeal of illicitly processed food products. Curry powder is a likely target for such adulteration, as it is a bright yellow mixture of several spices. Health hazards associated with

ingestion, and even improper handling of AO, include a high risk of several cancers, neural and liver toxicity, and even death. Despite bans on AO as a food additive, surveillance testing indicates its persistent use as an adulterant in foods and spices.

Misa (Metrohm Instant SERS Analyzer) achieves the rapid and sensitive detection of AO in curry powder in a simple assay format.


## INTRODUCTION

Misa is a versatile tool for the rapid and accurate detection of banned food colorants. This application

note details a facile extraction procedure for detecting AO in adulterated curry powder.

## REFERENCE MATERIAL AND LIBRARY CREATION

To establish a reference spectrum for AO, a pure standard in alkaline water (100 g/mL, pH 13) was analyzed using gold nanoparticles (Au NPs). The unique SERS spectrum shown in **Figure 1** can be used to create a library entry for AO.

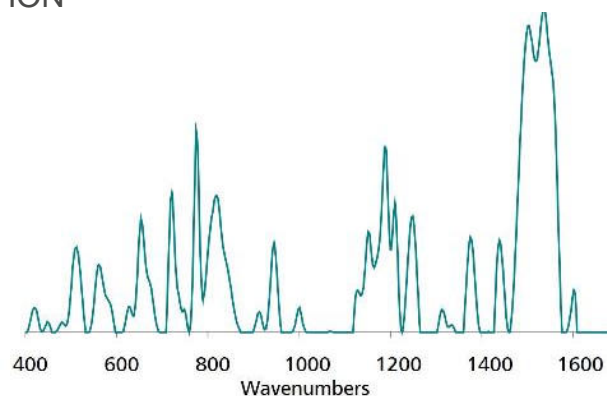

**Figure 1.** Standard SERS reference spectrum of Auramine O.

## EXPERIMENT

In a simulated test for AO in curry powder, solid AO was mixed with purchased curry powder to yield a concentration range of spiked test samples: 1000, 100, 10, 5, and 1 g/g. Liquid extraction of AO was performed by adding 1 mL of 0.1 mol/L NaOH to 100 mg of sample in a glass vial. This slurry was mixed and allowed to rest for 2 minutes. Ethyl acetate (EA, 1 mL) and NaCl (100 mg) were added to the vial, which was then inverted gently a few times (*do not shake vigorously*) to promote extraction of AO into the EA layer. After 10 minutes, 50 L of the top EA layer was added to a vial containing 400 L of Au NPs and 50 L of 0.5 mol/L NaCl. The vial was shaken to mix and immediately placed in the vial attachment on Misa for measurement.






**Table 1.** Experimental parameters

Instrument		Acquisition	
Firmware	0.9.33	Laser Power	5
Software	Misa Cal V1.0.15	Int. Time	10 s
Misa Vial Attachment	6.07505.040	Averages	10
ID Kit - Au NP	6.07506.440	Raster	ON

## RESULTS

Overlaid, baseline-corrected SERS spectra of basic EA extracts of curry powder spiked with varying concentrations of AO demonstrate reliable detection

down to 1 g/g (**Figure 2**). Note: Peaks in AO SERS spectra show solvent and pH-related shifts.

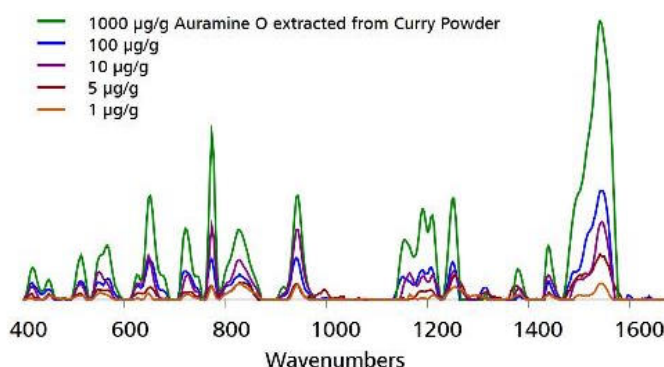**Figure 2.** Detection range of AO with Misa and Au NPs.

## FIELD TEST PROTOCOL

### Detection of Auramine O in the field

Using the large end of the scoop, add 3–4 scoops of sample to a 2 mL vial. Add NaOH solution to the vial until halfway full. Add 3–4 scoops solid NaCl, then fill vial to the top with ethyl acetate. Cap and invert the vial a few times to mix, but *do not shake the vial vigorously*. Let the sample rest for 5 minutes, as

distinct layers will form. Fill a *clean vial* halfway full with Au NPs. Using pipettes, add 2 drops each of the *top layer* of the sample solution and NaCl solution to Au NPs, cap and shake the vial gently to mix. Insert into vial attachment on Misa for measurement.

**Table 2.** Requirements for field test protocol

<b>ID Kit - Au NP</b>	<b>6.07506.440</b>
includes:	Gold nanoparticles (Au NP)
	Scoop
	Disposable pipettes
	2 mL glass vials
<b>Reagents</b>	
NaOH solution	0.4 g NaOH in 100 mL water
Solid NaCl	
Ethylene acetate	
NaCl solution	3 g NaCl in 100 mL water
<b>Test settings</b>	Use <b>ID Kit OP</b> on MISA

## CONCLUSION

The facile and sensitive detection of AO in adulterated curry powder is demonstrated using Misa. This analysis requires minimal user training and minimal consumables, making it an ideal analytical platform for on-site QC testing in food

manufacturing, shipping, and receiving facilities. Misa's portability and ease-of-use in trace detection of illicit colorants outperforms complex extraction and analysis procedures in a laboratory setting.

## CONTACT

Metrohm AG  
Ionenstrasse  
9100 Herisau

[info@metrohm.com](mailto:info@metrohm.com)

## CONFIGURATION



### MISA Advanced

Metrohm Instant SERS Analyzer (MISA) 是一款高性能、便携式分析系,可快速/定非法物、食品添加和微量食品染料。MISA 的特点是配了 Metrohm 的道光栅描 (ORS) 技的光。其空需求小和并且池寿命有所延,用于或移室用。MISA 提供各 1 激光附件,可活取。分析可通 BlueTooth 或 USB 接行。

MISA Advanced 套件是一个完整套件,其作用是用能用 Metrohms 米粒溶液和 P-SERS 条行 SERS 分析。MISA Advanced 套件包含了一个 MISA 小管附件、一个 P-SERS-附件、一个 ASTM 校正准件、一个 USB 迷、一个 USB 供元和用于行 MISA 器的 MISA Cal 件。随供了一个用来安全保管器和附件的固保箱。



### ID – Au NP

ID 套件 - Au NP 包含了 Mira/Misa 用使用体金溶液行 SERS 分析所需的件。套件包含了一个一次性抹刀、一个移液管、品小瓶和一个含金体的瓶子。