


OMNIS Titrator

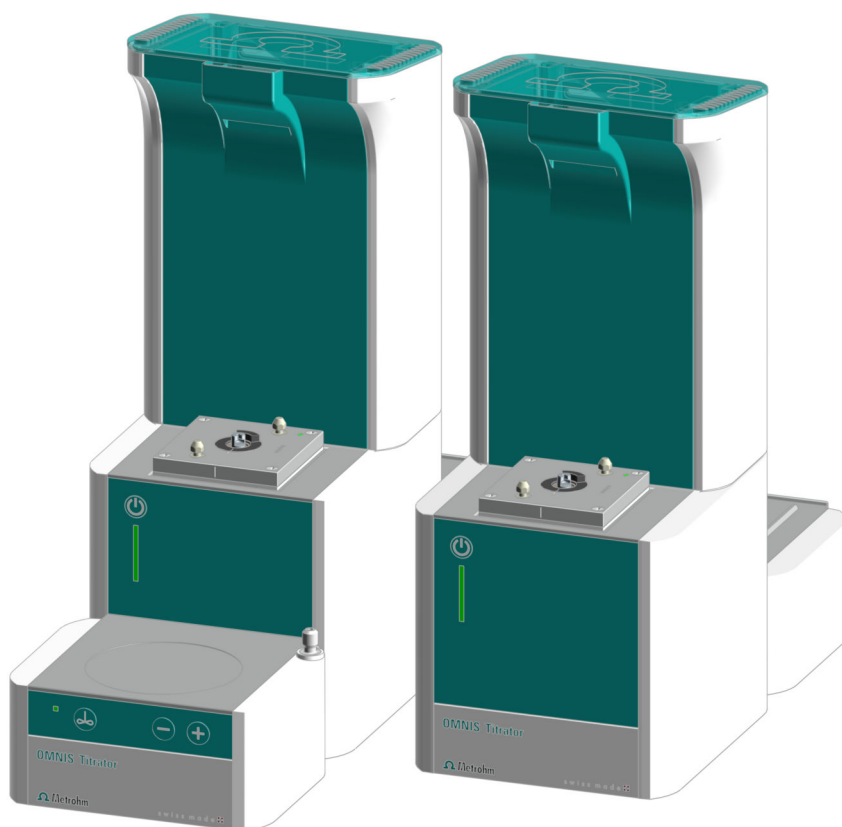

2.1001.0X20

Product manual

8.1001.8002EN / 2021-07-23



Metrohm AG
Ionenstrasse
CH-9100 Herisau
Switzerland
+41 71 353 85 85
info@metrohm.com
www.metrohm.com

OMNIS Titrator

2.1001.0X20

Product manual

8.1001.8002EN /
2021-07-23

This documentation is protected by copyright. All rights reserved.

This documentation has been prepared with great care. However, errors can never be entirely ruled out. Please send comments regarding possible errors to the address above.

Disclaimer

Deficiencies arising from circumstances that are not the responsibility of Metrohm, such as improper storage or improper use, etc., are expressly excluded from the warranty. Unauthorized modifications to the product (e.g. conversions or attachments) exclude any liability on the part of the manufacturer for resulting damage and its consequences. Instructions and notes in the Metrohm product documentation must be strictly followed. Otherwise, Metrohm's liability is excluded.

Table of contents

1	Overview	1
1.1	OMNIS Titrator – Product description	1
1.2	OMNIS Titrator – Product versions	1
1.3	Symbols and conventions	2
1.4	Further information	3
1.5	Accessories	3
2	Safety	4
2.1	Intended use	4
2.2	Responsibility of the operator	4
2.3	Requirements for operating personnel	5
2.4	Safety instructions	5
2.4.1	Danger from electrical potential	5
2.4.2	Danger from biological and chemical hazardous substances	5
2.4.3	Danger from highly flammable substances	6
2.4.4	Danger from leaking liquids	6
2.4.5	Danger during transport of the product	7
2.5	Design of warning messages	7
2.6	Meaning of warning signs	8
3	Functional description	10
3.1	OMNIS Titrator – Single-workplace titration system – Overview	10
3.2	OMNIS Titrator – Single-workplace KF titration system volumetric – Overview	12
3.3	OMNIS Titrator – Overview	13
3.3.1	Magnetic stirrer – Overview	15
3.3.2	Magnetic stirrer with accessories – Overview	16
3.3.3	Dosing unit – Overview	17
3.3.4	Bottle unit – Overview	21
3.3.5	Measuring Module Analog – Overview	24
3.3.6	Measuring Module Digital – Overview	25
3.3.7	Volumetric Karl Fischer titration cell – Overview	26
3.4	OMNIS Titrator – Function	27
3.4.1	Magnetic stirrer – Functional description	27
3.4.2	Dosing unit – Function	27
3.4.3	Bottle unit – Function	28
3.4.4	Measuring Module Analog – Functional description	29
3.4.5	Measuring Module Digital – Functional description	29

IV

8	Maintenance	77
8.1	Maintenance	77
8.2	Cleaning the product surface	77
8.3	Storing the cylinder unit OMNIS	78
8.4	Cleaning the cylinder unit OMNIS	80
8.5	Cleaning the OMNIS Liquid Adapter	82
9	Troubleshooting	85
9.1	Forcing a shutdown	85
10	Disposal	86
11	Technical specifications	87
11.1	Ambient conditions	87
11.2	OMNIS Titrator – Energy supply	87
11.3	Magnetic stirrer – Energy supply	87
11.4	Measuring module – Energy supply	88
11.5	OMNIS Titrator – Dimensions	88
11.6	Magnetic stirrer – Dimensions	88
11.7	Measuring module – Dimensions	89
11.8	OMNIS Titrator – Housing	89
11.9	Magnetic stirrer – Housing	89
11.10	Measuring module – Housing	90
11.11	OMNIS Titrator – Connectors specifications	90
11.12	Measuring Module Analog – Connectors specifications	91
11.13	Measuring Module Digital – Connector specifications	92
11.14	Display specifications	92
11.15	Measuring Module Analog – Display specifications	92
11.16	Measuring Module Digital – Display specifications	92
11.17	Measuring Module Analog – Measuring specifications ..	93
11.18	OMNIS Titrator – Liquid Handling specifications	94
11.19	Magnetic stirrer – Specifications	94

1 Overview

1.1 OMNIS Titrator – Product description

The OMNIS Titrator is the central instrument of an OMNIS titration system. The OMNIS Titrator is responsible for the energy supply to all of the modules in the titration system and for communication between the titration system and the OMNIS Software. The scope of the functions of the OMNIS Titrator is defined in the selected function license.

1.2 OMNIS Titrator – Product versions


The product is available in the following versions:

Table 1 Product versions

Article number	Designation	Magnetic stirrer	Function license
2.1001.0010	OMNIS Titrator	without	individual
2.1001.0020	OMNIS Titrator	with	individual
2.1001.0110	OMNIS Basic Titrator	without	Basic
2.1001.0120	OMNIS Basic Titrator	with	Basic
2.1001.0210	OMNIS Advanced Titrator	without	Advanced
2.1001.0220	OMNIS Advanced Titrator	with	Advanced
2.1001.0310	OMNIS Professional Titrator	without	Professional
2.1001.0320	OMNIS Professional Titrator	with	Professional

Information on function licenses is available on the [Metrohm website](#) or from the regional Metrohm representative.

The article number and serial number for identifying the product can be found on the type plate:

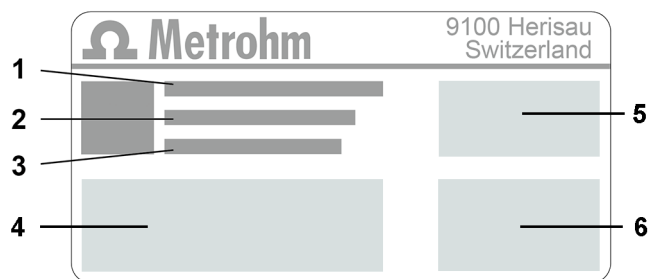

1	(01) = Article number in accordance with GS1 standard	2	(21) = Serial number
3	(240) = Metrohm article number	4	Certification
5	Certification	6	Technical specifications

1.3 Symbols and conventions

The following formatting may appear in the documentation:

(5-12)	Cross-reference to figure legend The first number refers to the figure number. The second number refers to the product part in the figure.
1	Instruction step Numbers indicate the order of the instructions steps.
Method	Names of parameters, menu items, tabs and dialog windows
File ► New	Menu path
[Continue]	Button or key

1.4 Further information

The Metrohm Knowledge Base <https://guide.metrohm.com> always provides the current version of this document. Further instructions, leaflets, release notes etc. may be available, depending on the product. You can directly access the required information or the associated PDF document using the full-text search function and filters.

1.5 Accessories

Up-to-date information on the scope of delivery and on optional accessories can be found on the Metrohm website. Download this information as follows:

Downloading the accessories list

- 1 Go to <https://www.metrohm.com>.
- 2 Enter the article number of the product (e.g. **2.1001.0010**) into the search field.
The search result is displayed.
- 3 Click on the product.
Detailed information regarding the product is shown on various tabs.
- 4 On the **Included parts** tab, click the link to download the PDF.
The PDF file with the accessories data is loaded.



NOTICE

Metrohm recommends downloading the accessories list from the Internet and keeping it for reference purposes.

2 Safety

2.1 Intended use

Metrohm products are used for the analysis and handling of chemicals.

Usage therefore requires the user to have basic knowledge and experience in handling chemicals. Knowledge regarding the application of fire prevention measures prescribed for laboratories is also mandatory.

Adherence to this technical documentation and compliance with the maintenance specifications make up an important part of intended use.

Any utilization in excess of or deviating from the intended use is regarded as misuse.

Specifications regarding the operating values and limit values of individual products are contained in the "Technical specifications" section, if relevant.

Exceeding and/or not observing the mentioned limit values during operation puts people and components at risk. The manufacturer assumes no liability for damage due to non-observance of these limit values.

The EU declaration of conformity loses its validity as soon as modifications are carried out on the products and/or the components.

2.2 Responsibility of the operator

The operator must ensure that basic regulations on occupational safety and accident prevention in chemical laboratories are observed. The operator has the following responsibilities:

- Instruct personnel in the safe handling of the product.
- Train personnel in the use of the product according to the user documentation (e.g. install, operate, clean, eliminate faults).
- Train staff on basic occupational safety and accident prevention regulations.
- Provide personal protective equipment (e.g. protective glasses, gloves).
- Provide suitable tools and equipment to carry out the work safely.

The product may be used only when it is in perfect condition. The following measures are required to ensure the safe operation of the product:

- Check the condition of the product before use.
- Remedy defects and malfunctions immediately.
- Maintain and clean the product regularly.

2.3 Requirements for operating personnel

Only qualified personnel may operate the product. Qualified personnel are persons who meet the following requirements:

- Basic regulations on occupational safety and accident prevention for chemical laboratories are known and complied with.
- Knowledge of handling hazardous chemicals is present. Personnel have the ability to recognize and avoid potential dangers.
- Knowledge regarding the application of fire prevention measures for laboratories is available.
- Safety-relevant information is communicated and understood. The personnel can operate the product safely.
- The user documentation has been read and understood. The personnel operate the product according to the instructions in the user documentation.

2.4 Safety instructions

2.4.1 Danger from electrical potential

Contact with electrical potential can cause serious injuries or death. To avoid danger from electrical potential, observe the following:

- Operate the product only if it is in perfect condition. The housing must also be intact.
- Only use the product with the covers fitted. If covers are damaged or missing, disconnect the product from the energy supply and contact the regional Metrohm service representative.
- Protect live components (e.g. power supply unit, power cord, connection sockets) against moisture.
- Always have maintenance work and repairs on electrical components carried out by a regional Metrohm service representative.
- Disconnect the product from the energy supply immediately if at least one of the following cases occurs:
 - The housing is damaged or open.
 - Live parts are damaged.
 - Moisture penetrates.

2.4.2 Danger from biological and chemical hazardous substances

Contact with biological hazardous substances may cause poisoning from toxins or infections from microorganisms. Contact with aggressive chemical substances may cause poisoning or chemical burns. To avoid danger from biological or chemical hazardous substances, observe the following:

- Label the product according to regulations if it is used for substances that have a potential for chemical hazards and are generally subject to the Hazardous Substances Ordinance.
- Wear personal protective equipment (e.g. protective glasses, gloves).
- Use exhaust equipment when working with vaporizing hazardous substances.
- Dispose of hazardous substances in accordance with regulations.
- Clean and disinfect contaminated surfaces.
- Only use detergents that do not cause any unwanted side reactions with the materials to be cleaned.
- Dispose of chemically contaminated materials (e.g. cleaning material) in accordance with regulations.
- Proceed as follows in case of a return shipment to Metrohm AG or a regional Metrohm representative:
 - Decontaminate the product or product component.
 - Remove the labeling for hazardous substances.
 - Create a declaration of decontamination and enclose it with the product.

2.4.3 Danger from highly flammable substances

Using highly flammable substances or gases may cause fires or explosions. To avoid danger from highly flammable substances, observe the following:

- Avoid ignition sources.
- Use protective grounding.
- Use exhaust equipment.

2.4.4 Danger from leaking liquids

Leaking liquids may cause injuries and may damage the product. To avoid danger from leaking liquids, observe the following:

- Check the product and its accessories for leakages and loose connections.
- Replace leaking parts and connecting elements without delay.
- Tighten loose connecting elements.
- Do not loosen tubing connections under pressure.
- Do not remove aspiration tubing under pressure.
- Carefully pull the ends of the tubing out of the containers.
- Carefully let liquids from tubing drain into suitable containers.
- Insert the buret tips completely into the containers.
- Remove and dispose of leaked liquids in accordance with regulations.
- If you suspect that liquid has penetrated the instrument, disconnect the instrument from the energy supply. Then have the instrument checked by a regional Metrohm service representative.

2.4.5 Danger during transport of the product

Chemical or biological substances may be spilled during the transport of the product. Parts of the product may fall down or may be damaged. There is a risk of injury from chemical or biological substances and pieces of broken glass. To ensure safe transport, observe the following:

- Remove loose parts (e.g. sample racks, sample vessels, bottles) before transport.
- Remove liquids.
- Lift and transport the product with both hands on the base plate.
- Lift and transport heavy products only according to instructions.

2.5 Design of warning messages

There are 4 hazard levels for warning messages. The following signal words are used for classifying the hazard levels in warning messages:

- **DANGER** indicates a hazardous situation which, if not avoided, will result in serious injury or death.
- **WARNING** indicates a hazardous situation which, if not avoided, could result in serious injury or death.
- **CAUTION** indicates a hazardous situation which, if not avoided, could result in minor or moderate injury.
- **NOTICE** indicates a hazardous situation which, if not avoided, could result in property damage.

Warning messages differ in design (color and warning sign) depending on the hazard level:



DANGER

Type and source of danger

Consequences when not observing the notice: An irreversible injury that may result in death is very probable.

- Measures to avoid the danger



WARNING

Type or source of danger

Consequences when not observing the notice: A serious injury that may result in death is probable.

- Measures to avoid the danger



CAUTION

Type or source of danger








Consequences when not observing the notice: A minor to moderate injury is probable.





- Measures to avoid the danger

2.6 Meaning of warning signs

This documentation uses the following warning signs:

Table 2 Warning sign according to ISO 7010

Warning sign	Meaning
	General warning sign
	Warning of electrical voltage
	Warning of hand injuries
	Warning of sharp object
	Warning of hot surface
	Warning of biological hazard
	Warning of toxic materials

Warning sign	Meaning
	Warning of flammable materials
	Warning of corrosive substances
	Warning of optical radiation
	Warning of laser beams

Depending on the intended use of the product, the corresponding warning sign stickers must be placed on the product.

3 Functional description

3.1 OMNIS Titrator – Single-workplace titration system – Overview

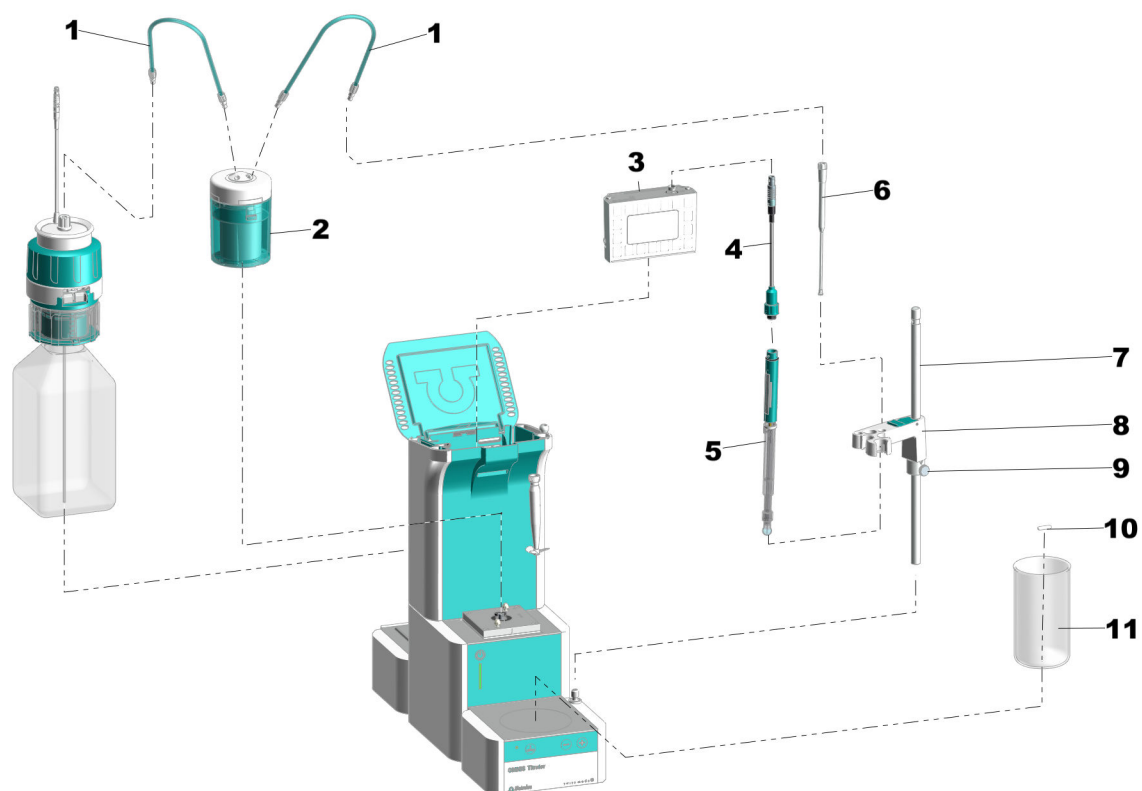


Figure 1 Single-workplace titration system – Front

1	FEP tubing	2	Cylinder unit
3	Measuring module	4	Electrode cable
5	Electrode	6	Antidiffusion tip
7	Support rod	8	Electrode holder
9	Clamping ring	10	Stirring bar
11	Sample beaker		



Figure 2 Single-workplace titration system – Rear

1	OMNIS Liquid Adapter	2	Bottle cap multi-use
3	Chemical bottle	4	Bottle holder

3.2 OMNIS Titrator – Single-workplace KF titration system volumetric – Overview

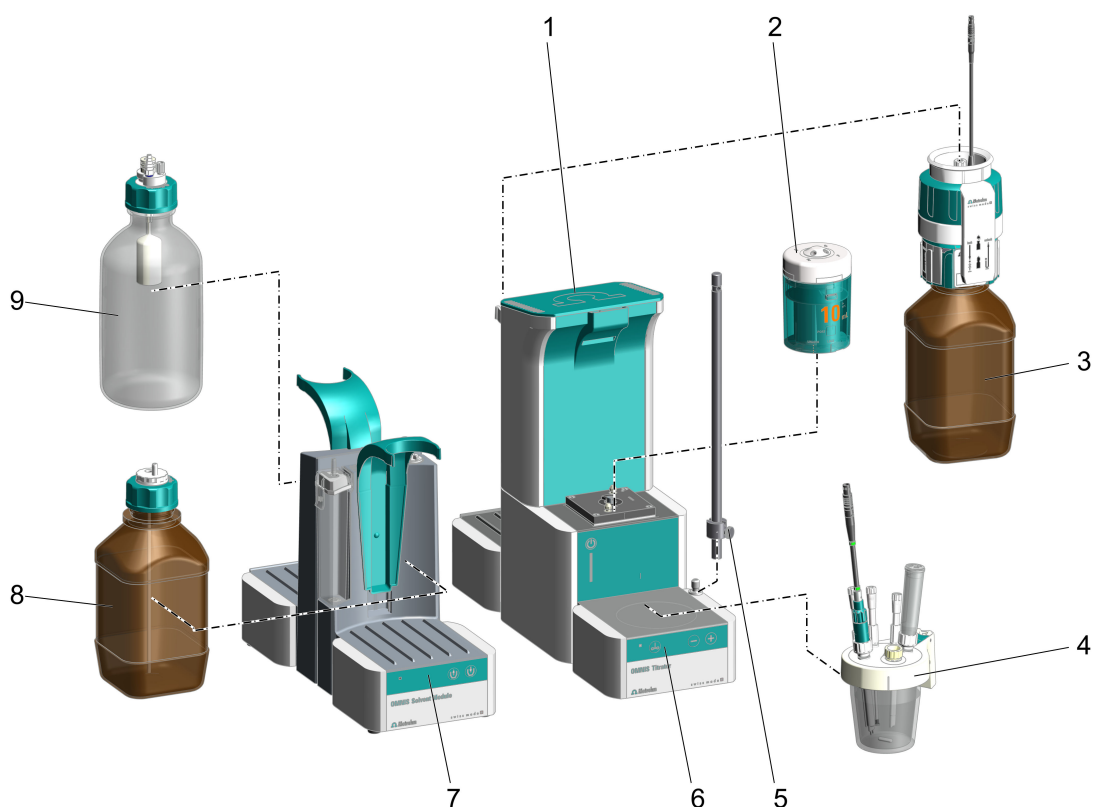


Figure 3 Single-workplace Karl Fischer titration system (volumetric) – Front

1	Measuring Module Analog, installed	2	Cylinder unit For the Karl Fischer titration, we recommend using the 10-mL cylinder units.
3	Chemical bottle (titrant) with OMNIS Liquid Adapter including bottle cap multi-use	4	OMNIS KF titration cell (volumetric)
5	Support rod with mounted clamping ring	6	OMNIS Titrator with magnetic stirrer
7	OMNIS Solvent Module	8	Chemical bottle (Solvent) with Siphon Breaker
9	Chemical bottle (Waste) with overflow protection		

3.3 OMNIS Titrator – Overview

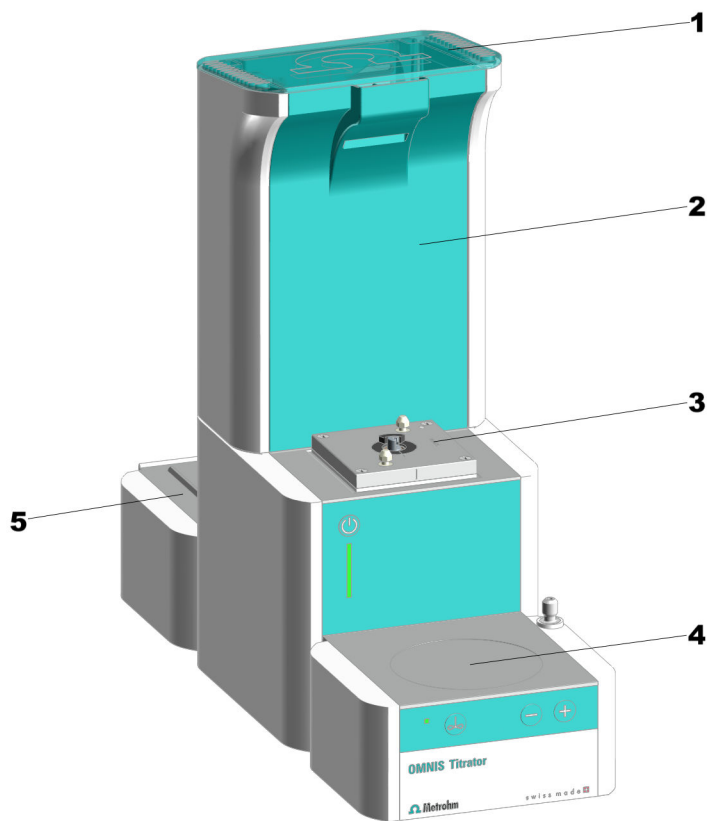


Figure 4 OMNIS Titrator – Front

1	Lid	2	Space for measuring modules
3	Dosing drive	4	Magnetic stirrer
5	Platform For chemical bottles		

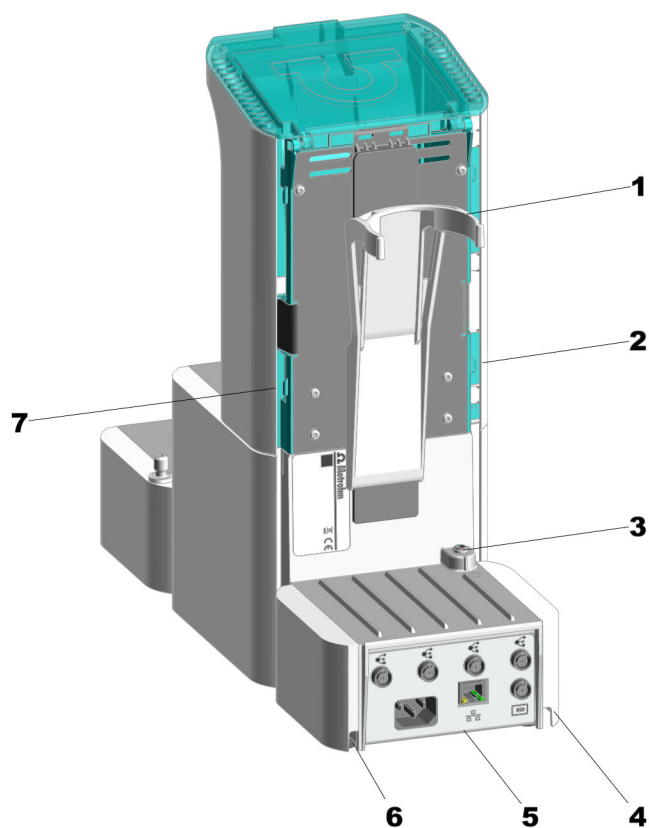


Figure 5 OMNIS Titrator – Rear

1 Bottle holder

2 Cable guide

On the side to the right. For Liquid Adapter cable

3 MSI connector

MSI = Metrohm Solution Identification. Connection socket for Liquid Adapter cable

4 Cable guide

Bottom right

5 Interfaces

6 Cable guide

Bottom left

7 Cable guide

On the side to the left. For rod stirrer cable or electrode cable

3.3.1 Magnetic stirrer – Overview

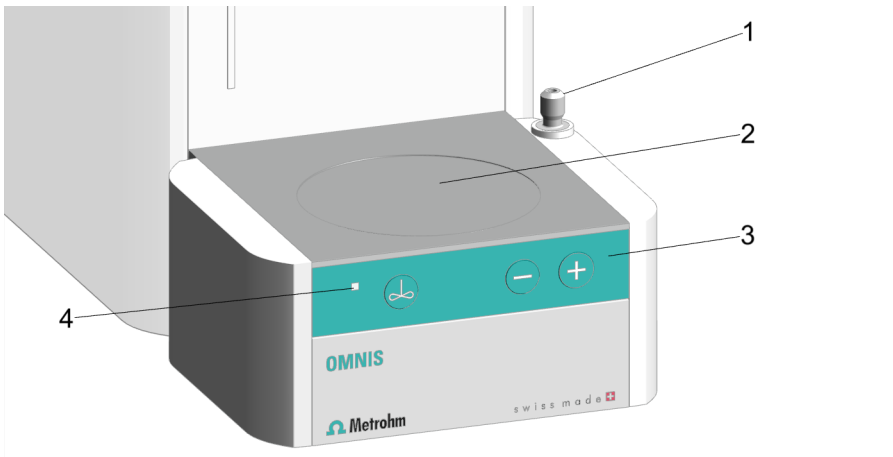


Figure 6 Magnetic stirrer – Overview

1	Stand attachment	2	Stirring area
3	Control bar	4	Status display LED. Multi-colored

3.3.2 Magnetic stirrer with accessories – Overview

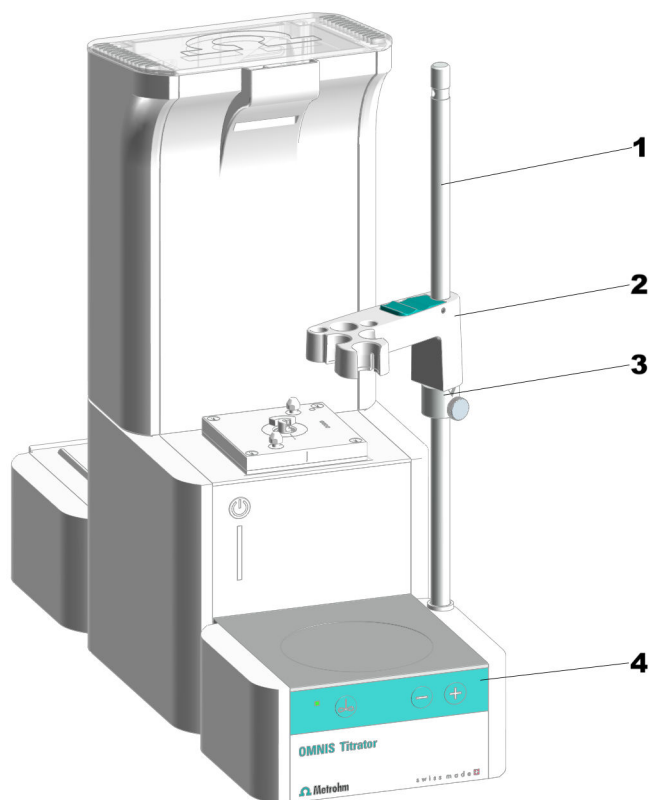


Figure 7 Magnetic stirrer with accessories – Overview

1	Support rod	2	Electrode holder
3	Clamping ring	4	Control bar



NOTICE

The support rod is grounded. The OMNIS Titrator and OMNIS Titration Module have a hole in the support rod that serves as earth contact for a banana plug (4 mm).

3.3.3 **Dosing unit – Overview**

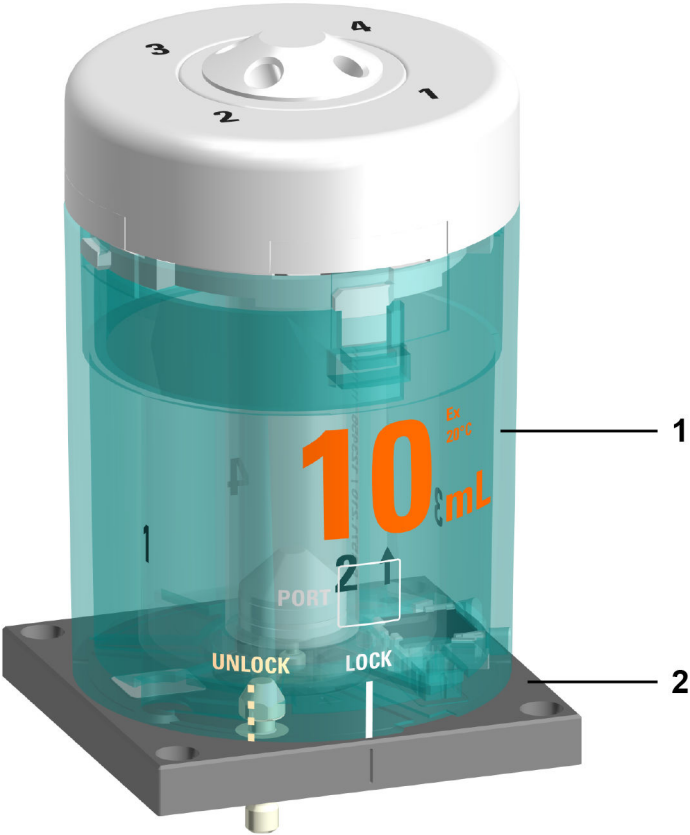


Figure 8 Dosing unit – Overview

- | | |
|-------------------------------------------------------------|------------------------------------------------------|
| 1 Cylinder unit
Available with various volumes | 2 Dosing drive
Not in scope of delivery |
|-------------------------------------------------------------|------------------------------------------------------|

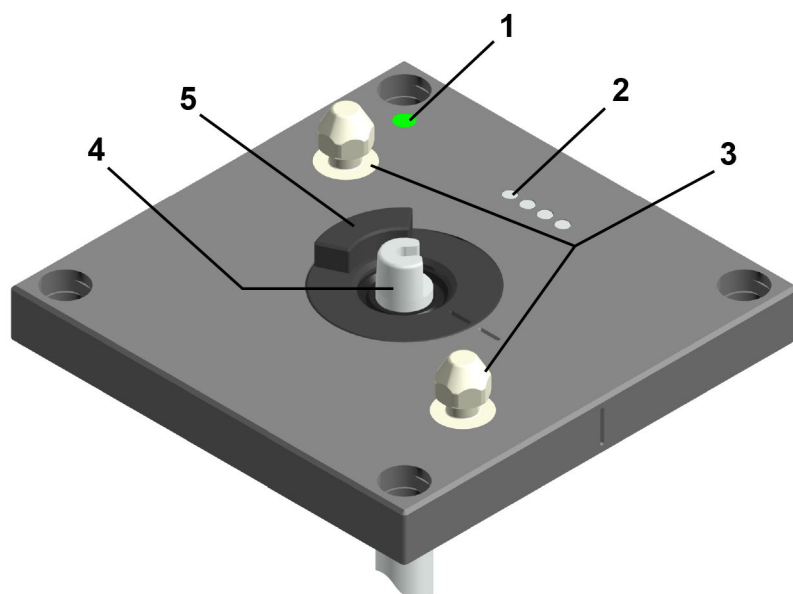


Figure 9 Dosing drive – Overview

- | | | | |
|----------|----------------------------------------------------|----------|-----------------------------------------------------------------|
| 1 | Status display
LED. Multi-colored | 2 | Contact pins
For communicating with the cylinder unit |
| 3 | Twistlocks
For locking the cylinder unit | 4 | Piston rod
For moving the dosing piston |
| 5 | Valve coupling | | |

3.3.3.1 **Cylinder unit OMNIS – Overview**

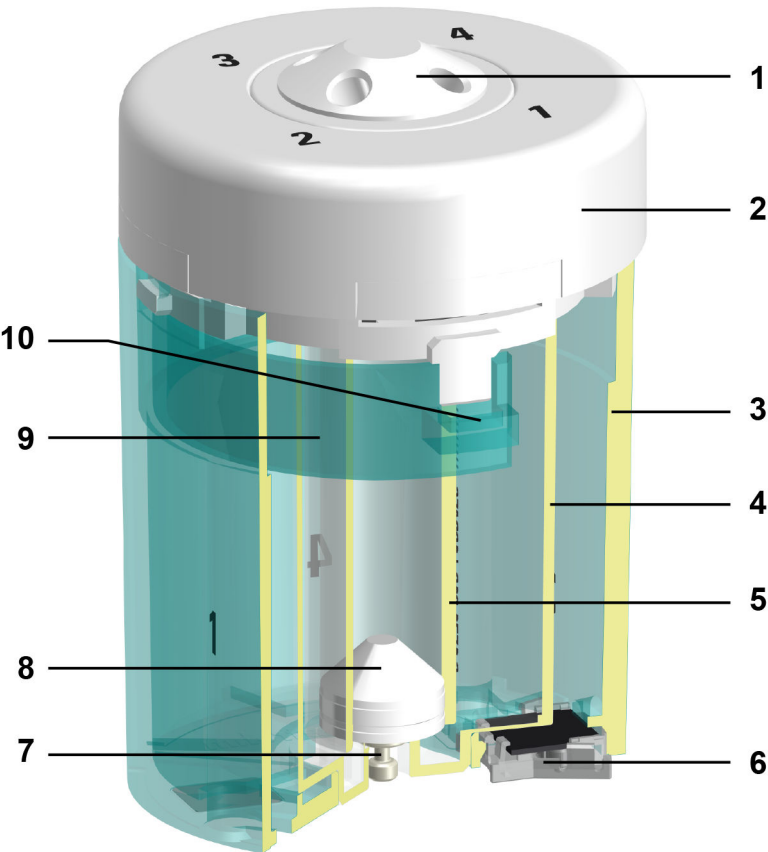


Figure 10 *Cylinder unit – Overview*

1	Distributor with 4 ports	2	Cylinder top piece
3	Cylinder housing	4	Centering tube
5	Cylinder	6	Data chip
7	Piston stopper	8	Dosing piston
9	Spring clip	10	Unlocking button

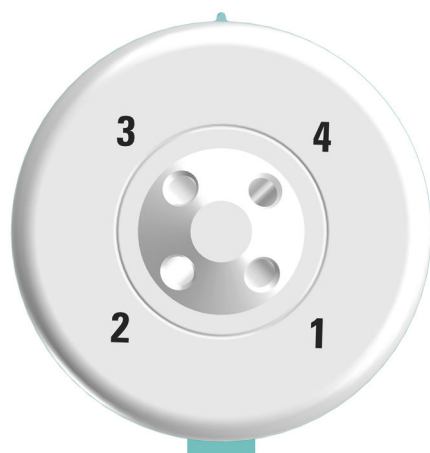


Figure 11 Cylinder unit – Overview from above

The following table shows the standard use of the 4 ports. The use of the ports can be changed in the OMNIS Software.

Port	Use	Connecting or sealing with
1	Dosing	Dosing tip
2	Filling the cylinder	Chemical bottle
3	Not used	Stopper
4	Not used	Stopper

3.3.4 Bottle unit – Overview

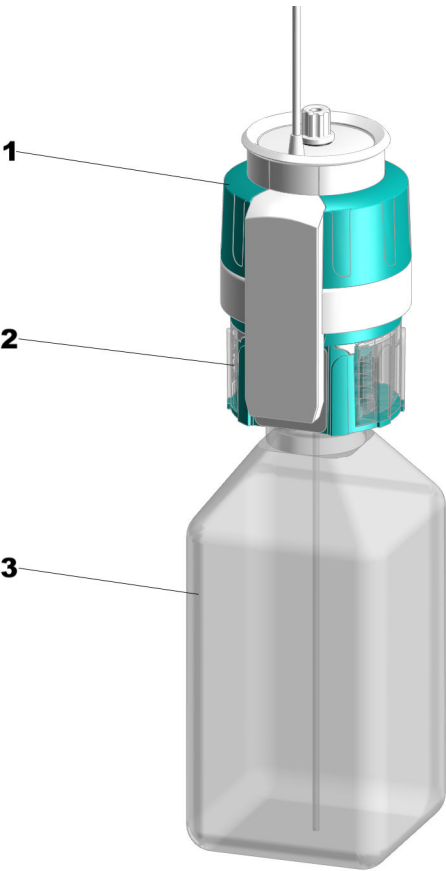


Figure 12 Bottle unit

1	OMNIS Liquid Adapter	2	Bottle cap multi-use
3	Chemical bottle		

3.3.4.1 OMNIS Liquid Adapter – Overview

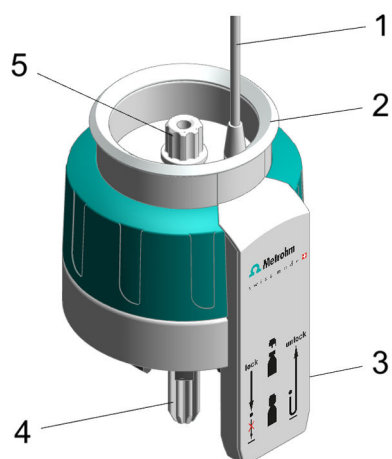


Figure 13 OMNIS Liquid Adapter – Parts

1	Cable	2	Status display
3	RFID reader	4	Aspiration tube Part of 6.01600.xxx
5	Tubing adapter Part of 6.01600.xxx		

3.3.4.2 **Bottle cap multi-use – Overview**

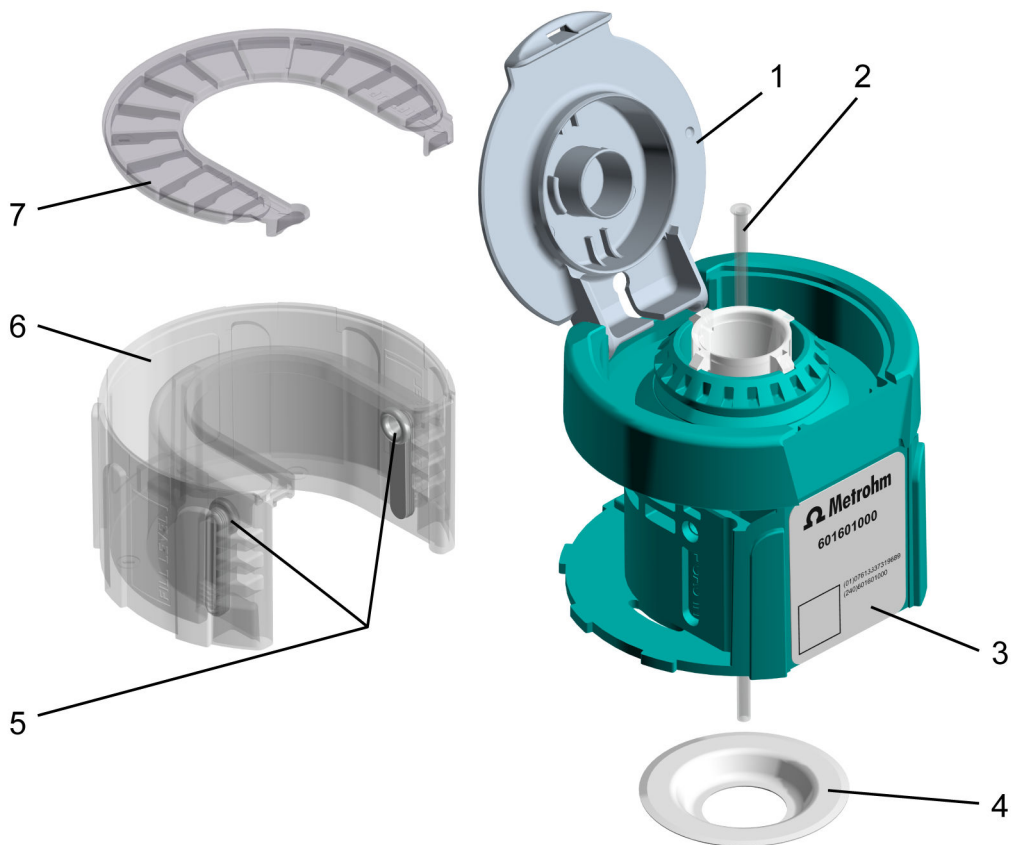


Figure 14 *Bottle cap multi-use, complete (6.01601.000)*

1 Flip-top lid	2 Aspiration tubing Aspiration tubing (6.1819.020)
3 RFID tag RFID chip for contact-free data transmission.	4 PTFE sealing ring PTFE sealing ring (6.02701.010)
5 Absorption cartridge insert (2 pieces)	6 Absorption cartridge housing
7 Absorption cartridge lid	5–7 Absorption cartridge, complete Absorption cartridge, complete (6.02701.000)

3.3.5 Measuring Module Analog – Overview

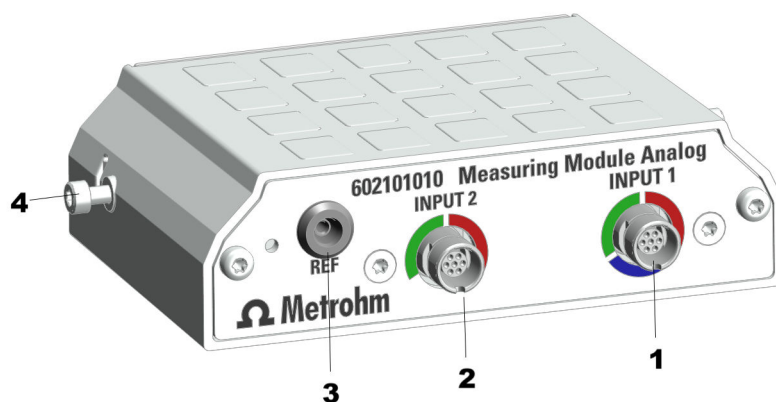


Figure 15 Measuring Module Analog – Overview

- | | |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------|
| <p>1 INPUT 1
Connection socket for potentiometric sensors (green coding), temperature sensors (red coding) and polarizable sensors (blue coding)</p> <hr/> | <p>2 INPUT 2
Connection socket for potentiometric sensors (green coding) and temperature sensors (red coding)</p> <hr/> |
| <p>3 REF
Connection socket for reference electrodes</p> | <p>4 Fastening screws
Fastening screws, left and right. These fasten the measuring module in the housing and ground the electronics.</p> |

3.3.6 Measuring Module Digital – Overview

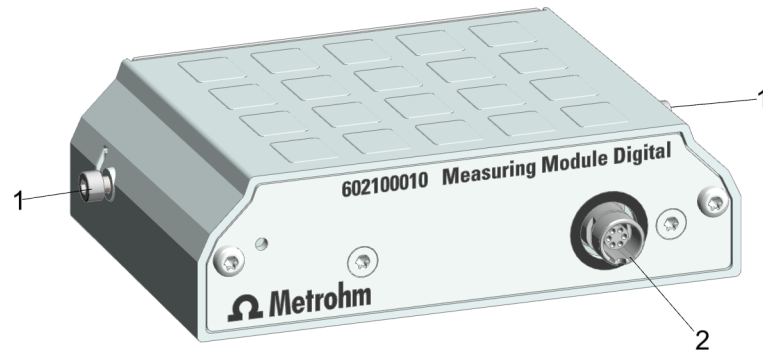


Figure 16 Measuring Module Digital – Overview

1 Fastening screws

Fastening screws, left and right. These fasten the measuring module in the housing and ground the electronics.

2 Connection socket

For dTodes

3.3.7 Volumetric Karl Fischer titration cell – Overview

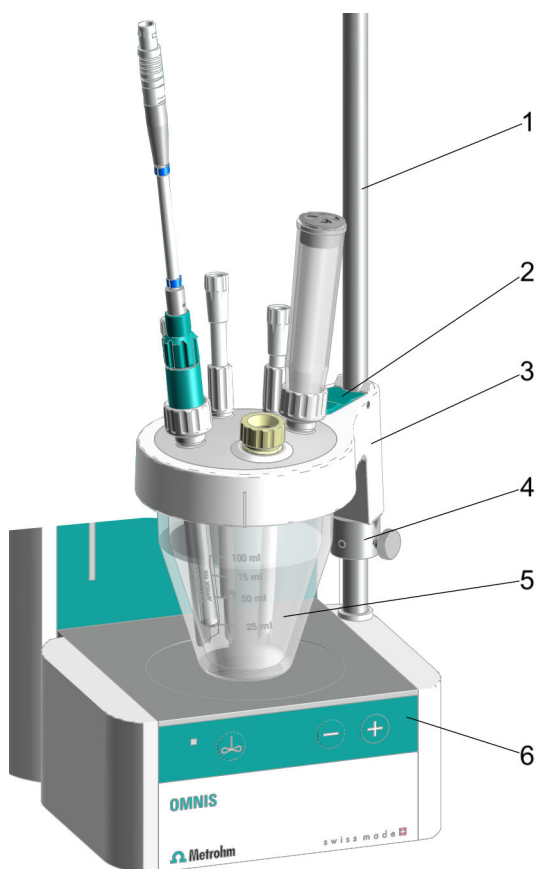


Figure 17 Installed volumetric Karl Fischer titration cell – Overview

1	Support rod	2	Locking lever
3	Titration vessel lid	4	Clamping ring
5	Titration vessel	6	Control bar for magnetic stirrer

3.4 OMNIS Titrator – Function

The OMNIS Titrator is connected to the energy supply and the Ethernet network. All of the other modules of the titration system are connected with the OMNIS Titrator.

The OMNIS Titrator is equipped with the following functional units:

- Connectors to the energy supply and the Ethernet network
- Interfaces for connecting additional modules
- Slots for a maximum of 2 measuring modules:
 - The Measuring Module Digital is used for connecting digital electrodes, the so-called dTodes.
 - The Measuring Module Analog is used for connecting analog electrodes.
- Dosing unit with an interchangeable cylinder unit
- 1 magnetic stirrer, depending on the product version
- Bottle unit with reagent detection
- Cable guides:

Cable guides are located at the rear and on the base of the titrator. The OMNIS Liquid Adapter cable can be guided to the connection socket in the cable guide on the right-hand side. The cable of a rod stirrer can be guided in the lateral and bottom cable guides on the left-hand side to the connection socket at the rear of the instrument.

The OMNIS Titrator with magnetic stirrer can also be used as a stand-alone instrument with manual sample change.

3.4.1 Magnetic stirrer – Functional description

The magnetic stirrer ensures that the sample is well mixed. The stirring rate can be adjusted depending on the amount and viscosity of the sample. The magnetic stirrer is operated via the instrument control bar or via the OMNIS Software.

3.4.2 Dosing unit – Function

Software control of the dosing unit can be used for accurate dosing of liquid volumes.

The dosing unit is comprised of the following units:

- Dosing drive
- Cylinder unit

The dosing drive is permanently installed in the housing of the instrument. The dosing drive is controlled via the OMNIS Software and is responsible for the accurate dosing of the solution.

If the cylinder unit is placed on top of the dosing drive, the dosing drive assumes responsibility for the following functions:

- **Raising and lowering the dosing piston:**

Solution is aspirated if the dosing piston is lowered. The cylinder fills up.

Solution is dosed if the dosing piston is raised. The cylinder empties.

- **Rotating the cylinder:**

The rotation of the cylinder controls which of the 4 ports the solution flows through.

The valve disk with an opening is located in the middle of the cylinder base.

The distributor disk with 4 openings corresponding to the 4 ports of the distributor is located at the bottom in the cylinder top piece.

The dosing drive rotates the cylinder by 90° stages so that the opening of the valve disk fits with an opening on the distributor disk. This results in a passage for the solution to the corresponding port of the distributor.

3.4.2.1 Cylinder unit OMNIS – Function

The cylinder unit is an accessory set for the dosing unit in OMNIS instruments with dosing function. The cylinder unit provides the liquid volume required for the analysis and can be equipped with different volumes. A distributor with 4 ports enables the filling and emptying of the cylinder.

The cylinder units 2 mL to 20 mL are used for dosing as well as titrating. The cylinder unit 50 mL is particularly suited for dosing (Liquid Transfer).

See also

Cylinder unit OMNIS – Overview (chapter 3.3.3.1, page 19)

3.4.3 Bottle unit – Function

The bottle unit contains the chemicals that are required for the analysis. The bottle unit in an OMNIS system is comprised of the following elements:

- Chemical bottle

- OMNIS bottle cap

The OMNIS bottle cap fits on a chemical bottle with GL 45 thread. Suitable adapters are available for chemical bottles with other threads.

The green bottle cap multi-use is equipped with a writeable RFID tag.

The RFID tag can be inscribed with information regarding the solution.

The red bottle cap single-use is equipped with an RFID tag. The RFID tag contains manufacturer's information regarding the composition and concentration of the solution.

- OMNIS Liquid Adapter

The OMNIS Liquid Adapter has an RFID tag reader that transfers the data between the RFID tag of the bottle cap and the OMNIS Software.

3.4.3.1 OMNIS Liquid Adapter – Functions

The basic principle

There is no fixed connection between the OMNIS Liquid Adapter and the chemical bottle that must be loosened when changing bottles. This means that you can transfer the OMNIS Liquid Adapter without difficulty from one bottle to another.

Liquid transfer

The liquid transfer is the responsibility of an aspiration tube installed in the OMNIS Liquid Adapter. When you place the OMNIS Liquid Adapter on the bottle, the aspiration tube presses on the aspiration tubing in the chemical bottle. A tightly sealed connection is made.

Contact-free data transmission

Information regarding the content of the chemical bottle is stored on an RFID tag on the bottle cap. The OMNIS Liquid Adapter has an RFID tag reader that reads this data.

The data cable fastened to the OMNIS Liquid Adapter transfers the information to the analyzer and to the software.

3.4.4 Measuring Module Analog – Functional description

The Measuring Module Analog contains the electronics necessary for the use of analog sensors. Analog electrodes and analog reference electrodes can be plugged into its connection sockets.

3.4.5 Measuring Module Digital – Functional description

The Measuring Module Digital is used as an interface for digital electrodes on an OMNIS Titrator or an OMNIS Titration Module.

The digital electrodes, the dTodes, can be plugged into its connection socket.

3.4.6 Volumetric Karl Fischer titration cell – Function

The volumetric Karl Fischer titration cell (KF titration cell) is a closed vessel for water content determination according to Karl Fischer that is fastened to the support rod of the magnetic stirrer.

The KF titration cell consists of:

- A titration vessel (different variants are available)
- A titration vessel lid (different variants are available, with and without homogenizer (Polytron PT 1300 D))

Seepage of moisture into the KF titration cell is prevented with seals and with the adsorber tube, which is filled with molecular sieve.

3.5 OMNIS main instrument – Indicators and controls

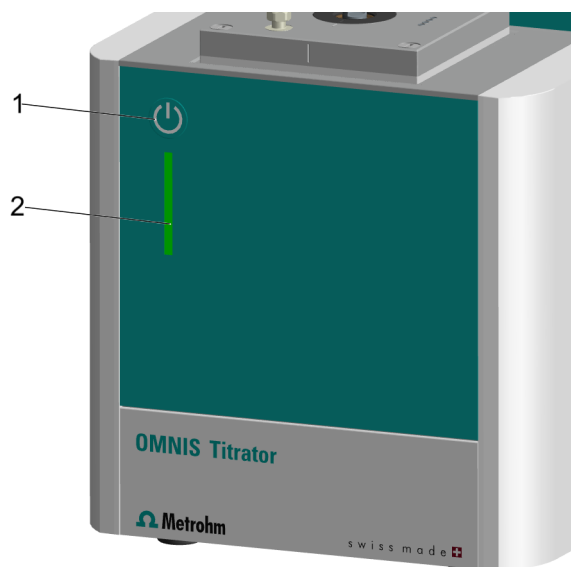


Figure 18 OMNIS main instrument – Indicators and controls

- | | | | |
|----------|----------------------|----------|----------------------------------------|
| 1 | On/off switch | 2 | Status display
Multi-colored |
|----------|----------------------|----------|----------------------------------------|

Indicators

The status of the main instrument is displayed with the status display (18-2) using different colors *System – Signals* (see chapter 3.6, page 31).

Controls







The on/off switch (**18-1**) is used for the hardware-side operation of the OMNIS main instrument.

Table 3 Behavior of the on/off switch

Pressure duration	Acoustic signal	Function on the OMNIS main instrument
Short pressing (1 to 5 s)	Beep after 1 s	Switch on the instrument. Shut down the instrument.
Very long pressing (> 10 s)	Continuous beep after 8 s	Force a shutdown.

3.6 System – Signals

System components with status indicators show their operating status with colors and/or flashing patterns. The meaning of the colors and flashing patterns is explained in the following table.

Visual signal		Meaning
	LED lights up yellow.	System start or initialization
	LED flashes yellow (slowly).	Ready for connection setup or locking
	LED flashes yellow (fast).	Connection setup started or locking underway
	LED lights up green.	Ready for operation
	LED flashes green (slowly).	In operation
	LED flashes red (fast).	Malfunction or error

Some system components only use part of the explained flashing patterns.

3.7 OMNIS Titrator – Interfaces

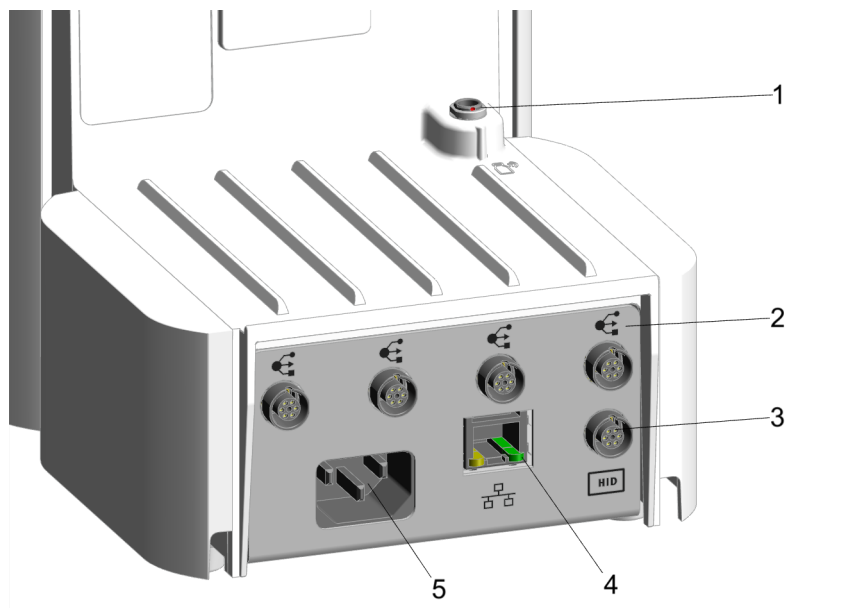


Figure 19 OMNIS Titrator – Interfaces and connectors

- | | |
|---------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------|
| <p>1 MSI connector
MSI = Metrohm Solution Identification. Connection socket for Liquid Adapter cable</p> <hr/> | <p>2 MDL connectors
MDL = Metrohm Device Link. Connection socket for connecting cable between OMNIS products</p> <hr/> |
| <p>3 HID connector
HID = Human Interactive Device. Connection socket for external operating units</p> <hr/> | <p>4 LAN connector
LAN = Local Area Network. Connection socket for a connecting cable to the local network</p> <hr/> |
| <p>5 Power socket
Connection socket for the energy supply</p> | |

3.8 Measuring Module Analog – Interfaces

The Measuring Module Analog has 3 measuring inputs for analog electrodes.

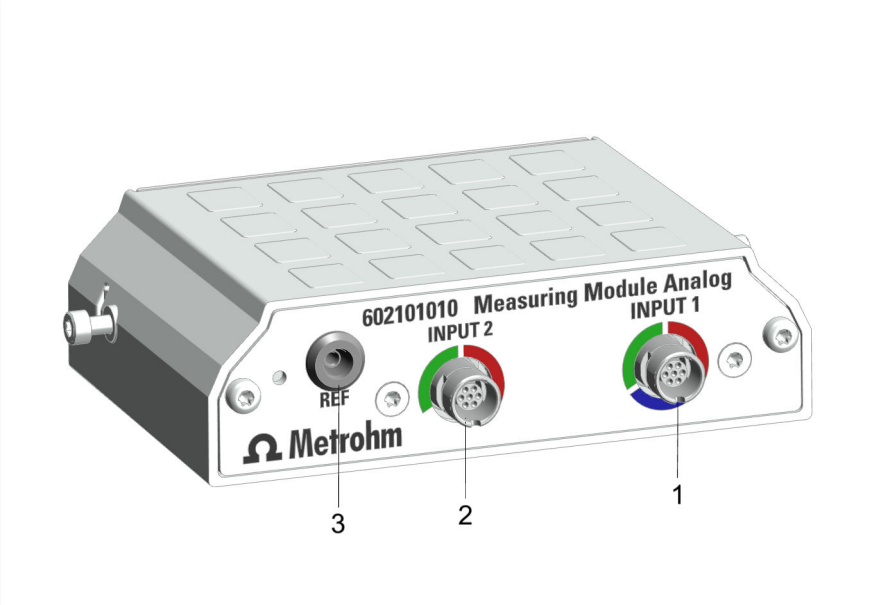


Figure 20 Measuring Module Analog – Measuring inputs

1	INPUT 1	2	INPUT 2
3	REF		

Measuring inputs INPUT 1 and INPUT 2

The measuring inputs **INPUT 1** and **INPUT 2** are marked with colored circle segments. The markings indicate that only certain types of electrode cables may be plugged into the connection socket:

Table 4 Meaning of the colors

Red	The connector supports temperature sensors.
Blue	The connector supports polarized sensors.
Green	The connector supports potentiometric sensors.

REF measuring input

Reference electrodes can be plugged into the **REF** measuring input.

3.9 Measuring Module Digital – Interfaces

The Measuring Module Digital has 1 connection socket for a digital electrode.





Figure 21 Measuring Module Digital – Interfaces and connectors

1 Connection socket for digital electrodes

Connection socket

The black marking around the connection socket indicates that only a cable of a digital electrode may be plugged in here.

1

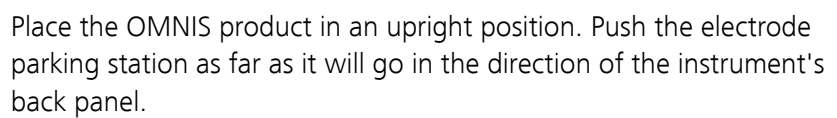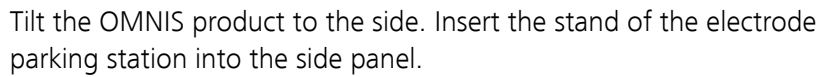


■■■■■■■ 37

- Titration tip (6.1543.200)

1 Aligning the cylinder unit


Rotate the cylinder unit until the marking with the label **UNLOCK** is in line with the marking on the dosing drive.

2

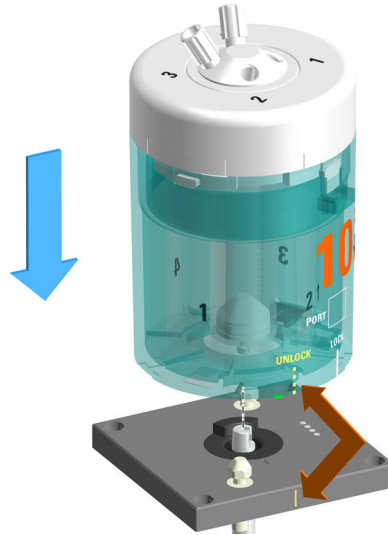

Set the cylinder unit down onto the two twistlocks straight from above.

3 Locking the cylinder unit

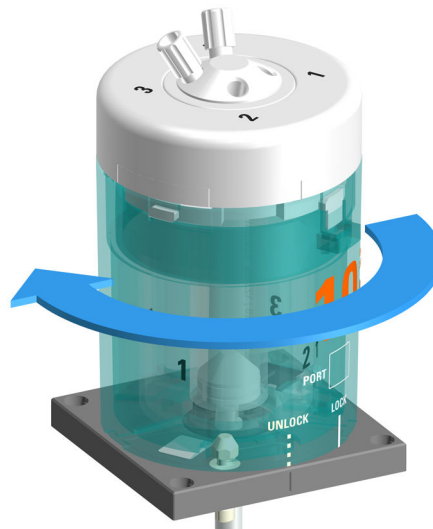

Rotate the cylinder unit to the left until it stops.
Use the marking with the label **LOCK** as a guide.

4 Mounting the tubings



Tighten an FEP tubing (6.1805.100) into port 1.

This FEP tubing is used as dosing tubing. Tighten the other end to the titration tip (6.1543.200).

5 Tighten the other FEP tubing (6.1805.100) into port 2.

This FEP tubing is used as filling tubing. Tighten the other end to the OMNIS Liquid Adapter.

6 Firmly tighten the tubing with the wrench (6.2739.000).

See also


Cylinder unit OMNIS – Overview (chapter 3.3.3.1, page 19)

5.5 Magnetic stirrer – Installing the accessories


Mounting the electrode holder

Accessories



- Support rod, 30 cm (6.2016.050)
- Clamping ring, 10 mm (6.2013.010)
- Electrode holder (6.02005.000)

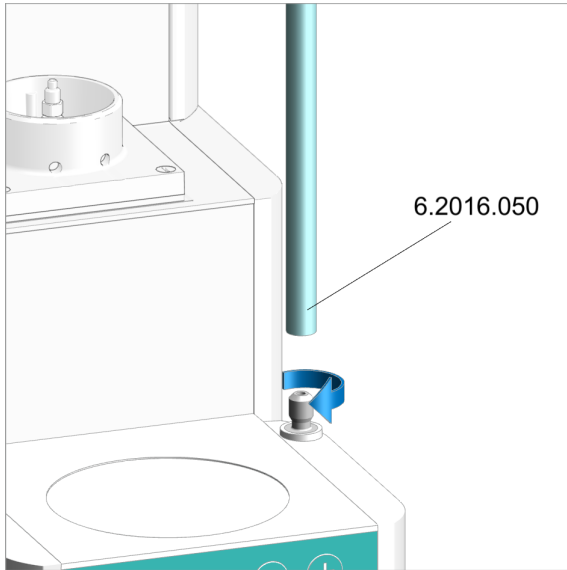**Mounting the support rod**

Screw the support rod onto the stand attachment .

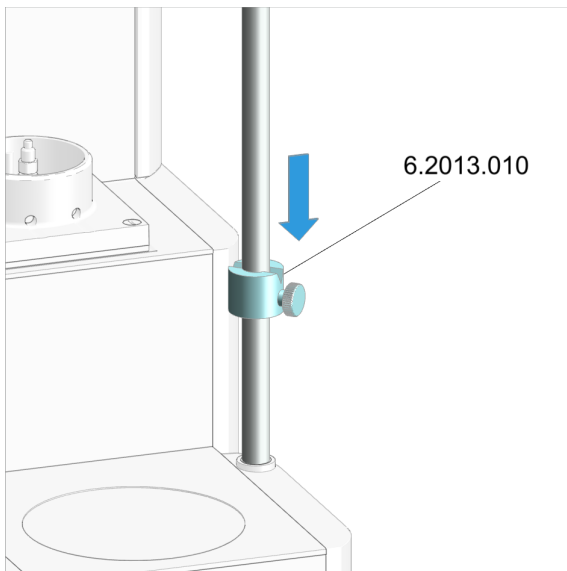**Mounting the clamping ring**

Push the clamping ring over the support rod with the indent facing upward.

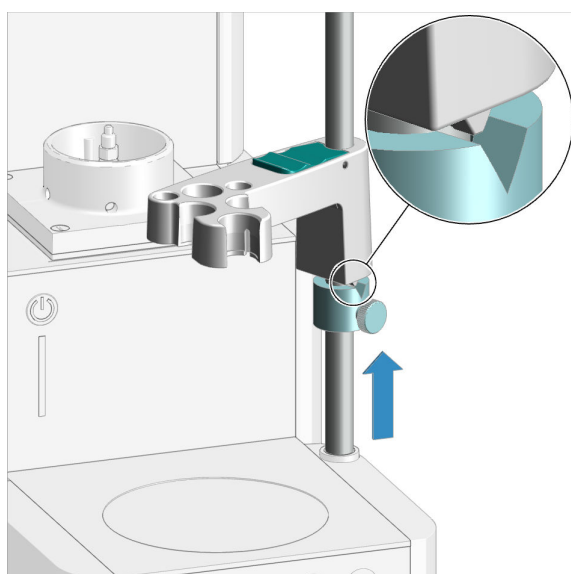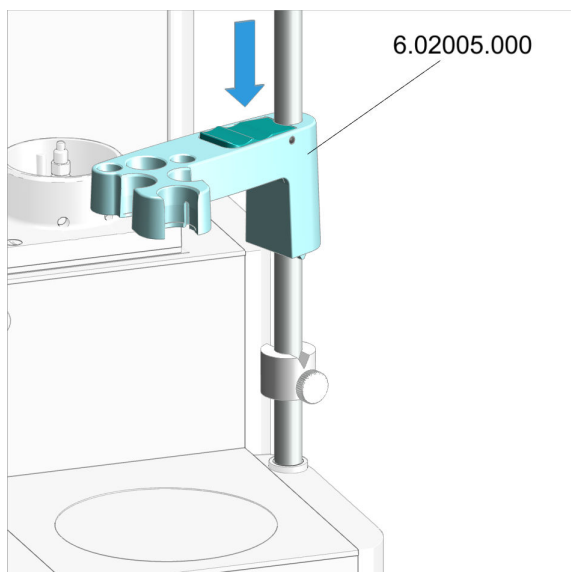

Mounting the electrode holder

1. Press the green locking lever on the electrode holder.
2. Push the electrode holder over the support rod.
3. To fix in place, release the green locking lever at the desired height.

The electrode holder is fixed in place.




NOTICE

The clamping ring is used as the lower stop for the electrode holder. The clamping ring prevents the electrode holder with the mounted electrode from being lowered too far.

1. Push the clamping ring under the electrode holder.
2. Rotate the clamping ring in such a way that the wedge on the electrode holder fits in the indent in the clamping ring.
3. Fix the clamping ring with the knurled screw in place at the desired height.

Preparing the magnetic stirrer for operation

Placing beaker and stirring bar

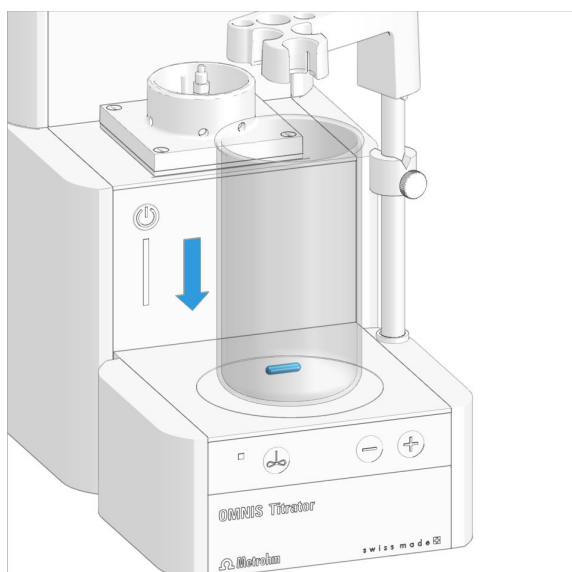

1. Place a 16 mm (6.1903.020) or 25 mm (6.1903.030) PTFE stirring bar in the sample beaker.
2. Place the sample beaker on the contact surface of the magnetic stirrer.

5.6 OMNIS Karl Fischer products – Replacing the adsorber material




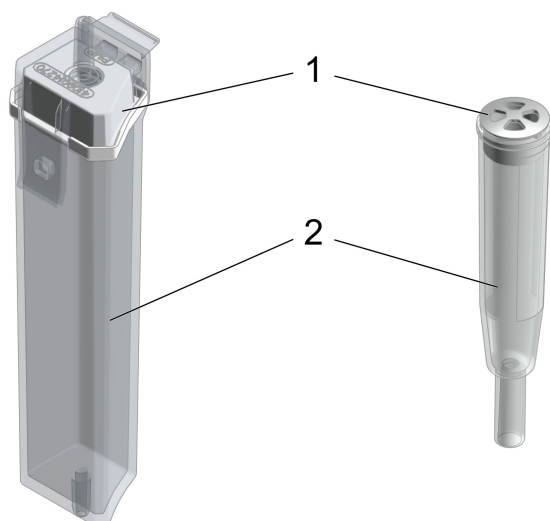


Figure 22 Adsorption cartridge of the OMNIS Solvent Module and adsorber tube for the KF titration cell

1 Lid

2 Housing

Replacing the adsorber material of the adsorption cartridge or adsorber tube

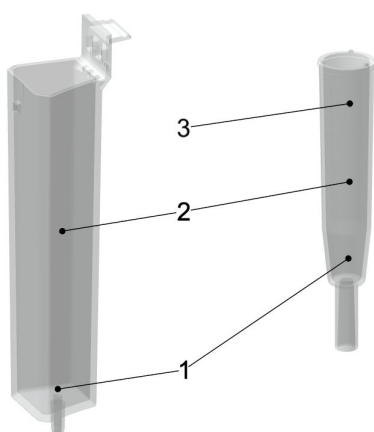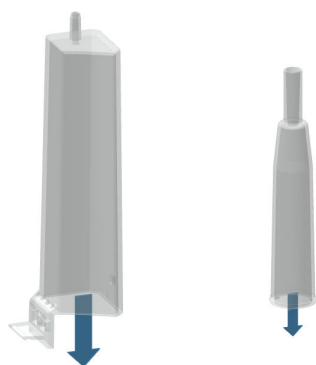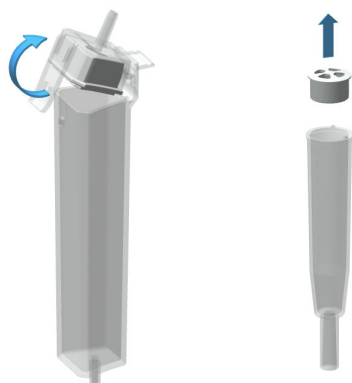

1. Removing the lid

Adsorption cartridge: Unlatch and remove the lid including the seal from the housing.

Adsorber tube: Remove the lid by pulling it out of the housing.

2. Removing the adsorber material

Remove the entire content.

This step is not necessary if the housing is empty.




NOTICE

The molecular sieve can be regenerated at 300 °C in the drying oven, see <https://www.metrohm.com/en/support-und-service/faq-kft/>.

3. Filling the housing with adsorber material

1. Place a cotton plug loosely at the base of the housing. Do not pack the cotton too tightly as sufficient gas flow has to be possible.
2. Fill the housing with molecular sieve to approx. 1 cm under the housing edge.
3. **Adsorber tube:** Place a small cotton plug on the molecular sieve. Do not pack the cotton too tightly as sufficient gas flow has to be possible.

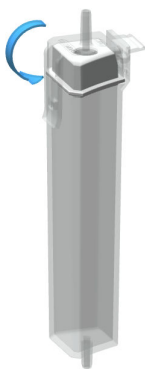

4. Sealing the housing with the lid



NOTICE

Make sure that the sealing surface between the housing and the lid is clean and dry and that there are no residuals of filling material!

Adsorption cartridge: Hook the lid including the seal into the housing side and close it by clicking it into place.

Adsorber tube: Seal the housing with the lid.



NOTICE

At average humidity, we recommend replacing the adsorber material approx. every 6 weeks.

An increase in drift indicates that the leak-tightness of the KF titration cell should be inspected and that the molecular sieve should possibly be replaced.

Hint:

Write the date on the adsorption cartridge when you replace the molecular sieve. With this, you always know when it was last filled or replaced.

5.7 OMNIS system – Mounting the electrode



CAUTION

Risk of cuts from sharp edges

Cut injuries from pieces of broken glass and/or sharp edges.

- Handle glass parts (e.g. electrode, sample beaker) with care.
- Only use undamaged glass parts.
- Dispose of damaged glass parts immediately.




Inserting the electrode in the electrode holder




Insert the electrode from above into the front opening of the electrode holder.




Push the green upper part of the electrode downwards until it stops. The green part of the electrode is flush with the lower edge of the electrode holder.

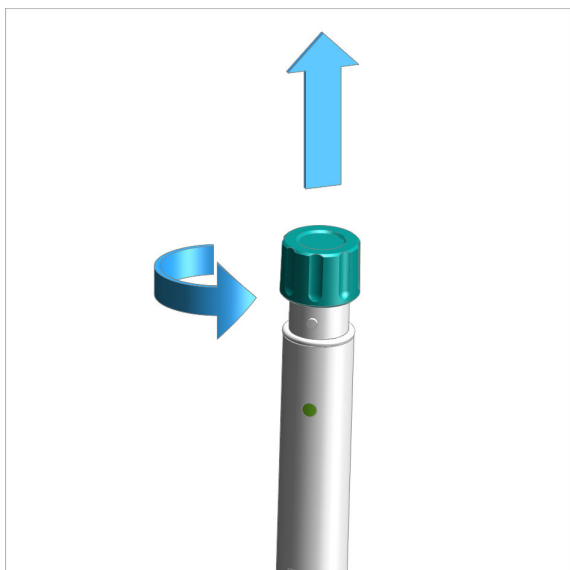

Unscrew and remove the protective cap from the electrode.

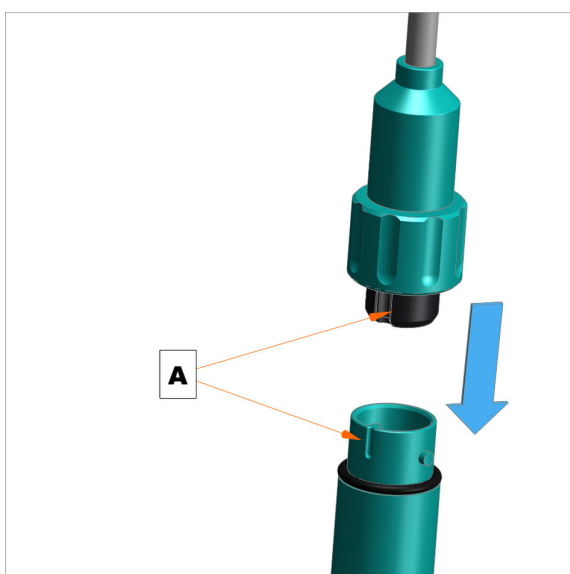

Plug in the plug of the electrode cable.
Observe the orientation (A).

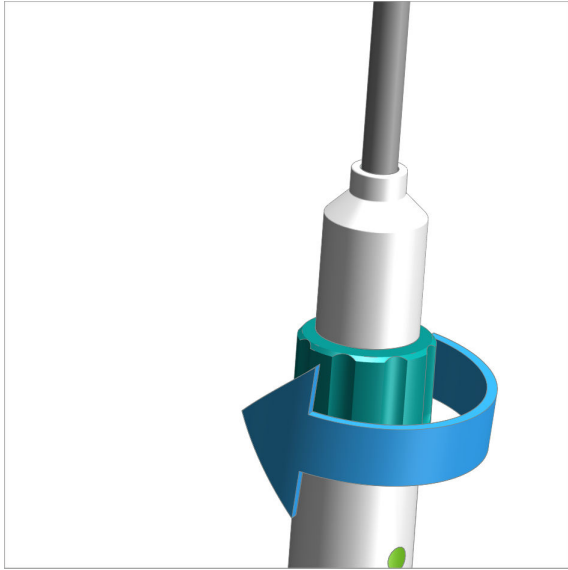


Screw the electrode cable tight.

5.8 Mounting the measuring module


The measuring module is supplied with mounted fastening screws. These fastening screws are used to secure the measuring module in the instrument to ensure trouble-free operation.

1 Opening the lid

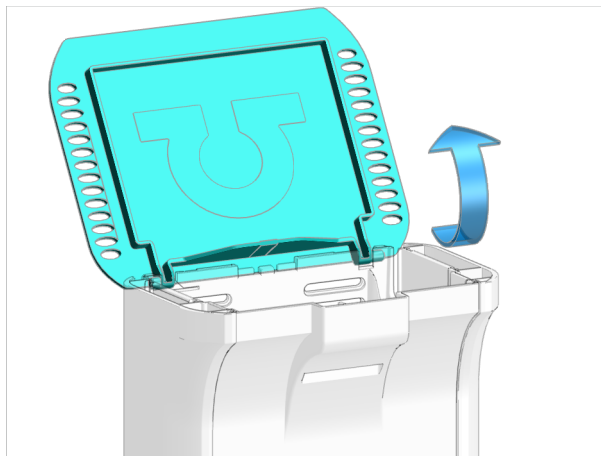

- Open the lid.


2 Removing the side parts

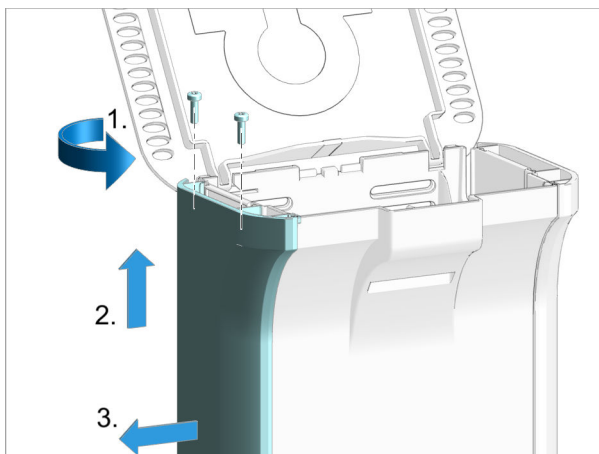

Carry out the following steps on **both sides** of the instrument.


- Remove the 2 screws from above using the hex key.
- Push the side covering upwards until it can be removed from the side.
- Remove the side covering to the side.

3 Removing the fastening screws

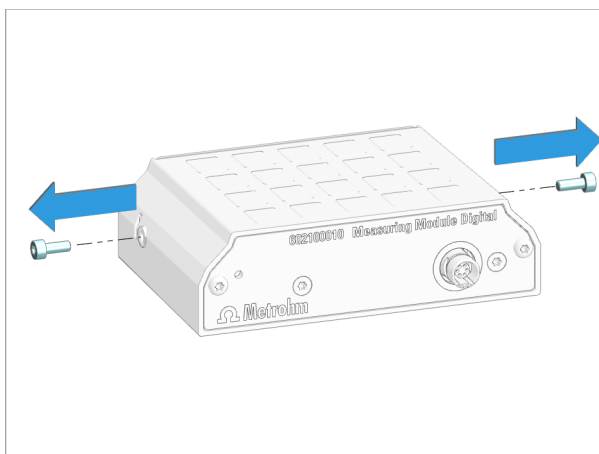

- Remove the 2 fastening screws from the measuring module using the hex key.

6 Mounting the side parts

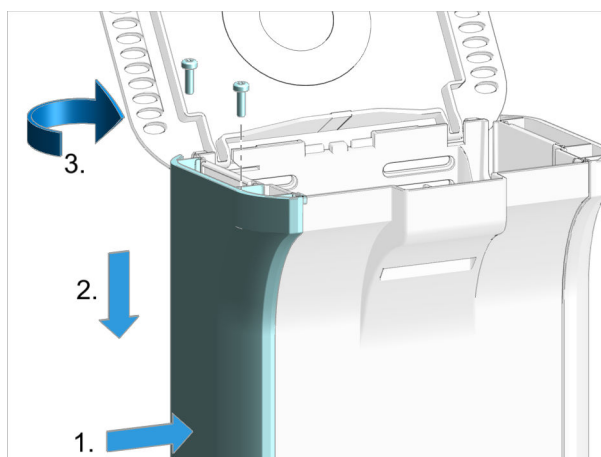

Carry out the following steps on **both sides** of the instrument.

- Position the side covering from the side in an elevated position.
- Insert the side covering into the guide rail and push it downwards.
- Insert the 2 screws and tighten them from above using the hex key.

5.9 OMNIS Titrator – Mounting the volumetric Karl Fischer titration cell

Setting up the volumetric KF titration cell

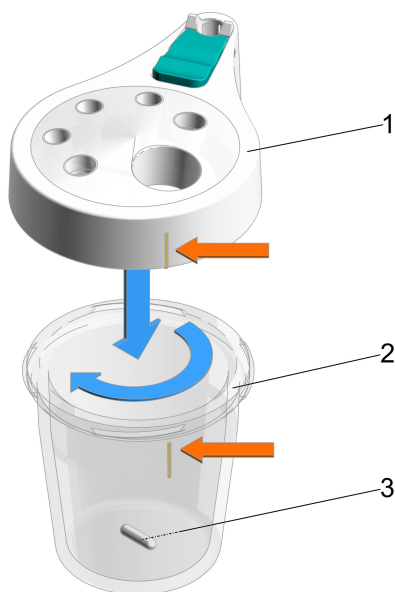


Figure 23 Preparing the volumetric KF titration cell

- | | |
|-----------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------|
| <p>1 KF titration vessel lid
(6.01405.010) or
for usage with homogenizer (6.01405.040)</p> | <p>2 KF titration vessel
20 – 90 mL (6.01406.220) or
50 – 150 mL (6.01406.250)</p> |
| <p>3 Stirring bar (without homogenizer!)
16 mm (6.1903.020) or
25 mm (6.1903.030)</p> | |

Assembling the volumetric KF titration cell

- 1 Select the KF titration vessel (23-2) of the required size and place a suitable stirring bar (23-3) inside.
- 2 Tighten the corresponding KF titration vessel lid (23-1) to the KF titration vessel.



NOTICE

When doing so, ensure that the color marking on the KF titration vessel aligns with the raised marking on the KF titration vessel lid.

The scale of the volumetric KF titration cell will thus face towards you.

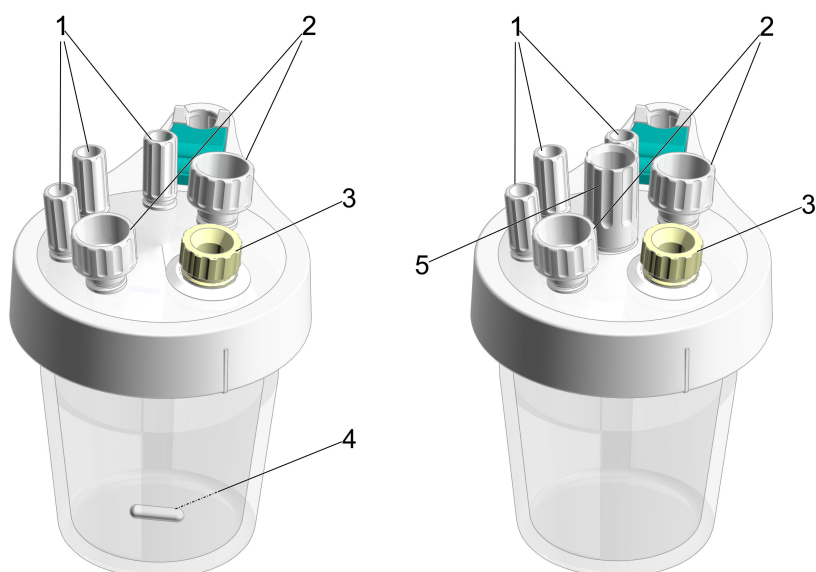


Figure 24 Volumetric KF titration cell (6.01405.010) – Volumetric KF titration cell for use with homogenizer (6.01405.040)

- | | |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------|
| <p>1 M10 screw nipple
(6.02709.010)</p> | <p>2 M12 screw nipple
(6.02709.030)</p> |
| <p>3 Septum stopper (or spoon for paste)</p> | <p>4 Stirring bar (without homogenizer!)
16 mm (6.1903.020) or
25 mm (6.1903.030)</p> |
| <p>5 Guide sleeve for homogenizer
Guide sleeve for Polytron PT 1300 D
(6.02709.050), further information see section (see "Optional: Mounting the volumetric KF titration cell for use with homogenizer on the OMNIS product (e.g. OMNIS Titrator)", page 61)</p> | |

Preparing the volumetric KF titration cell

- 1** Insert the 3 screw nipples (24-**1**) in the M10 openings of the KF titration vessel lid.
- 2** Insert the 2 screw nipples (24-**2**) in the M12 openings of the KF titration vessel lid.

Equipping the volumetric KF titration cell

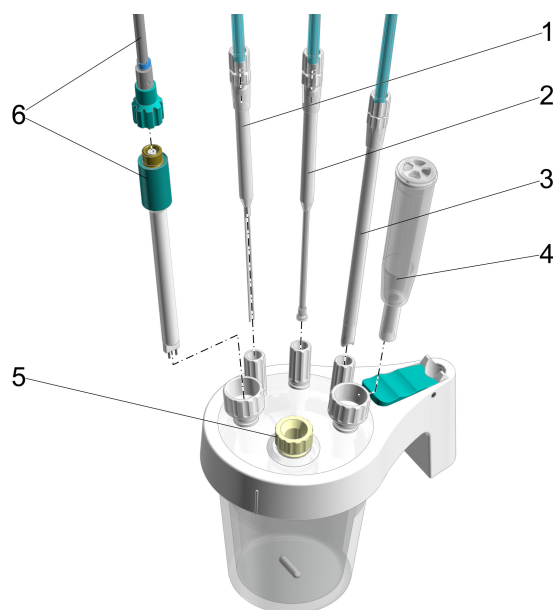


Figure 25 Volumetric KF titration cell – Equipping

1 Dosing tip (6.1543.110) with M8 tubing (6.1805.200)	2 Buret tip (6.01543.120) with M6 tubing (6.1805.100)
3 Aspiration tip (6.01543.000) with M8 tubing (6.1805.200)	4 Adsorber tube (6.01406.010)
5 Septum stopper with inserted septum (6.02709.020), stopper (6.02709.010) or OMNIS spoon for paste (6.02711.000)	6 Double Pt electrode (6.0338.100) with cable (6.02104.040)

Prerequisites:

- The volumetric KF titration cell is set up (see *"Setting up the volumetric KF titration cell"*, page 52).
- The adsorber tube with lid (25-4) is filled with fresh molecular sieve, see OMNIS Karl Fischer products – *Replacing the adsorber material* (see chapter 5.6, page 43).

- 1 Insert the dosing tip (25-**1**) in the M10 screw nipple on the left (24-**1**) and screw it tight.

The dosing tip should be located just above the stirring bar, but should not impede it.
- 2 Insert an M8 tubing in the M8 connector of the dosing tip (25-**1**) and screw it tight.

- 3** Insert the buret tip (25-2) of the cylinder unit in the M10 screw nipple in the middle (24-1) and screw it tight.

The antidiffusion valve of the buret tip should be located just above the stirring bar, but should not impede it.

- 4** Insert the M6 tubing in the M6 connector of the buret tip (25-2) and screw it tight.

- 5** Insert the aspiration tip (25-3) in the M10 screw nipple on the right (24-1) and screw it tight.

When solvent is aspirated, the end of the aspiration tip must touch the vessel base, but it must not inhibit the action of the stirring bar.

The aspiration tip can, if needed, be pulled out of the solvent.

- 6** Insert an M8 tubing in the M8 connector of the aspiration tip (25-3) and screw it tight.

- 7** Insert the double Pt electrode (25-6) in the M12 screw nipple on the left (24-2) and then tighten the screw nipple until it seals.

- 8** Screw an electrode cable with blue coding tightly onto the electrode (25-6).

- 9** Insert the adsorber tube (25-4) in the M12 screw nipple on the right (24-2) and then tighten the screw nipple until it seals.

- 10** Introduce the septum stopper (with septum inserted) into the front opening (24-3) of the titration vessel lid.

Select a different insert if required:

- Stopper
- Spoon for paste

Mounting the volumetric KF titration cell on the OMNIS product (e.g. OMNIS Titrator)

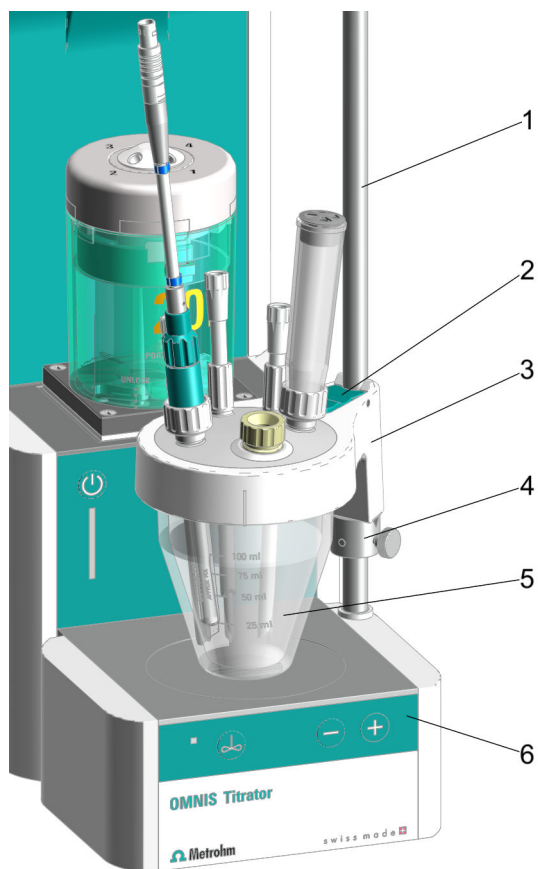


Figure 26 Mounting the volumetric KF titration cell on the OMNIS product – Overview

- | | | | |
|----------|----------------------------------------------------------------------------------------|----------|-----------------------------------------------------|
| 1 | Support rod | 2 | Locking lever on the KF titration vessel lid |
| 3 | KF titration vessel lid
(6.01405.010) | 4 | Clamping ring |
| 5 | KF titration vessel
20 – 90 mL (6.01406.220) or
50 – 150 mL (6.01406.250) | 6 | Magnetic stirrer |



NOTICE

If you want to mount the volumetric KF titration cell on the OMNIS product for use with a homogenizer, you can find additional information under (see "*Optional: Mounting the volumetric KF titration cell for use with homogenizer on the OMNIS product (e.g. OMNIS Titrator)*", page 61).

Prerequisites:

- The support rod is mounted on the OMNIS product with the clamping ring, see *Magnetic stirrer – Installing the accessories (see chapter 5.5, page 40)*.
The clamping ring is used as the lower stop for the KF titration vessel lid. Thus, the clamping ring ensures that the KF titration cell is always positioned at the same height and exactly in the center on the magnetic stirrer.
- The OMNIS titration system is connected to the cylinder unit and all parts are fully operational, see *OMNIS titration system – Mounting the bottle unit (see chapter 5.10, page 63)* and *Attaching the cylinder unit OMNIS (see chapter 5.4, page 38)*.
- The bottle caps are fully equipped and connected to the OMNIS Solvent Module, see .
- The volumetric KF titration cell is fully equipped, see steps *OMNIS Titrator – Mounting the volumetric Karl Fischer titration cell (see chapter 5.9, page 51)*.

1 Press the green locking lever (26-2) on the KF titration vessel lid (26-3).

2 Push the KF titration cell, comprising of (26-3) and (26-5), over the support rod (26-1).

3 Push the KF titration cell down so that it is positioned approx. 1 mm above the magnetic stirrer (26-6) and position it in the center of the magnetic stirrer.

To fix in place, release the green locking lever.

4 Push the clamping ring (26-4) under the KF titration vessel lid.

Rotate the clamping ring in such a way that the wedge on the KF titration vessel lid fits in the indent in the clamping ring.

5 Fix the clamping ring in place at the desired position with the knurled screw.

The position of the KF titration cell is now fixed by the clamping ring.

Mounting the joining elements for the volumetric KF titration cell on the OMNIS product (e.g. OMNIS Titrator) and OMNIS Solvent Module

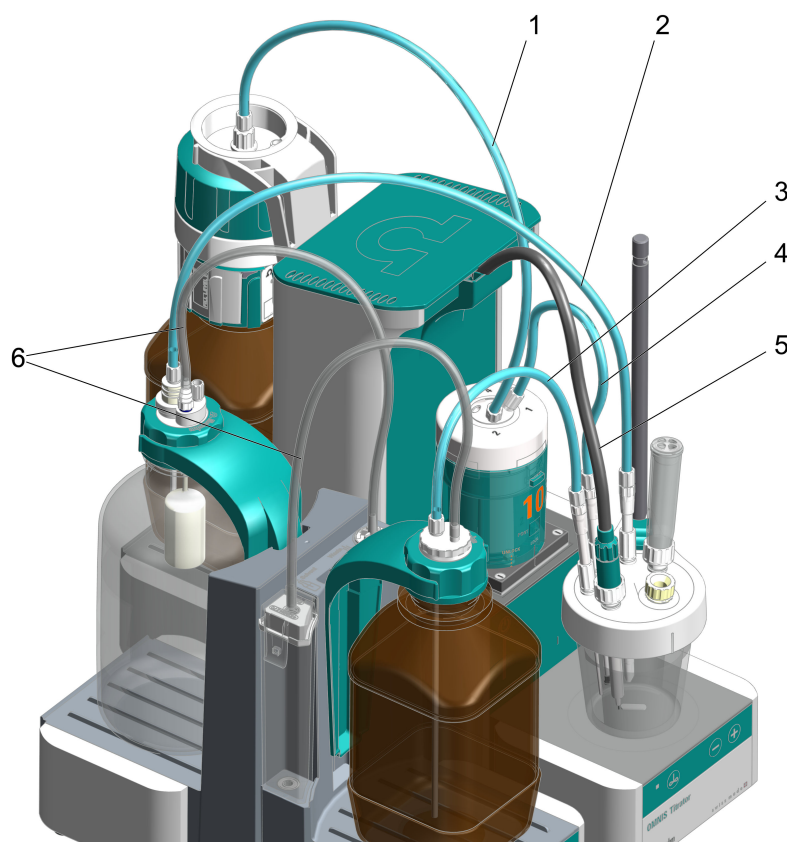


Figure 27 Connecting the KF titration cell with the OMNIS product and OMNIS Solvent Module

- | | | | |
|----------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 1 | M6 PTFE tubing from the fill port of the cylinder unit to the titrant bottle
(6.1805.100) | 2 | M8 PTFE tubing between aspiration tip and waste bottle (Waste)
Aspiration tip (6.01543.000) with M8 PTFE tubing (6.1805.200) leading to waste bottle (Waste) |
| 3 | M8 PTFE tubing between dosing tip and reagent bottle (Solvent)
Dosing tip (6.1543.110) with M8 PTFE tubing (6.1805.200) leading to the reagent bottle (Solvent) | 4 | M6 PTFE tubing between buret tip and dosing port 1 of cylinder unit
For adding the titrant via the buret tip (6.1543.200) with M6 PTFE tubing (6.1805.100) into the KF titration cell |
| 5 | Electrode cable leading to measuring module
Double Pt electrode (6.0338.100) with electrode cable (6.02104.040) leading to Measuring Module Analog | 6 | PVC tubings from the bottles to the OMNIS Solvent Module
(6.01804.210) |

Prerequisites:

- The adsorption cartridge is filled with molecular sieve, sealed tightly and mounted to the OMNIS Solvent Module, see .
- The Siphon Breaker and the bottle cap are fully equipped and connected to the OMNIS Solvent Module, see and *OMNIS titration system – Mounting the bottle unit (see chapter 5.10, page 63)*.
- The volumetric KF titration cell is fully equipped, see steps (see *"Equipping the volumetric KF titration cell", page 54*).

1 Connecting the volumetric KF titration cell with the titrant

Insert the M6 PTFE tubing from the buret tip (27-4) into the corresponding connector of the cylinder unit and screw it tight, see *Attaching the cylinder unit OMNIS (see chapter 5.4, page 38)*.

2 Connecting the volumetric KF titration cell with the bottles

Plug the M8 PTFE tubing from the dosing tip (27-3) onto the M8 connector of the Siphon Breaker of the reagent bottle (Solvent) and screw it tight.

- 3** Plug the M8 PTFE tubing from the aspiration tip (27-2) onto the olive of the ground-joint stopper SGJ 14/M8 on the GL 45 bottle cap of the waste bottle (Waste) and screw it tight.

4 Connecting the double Pt electrode with the electrode cable

Screw an electrode cable (27-5) with blue coding tightly onto the double Pt electrode.

Connecting a double Pt electrode to the analog measuring module

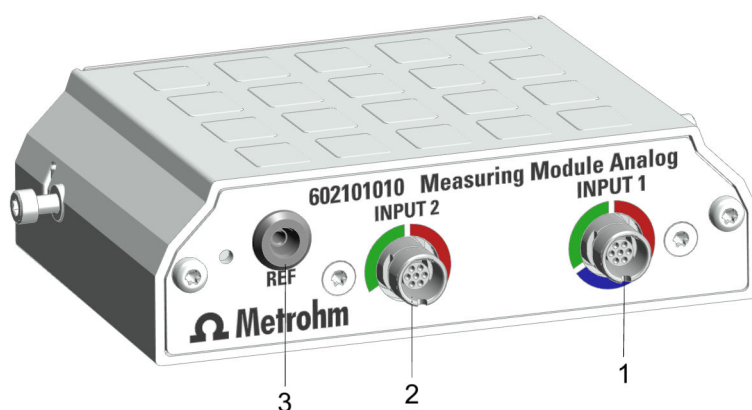


Figure 28 Measuring Module Analog – Interfaces

- | | | | |
|---|---------------------------|---|---------------------------|
| 1 | Connection socket INPUT 1 | 2 | Connection socket INPUT 2 |
| 3 | Connection socket REF | | |

Prerequisites:

An electrode cable with blue coding is screwed tightly onto the double Pt electrode.

- 1** If applicable, mount an analog measuring module in the OMNIS product, see *Mounting the measuring module* (see chapter 5.8, page 48).
- 2** Connect the electrode cable with blue coding (27-**1**) to the socket **INPUT 1** (see "Measuring inputs INPUT 1 and INPUT 2", page 33), see *Plugging in the sensor* (see chapter 5.11, page 64).



NOTICE

Note that only the socket **INPUT 1** can be used for KF titration, see blue marking!

Optional: Mounting the volumetric KF titration cell for use with homogenizer on the OMNIS product (e.g. OMNIS Titrator)

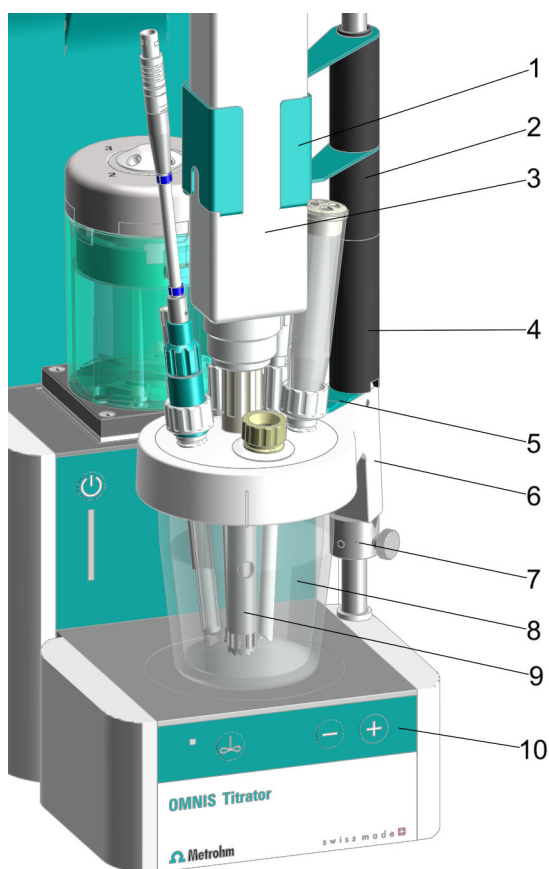


Figure 29 Mounting the volumetric KF titration cell for use with homogenizer on the OMNIS product – Overview

- | | |
|--------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|
| <p>1 Holder for homogenizer
Holder for Polytron PT 1300 D
(6.02008.010)</p> | <p>2 Spacer 35 mm</p> |
| <p>3 Polytron PT 1300 D
(2.1360.100) with dispersing aggregate</p> | <p>4 Spacer 65 mm</p> |
| <p>5 Locking lever on the KF titration vessel lid</p> | <p>6 KF titration vessel lid for use with homogenizer
(6.01405.040)</p> |

- 6** Connect the homogenizer handset (29-3) to the control instrument of the homogenizer (Polytron PT 1300 D).
- 7** Connect the homogenizer (Polytron PT 1300 D) (29-3) with the RS-232 cable to the PC.

Hint:

We recommend using the dispersing aggregates as follows:

- **Dispersing aggregate 125 mm**
 - Applications with viscous samples
 - Samples with a diameter that is smaller than the diameter of the aggregate
 - Powders and salts that are difficult to dissolve
- **Dispersing aggregate 157 mm**
 - Applications with solid samples
 - Samples with a diameter that is bigger than the diameter of the aggregate

5.10 OMNIS titration system – Mounting the bottle unit

The bottle unit in an OMNIS system is comprised of the following elements:

- Chemical bottle
- OMNIS bottle cap
- OMNIS Liquid Adapter

Several chemical manufacturers offer chemical bottles with an OMNIS bottle cap single-use. An OMNIS bottle cap multi-use is available for other commercially available chemical bottles. If the chemical bottle does not feature a red OMNIS bottle cap, replace the original lid of the chemical bottle with a bottle cap multi-use.


Mounting the bottle unit

- 1** Assemble the OMNIS Liquid Adapter.
- 2** Mount and connect the OMNIS Liquid Adapter.
- 3** If the chemical bottle is not sealed with a red OMNIS bottle cap single-use,
 - prepare an OMNIS bottle cap multi-use.
 - Remove the original lid from the chemical bottle.
 - Screw the OMNIS bottle cap multi-use onto the bottle.

- 4 ▪ Lock the OMNIS Liquid Adapter with the chemical bottle.
- Place the chemical bottle on the platform.

5.11 Plugging in the sensor

1 Plugging in the electrode cable

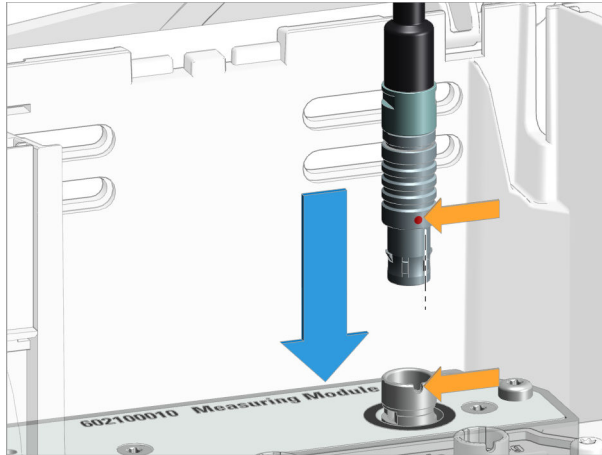


NOTICE

It should be easy to plug in the plug.

- Do not apply force if the plug cannot be inserted easily.
 - Rotate the plug to the right or left using light pressure until it latches in the socket.
-
- Align the red dot on the plug with the groove on the connection socket.
 - Plug in the plug until you can feel it snap in.

2 Guiding out the cable

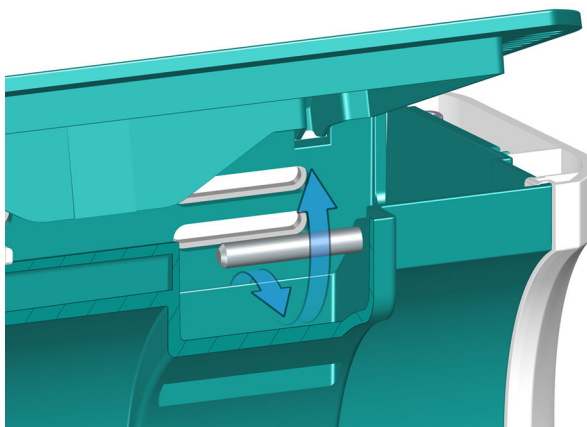

- Guide the cable out going under the bar.
- Close the lid.

See also

OMNIS Titrator – Interfaces (chapter 3.7, page 32)

OMNIS Titrator – Interfaces (chapter 3.7, page 32)

Measuring Module Analog – Interfaces (chapter 3.8, page 33)

Measuring Module Digital – Interfaces (chapter 3.9, page 34)

5.12 Plugging in the power cord



WARNING


Electrical potential

Contact with electrical potential can cause serious injuries or death.

- Operate the product only if it is in perfect condition. The housing must also be intact.
- Only use the product with the covers fitted.
- Protect live components (e.g. power supply unit, power cord, connection sockets) against moisture.
- Always have maintenance work and repairs on electrical components carried out by a regional Metrohm service representative.

Required accessories:

- Power cord:
 - Length: max. 2 m
 - Number of conductors: 3, with protective ground
 - Conductor cross-section: 3x min. 1.0 mm² / 18 AWG
- Instrument plug:
 - IEC 60320, type C13, 10 A
- Power plug:
 - 6.2122.XX0 (according to customer requirement), min. 10 A

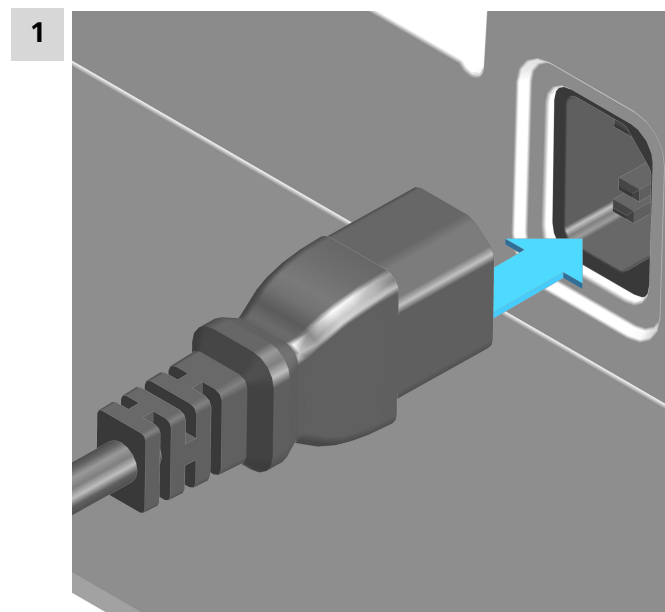

- Plug the power cord into the product's power socket. Only use permitted power cords.
- Connect the power cord to the energy supply.
- To power down the product, unplug the power cord from the energy supply.

6 Start-up

6.1 Initial start-up by Metrohm

As a basic rule, the initial start-up of the system is carried out by the regional Metrohm service representative.

Notes on handling



CAUTION

Piston wear

Solutions of solids (e.g. salts or hydroxides) increase the wear of the dosing piston, which may lead to leakage.

- Fill the cylinder with solution and move to the exchange position after each titration/dosing.

If no continuous sample throughput is ensured, fill the cylinder with solution and move to the exchange position, especially when using:

- concentrated solutions that tend to crystallize
- EDTA solutions, ultrapure solvents and ultrapure water
- organic solvents
- alkaline (e.g. KOH or isopropyl alcohol), corrosive or high-concentration reagents

The cylinder unit does not automatically move into the exchange position. In order to automatically move into exchange position after each titration/dosing, insert the **FILL** and **VALVE POS** commands in the method, see [software help](#).

For storing the cylinder unit over an extended period of time *Storing the cylinder unit OMNIS* (see chapter 8.3, page 78).

Usage of dosing tubing with antidiffusion valve

If using an antidiffusion valve, the maximum dosing rate is 150 mL/min.

The dosing rate can be saved on the memory chip of the cylinder unit: Enter the dosing rate in the OMNIS Software in **Properties ► Specific data**.

Usage of dosing tubing without antidiffusion valve

For usage without antidiffusion valve, do not immerse the dosing tubing in the sample solution.

There is a risk of back diffusion of the sample solution from the vessel into the tubing due to open tubing ends.



NOTICE

The cylinder unit and its parts cannot be autoclaved. The sterility of a germ-free solution cannot be guaranteed.

7.2.1 Attaching the cylinder unit OMNIS



NOTICE

Default settings for the ports 1 and 2

Port 1 is defined as dosing port and port 2 as fill port in the data chip default settings of the cylinder unit. The following instructions describe the default setting.

If the ports should be used differently from the default setting, adjust the ports in the OMNIS Software in **Properties ► Specific data**.

Preparing for attaching

- 1 Open the **Manual control** of the dosing unit in the OMNIS Software, see [software help](#).
- 2 Start the **Exchange position** function.

Attaching the cylinder unit



NOTICE

These instructions describe the default installation as defined in the OMNIS Software.

Prerequisite:

- Dosing drive: Valve coupling and piston rod are in the exchange position (port 2 is set).
- Cylinder unit: The piston stopper is flush with the base of the cylinder housing. The centering tube is in the correct position.


Required accessories:

- Wrench (6.2739.000)
- 2 FEP tubings (6.1805.100)
- Titration tip (6.1543.200)

- ## 1 Aligning the cylinder unit


Rotate the cylinder unit until the marking with the label **UNLOCK** is in line with the marking on the dosing drive.

2

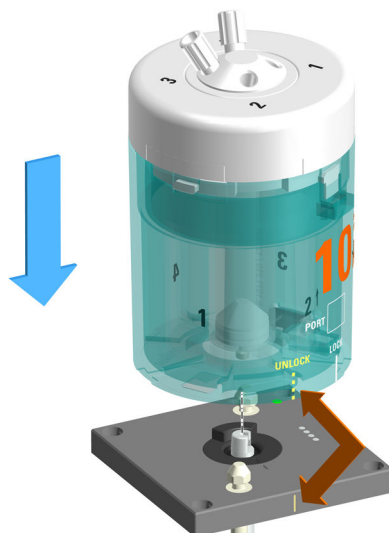


Set the cylinder unit down onto the two twistlocks straight from above.

3 Locking the cylinder unit

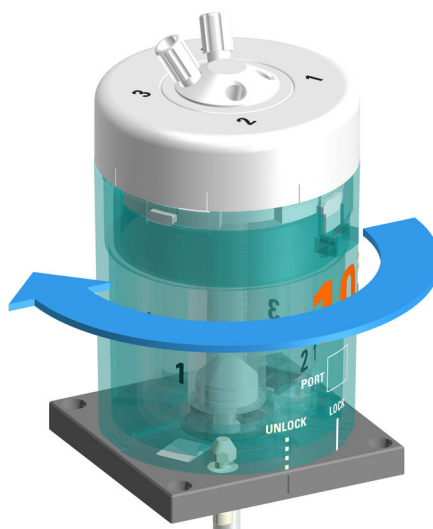

Rotate the cylinder unit to the left until it stops.
Use the marking with the label **LOCK** as a guide.

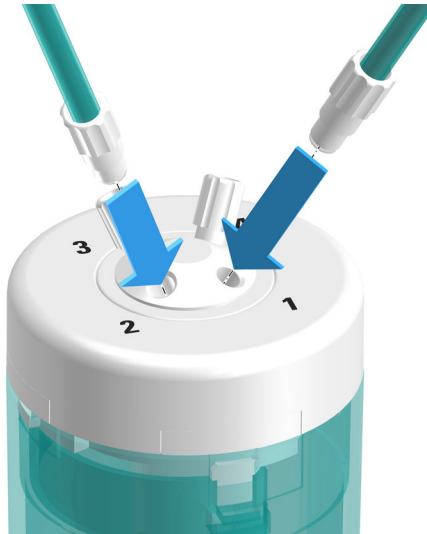

Tighten an FEP tubing (6.1805.100) into port 1.

This FEP tubing is used as dosing tubing. Tighten the other end to the titration tip (6.1543.200).

5 Tighten the other FEP tubing (6.1805.100) into port 2.

This FEP tubing is used as filling tubing. Tighten the other end to the OMNIS Liquid Adapter.

6 Firmly tighten the tubing with the wrench (6.2739.000).

See also

Cylinder unit OMNIS – Overview (chapter 3.3.3.1, page 19)

7.2.2 Removing the cylinder unit OMNIS

Preparing the removal

- 1 Open the **Manual control** of the dosing unit in the OMNIS Software, see [software help](#).
- 2 Start the **Empty** function.
- 3 Start the **Exchange position** function.

Removing the cylinder unit

Prerequisite:

- Dosing drive: Valve coupling and piston rod are in the exchange position (port 2 is set).
- Cylinder unit: The piston stopper is flush with the base of the cylinder housing. The centering tube is in the correct position.

1 Removing the tubing




Unscrew the dosing tubing and the filling tubing.

2 Unlocking the cylinder unit



Rotate the cylinder unit to the right as far as the **UNLOCK** position.

3 Raising the cylinder unit

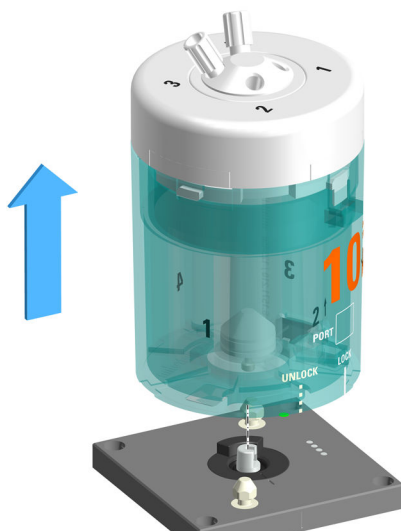

Raise the cylinder unit straight upwards.

See also

Cylinder unit OMNIS – Overview (chapter 3.3.3.1, page 19)

7.3 Magnetic stirrer – Operation

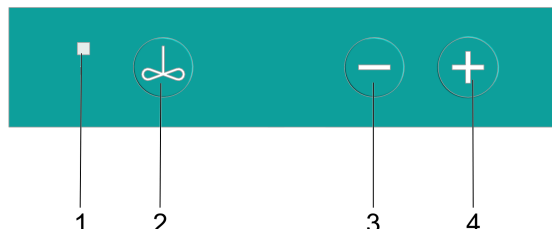


Figure 30 Magnetic stirrer – Control bar

<p>1 Status display Multi-colored</p>	<p>2 On/Off <i>Switching the magnetic stirrer on and off (see chapter 7.3.1, page 75)</i></p>
<p>3 Reduce stirring rate <i>Setting the magnetic stirrer (see chapter 7.3.2, page 76)</i></p>	<p>4 Increase stirring rate <i>Setting the magnetic stirrer (see chapter 7.3.2, page 76)</i></p>

Other functions in the software

The following functions can be executed only with the OMNIS Software (see [OMNIS Help](#)):

- **Deactivate keys**
The magnetic stirrer can be operated only via the software.
- **Switch over the keys for the rod stirrer**
The keys on the magnetic stirrer operate the rod stirrer.
- **Set the stirring direction**

7.3.1 Switching the magnetic stirrer on and off

1 Switching the magnetic stirrer on

Press the  key.

The magnetic stirrer stirs with the most recently used stirring rate.

2 Switching the magnetic stirrer off

Press the  key once again.

The magnetic stirrer stops.



NOTICE

If the magnetic stirrer is running at a high stirring rate, reduce the stirring rate before switching it off.

As an alternative, switch the magnetic stirrer on and off in the OMNIS Software under *Manual control*.

7.3.2 Setting the magnetic stirrer

The stirring rate can be adjusted in 15 steps.

Prerequisite:

The magnetic stirrer is switched on.

1 Increasing the stirring rate in steps

Press the  key.

Each pressing of the key increases the stirring rate by 1 step. The current stirring rate appears in the OMNIS Software under **Manual control**.

2 Reducing the stirring rate

Press the  key.

Each pressing of the key reduces the stirring rate by 1 step. The current stirring rate appears in the OMNIS Software under **Manual control**.

Alternatively, the stirring rate can also be set in the OMNIS Software under *Manual control*.



NOTICE

The stirring direction can be exclusively set in the OMNIS Software under **Manual control**.

8 Maintenance

8.1 Maintenance

Regularly perform maintenance work on the product to prevent malfunctions and to ensure a long service life.

- Metrohm recommends having the products maintained by specialist personnel of Metrohm AG as part of an annual service. Shorter maintenance intervals may be necessary if you frequently work with caustic and corrosive chemicals.
- Only perform maintenance work that is described in this instruction. Contact your regional Metrohm service representative for further maintenance work and repairs. The regional Metrohm service representative offers every form of technical advice for maintenance and service of all Metrohm products.
- Only use spare parts that meet the technical requirements of the manufacturer. Original spare parts always meet these requirements.

8.2 Cleaning the product surface

Regularly clean the product to prevent malfunctions and to ensure a long service life.

- Remove spilled chemicals immediately.
- Protect plug connections against contamination.



WARNING

Chemische Gefahrstoffe

Der Kontakt mit aggressiven chemischen Stoffen kann Vergiftungen oder Verätzungen verursachen.

- Persönliche Schutzausrüstung (z. B. Schutzbrille, Handschuhe) tragen.
- Absaugeinrichtung bei Arbeiten mit verdampfenden Gefahrstoffen verwenden.
- Verunreinigte Oberflächen reinigen.
- Nur Reinigungsmittel verwenden, die mit den zu reinigenden Materialien keine unerwünschten Nebenreaktionen auslösen.
- Chemisch verunreinigte Materialien (z. B. Reinigungsmaterial) vorschriftsmässig entsorgen.



WARNING

Electrical potential

Contact with electrical potential can cause serious injuries or death.

- Operate the product only if it is in perfect condition. The housing must also be intact.
- Only use the product with the covers fitted.
- Protect live components (e.g. power supply unit, power cord, connection sockets) against moisture.
- Always have maintenance work and repairs on electrical components carried out by a regional Metrohm service representative.

Prerequisite:

- The product is switched off and disconnected from the energy supply.

Required accessories:

- Cleaning cloth (soft, lint-free)
- Water or ethanol

- 1 Clean the surface with a damp cloth. Remove persistent contamination with ethanol.
- 2 Wipe the surface with a dry cloth.
- 3 Clean the connectors with a dry cloth.

8.3 Storing the cylinder unit OMNIS



NOTICE

If the cylinder unit will not be used for a prolonged period of time, rinse and fill the cylinder with deionized water to prevent the valve disk and distributor disk from sticking together, especially when using:

- concentrated solutions that tend to crystallize
- EDTA solutions, ultrapure solvents and ultrapure water
- organic solvents
- alkaline (e.g. KOH in isopropanol), corrosive or high-concentration reagents



- 1 Soak the filling tubing in a bottle with rinsing liquid.
- 2 Open the **Manual control** of the dosing unit in the OMNIS Software, see [software help](#).
- 3 Execute the **Preparing** function 2 to 3 times with rinsing solution.
- 4 If the cylinder unit is to be stored empty,
 - remove the filling tubing from the bottle with rinsing solution and
 - start the **Empty** function.
- 5 Start the **Exchange position** function.
- 6 Store the cylinder unit at ambient temperature and protect it from direct sunlight.



To automatically rinse the cylinder unit, download the method for automatically rinsing the cylinder unit as a template or create it yourself .



The cylinder unit has been removed from the dosing drive. *Removing the cylinder unit OMNIS (see chapter 7.2.2, page 72)*

1 Cleaning the cylinder housing




NOTICE

The cylinder housing is not dishwasher proof.

Clean the cylinder housing with lukewarm water and dishwashing detergent.

- 2 If the cylinder top piece is stuck, place the cylinder unit with the cylinder top piece facing down in warm water (possibly with a little dishwashing detergent) for at least 30 minutes.

3 Cleaning the electrical contacts of the cylinder unit

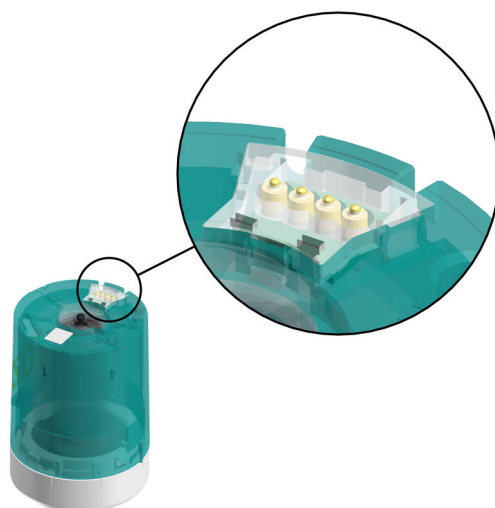


If the electrical contacts are only lightly contaminated, clean the electrical contacts with a cloth moistened with water.

- 4 If the electrical contacts are heavily contaminated,
 - add dishwashing detergent or ethanol to the moist cloth and clean the electrical contacts or
 - clean the electrical contacts in an ultrasonic bath with a little dishwashing detergent or ethanol.

Do not exceed 50 °C when drying. Use compressed air if necessary.


5 Cleaning the electrical contacts of the dosing drive

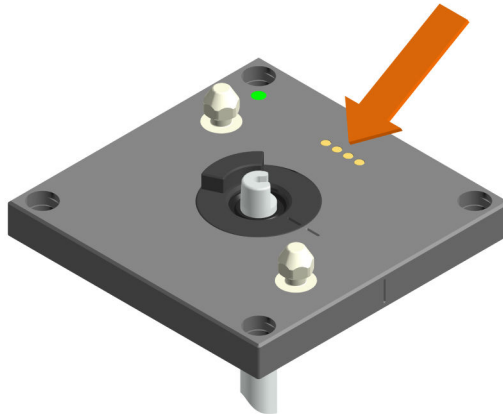

- If the electrical contacts are only lightly contaminated, clean the electrical contacts with a cloth moistened with water.
- If the electrical contacts are heavily contaminated, add dishwashing detergent or ethanol to the moist cloth and clean the electrical contacts.

8.5 Cleaning the OMNIS Liquid Adapter

Cleaning the aspiration tube

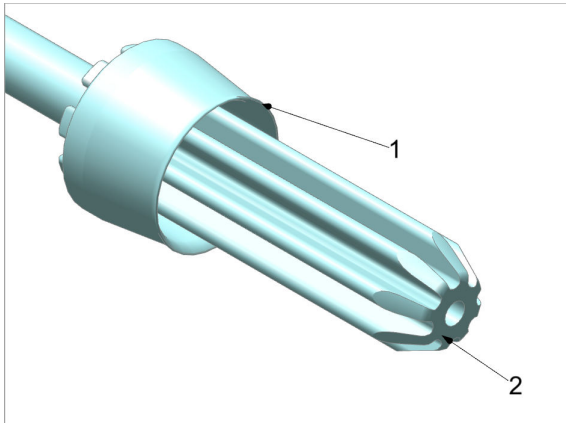

1. Rinse the aspiration tube thoroughly under running water.
Rub dry with a lint-free, damp cloth.
2. Check whether the sealing lip **(1)** and the sealing surface **(2)** are clean and intact.

If the sealing lip or the sealing surface of the aspiration tube can no longer be cleaned, mount a new aspiration tube.

Cleaning the OMNIS Liquid Adapter



CAUTION

Instrument damage through inward seepage of liquid


Property damage to the instrument or malfunction through the inward seepage of liquids (e.g. when cleaning).

The instrument is not resistant to splash water. Water can seep into the interior during cleaning and cause damage (e.g. to the electronics).

- Do not clean the instrument under running water.
- Do not use a wash bottle to clean the instrument.
- Only wipe the instrument thoroughly with a damp cloth.



Wipe the OMNIS Liquid Adapter thoroughly from the outside with a damp cloth.

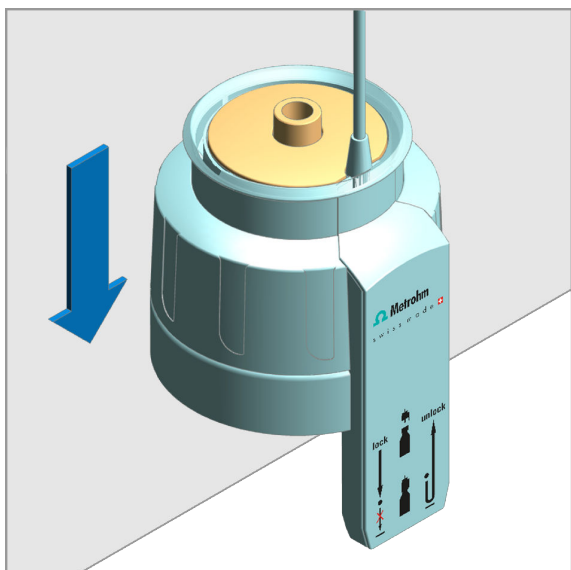

1. In order to access the lowered part better, set the OMNIS Liquid Adapter on a table edge. Push the OMNIS Liquid Adapter down and keep it pushed down. The lowered part rises.
2. Thoroughly wipe the surface and the slot for the aspiration tube with a damp cloth.
3. If the interior of the slot is contaminated, wipe it with a moistened cotton bud.
4. Release the housing of the OMNIS Liquid Adapter.

Contamination with organic substances

If the OMNIS Liquid Adapter is contaminated with organic substances, clean it with ethanol, methanol and/or isopropanol.



NOTICE

Do not use solvents that contain acetone to clean the OMNIS Liquid Adapter. Acetone will corrode the labels on the OMNIS Liquid Adapter.

9 Troubleshooting

Messages on malfunctions and errors are displayed in the control software or in the embedded software (e.g. on the display of an instrument) and contain the following information:

- Descriptions of causes of malfunctions (e.g. jammed drive)
- Descriptions of problems with the control (e.g. missing or invalid parameter)
- Information on how to solve the problem

System components with status display elements also indicate malfunctions and errors with a red flashing LED.

Troubleshooting on the product is often only possible with the control software or the embedded software (e.g. initializing, moving to a defined position).

See also


System – Signals (chapter 3.6, page 31)

9.1 Forcing a shutdown

Prerequisite:

The OMNIS main instrument cannot be switched off.

1

Press the on/off switch  for 8 seconds until the acoustic signal is heard in short intervals.

The acoustic signal sounds for 2 seconds. The status display goes out and the OMNIS main instrument is switched off.

11.1 Ambient conditions

Storage	+5 to +45 °C
----------------	--------------

Nominal voltage range	100 VAC - 240 VAC
------------------------------	-------------------

Power consumption max. 100 W

Fuse protection		
<i>Internal fuse</i>	4 ATH	cannot be replaced by the user

Nominal voltage	24 VDC	internal
------------------------	--------	----------

11.4 Measuring module – Energy supply

Power consumption	max. 0.6 W	-
Energy transmission	-	inductive coupling

11.5 OMNIS Titrator – Dimensions

Measurements

<i>Width</i>	142 mm
<i>Height</i>	358 mm
<i>Depth</i>	
Without magnetic stirrer	284 mm
With magnetic stirrer	400 mm

Weight

Type	
Without magnetic stirrer	4.4 kg
With magnetic stirrer	5.1 kg

11.6 Magnetic stirrer – Dimensions

Measurements

<i>Width</i>	142 mm
<i>Height</i>	70 mm
<i>Depth</i>	116 mm

Weight 700 g

11.7 Measuring module – Dimensions

Measurements

<i>Width</i>	105 mm
<i>Height</i>	31 mm
<i>Depth</i>	72 mm

Weight	approx. 420 g
---------------	---------------

11.8 OMNIS Titrator – Housing

Materials

<i>Lid</i>	PET	poly(ethylene terephthalate)
<i>Back panel</i>	AW-5754 H12/H22	aluminum, coated
<i>Base</i>	1.4301	high-grade steel
<i>Enclosure</i>	PBT	poly(butylene terephthalate)
<i>Front foils</i>	PET	poly(ethylene terephthalate), mat

IP degree of protection	IP 40
--------------------------------	-------

11.9 Magnetic stirrer – Housing

Materials

<i>Cover</i>	PBT	poly(butylene terephthalate)
<i>Base</i>		chromium steel
<i>Enclosure</i>	PBT	poly(butylene terephthalate)
<i>Front foils</i>	PET	poly(butylene terephthalate), mat

IP degree of protection	IP 40
--------------------------------	-------



<i>Cable type</i>	(min. FFTP)	shielded
<i>Cable length</i>	max. 10 m	from Metrohm accessories (link)
Measuring module		2 slots
<i>Power output</i>	max. 0.6 W	per measuring module
<i>Energy transmission</i>		inductive coupling
<i>Data transmission</i>		optical
Contacts	4	contact surfaces for cylinder unit

11.12 Measuring Module Analog – Connectors specifications

Measuring inputs

INPUT 1

Socket		round plug 7-pin, size 0, 45°
Potentiometric	pH, ISE, Redox	measuring input for potentiometric electrodes
Temperature	Temp.	measuring input for temperature sensors of the Pt1000 or NTC type for automatic temperature compensation
Polarizer	Pol.	measuring input for polarizable electrodes

INPUT 2

Socket		round plug 7-pin, size 0, 45°
Potentiometric	pH, ISE, Redox	measuring input for potentiometric electrodes

Temperature	Temp.	measuring input for temperature sensors of the Pt1000 or NTC type for automatic temperature compensation
REF		reference potential
Type	2 mm	
(INPUT 1 - INPUT 2)	pH, ISE, Redox	potentiometric differential measurement, with respect to REF

11.13 Measuring Module Digital – Connector specifications

Connection socket for a digital electrode

Type	round plug 6-pin, size 0, 60°
------	----------------------------------

11.14 Display specifications

Status display	LED	multi-colored
-----------------------	-----	---------------

11.15 Measuring Module Analog – Display specifications

Status display	LED	green-red
-----------------------	-----	-----------

11.16 Measuring Module Digital – Display specifications

Status display	LED	green
----------------	-----	-------

11.17 Measuring Module Analog – Measuring specifications

Potentiometric

<i>Measuring range</i>	–2,400 to +2,400 mV	
<i>Resolution</i>	1.56 μ V	
<i>Measuring accuracy</i>	± 0.5 mV	in the measuring range –2,000 mV to +2,000 mV
<i>Input resistance</i>	$\geq 1 \cdot 10^{12} \Omega$	
<i>Offset current</i>	$\leq \pm 1 \cdot 10^{-12}$ A	

Temperature

<i>Pt1000</i>		
Measuring range	–150 to +250 °C	
Resolution	approx. 0.002 °C	
Measuring accuracy	± 0.4 °C	in the measuring range –20.0 to +150.0 °C

Polarizer

<i>I_{pol} DC</i>		
Polarization current	–200.0 to +200.0 μ A	adjustable in 0.5 μ A steps
Measuring range	–2,400 to +2,400 mV	
Measuring resolution	0.1 mV	

Potentiometric differential measurement

<i>Measuring range</i>	–2,400 to +2,400 mV	
<i>Measuring resolution</i>	1.56 μ V	
<i>Measuring accuracy</i>	± 1.0 mV	in the measuring range –2,000 mV to +2,000 mV

Reference conditions

min. 30 minutes in operation

annual

applies for all measuring ranges without sensor error, under reference conditions, measuring interval 100 ms

11.18 OMNIS Titrator – Liquid Handling specifications

Dosing drive

<i>Dosing resolution</i>	100,000	steps per cylinder volume
<i>Dosing accuracy</i>	0.01%	typically, of the dosed volume with respect to the cylinder volume

11.19 Magnetic stirrer – Specifications

OMNIS product version		with attached magnetic stirrer
Adjustment range for rotational speed	+1 to +15	Rotation in counter-clockwise direction (seen from above)
	−1 to −15	Rotation in clockwise direction (seen from above)
Rotational speed change per step	120 rpm	
Maximum rotational speed	1,800 rpm	
Stirring bar lengths	8, 12, 16, 25, 30 mm	