

Aquatrode plus with plug-in head U



NEW:

The proven Aquatrode plus for water analysis – now also with removable electrode cable incl. temperature sensor Pt1000


The advantages at a glance

The Aquatrode plus with plug-in head U was specially developed for measurements and titrations in low-ion and poorly buffered samples, such as surface waters and deionized water.

- **Special membrane glass** with low resistance for measurement in low-ion solutions. This guarantees very short response times even in such difficult samples (fast measurement or end point titration)
- **Fixed ground-joint diaphragm** with low resistance and very low, clearly defined electrolyte flow for fast, very accurate and reproducible pH measurement / titration
- **Double-junction design** with maintenance-free inner electrolyte (gel) enables the bridge electrolyte to be adapted to the measurement conditions
- **Shielding of the inner tube** prevents interference due to electrostatic effects
- In addition to these advantages, the Aquatrode plus now offers **easy storage / fast electrode replacement / complete flexibility** due to its removable cable

Ordering information

6.0257.600	Aquatrode plus with Pt1000 (plug-in head U)
6.2104.600	Electrode cable, plug-in head U / plug F, 2 x B, 1 m
6.2104.610	Electrode cable, plug-in head U / plug F, 2 x B, 2 m

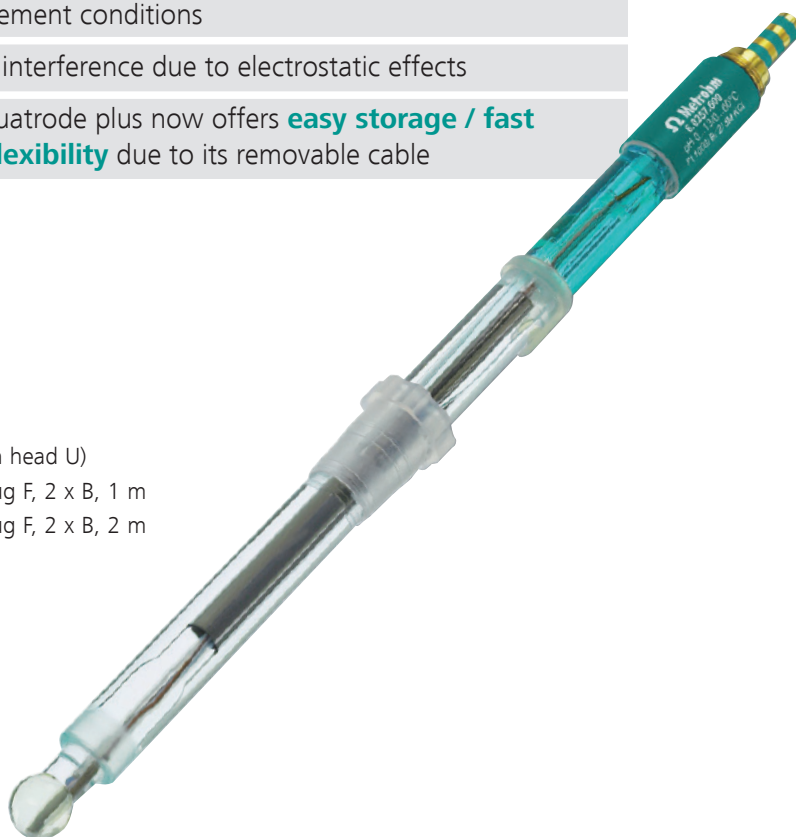

www.metrohm.com

 **Metrohm**