




Ecotrode Gel

The proven Ecotrode now maintenance-free as gel electrode with aging indicator

HIGHLIGHTS

- Ideal for routine measurements in similar kinds of samples
- Maintenance-free due to gel electrolyte
- Enhanced reliability due to ageing indicator

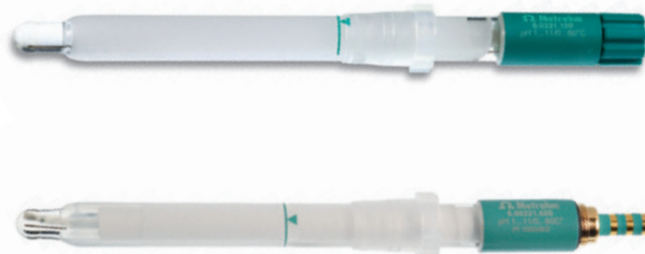

Enhanced reliability for your routine pH measurements

The maintenance-free pH electrode features a gel electrolyte with an integrated temperature sensor (Pt1000) and is ideally suited for routine measurements in similar kinds of samples. The electrode is stored in $c(\text{KCl})=\text{sat.}$ (6.2308.000) and is not suitable for solutions with low ion concentration.

No refilling / monitoring of the electrolyte necessary

An indicator shows in due time when the electrode should be replaced: With ongoing usage the white electrolyte gel loses its color and becomes greyish. As long as there is still some white electrolyte shown below the marker, the electrode is good for use and accurate results are guaranteed.



Everything ok: With ongoing usage, the white electrolyte gel loses its color towards the green marker.



ORDERING INFORMATION

Instrument	
6.0221.100	Ecotrode Gel
6.00221.600	Ecotrode Gel (with Pt1000, and cable plug-in head U)
Options	
6.2104.020	Electrode cable for plug-in head G / plug F, 1 m)
6.2104.030	Electrode cable for plug-in head G / plug F, 2 m)
6.2104.600	Electrode cable for plug-in head U / plug F, 2 x B, 1m)
6.2104.610	Electrode cable for plug-in head U / plug F, 2 x B, 2m)
6.2308.000	Electrolyte KCl (saturated), 250 mL

