

Audit Report: Categorisation acc. to GAMP

## MagIC Net 2.2



Metrohm Ltd.  
CH-9101 Herisau  
Switzerland  
Tel. +41 71 353 85 85  
Fax +41 71 353 89 01

[www.metrohm.com](http://www.metrohm.com)

**Date:** 17.05.2010  
**Author:** Sieghard Wagner, mech. Engineer (grad.), Chemengineering Business Design GmbH  
**Goal:** Classification of *MagIC Net 2.2* in software categories according to GAMP

**Description:** *MagIC Net* is a control and database software for ion chromatography devices, dosing devices and sample processors for creating methods, carrying out determinations, data evaluation and generating reports. *MagIC Net* is the successor of IC Net.  
*MagIC Net* was developed and manufactured by Metrohm Ltd. in accordance with the requirements of ISO 9001 quality system regarding design, manufacturing and maintenance.

**Categorisation:** *MagIC Net* is a standard - software according to GAMP category 3.

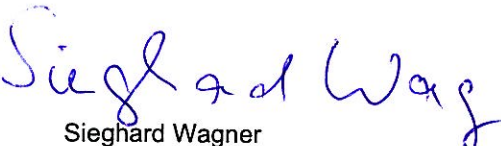
**Justification:**

The software configuration is limited to:

- Adaptation of the system's runtime environment, e.g.:
  - Configuration of the connected instruments,
  - Installation of databases
  - Definition of users and user groups (with pre-defined privileges)
  - Configuration of security settings.

However, no structural modifications or customisations of the software concerning user-specific procedures are done.

- Creation of *methods*:  
*Methods* are consisting of a *device arrangement*, a *time program* and *evaluation parameters*. As part of the regular system operation, methods are adapted to specific analytical procedures on a case by case bases. The resulting method settings have to be approved accordingly. That applies also and in particular to calculations and reports.

  
Sieghard Wagner