

Release Notes tiamo 2.5 Patch 1

1 Purpose

This documents describes the differences between the new software version **tiamo 2.5 Patch 1 (Build 116-1)** and its predecessor **tiamo 2.5 (Build 116)**.

2 New features

New instrument firmware

- New firmware version for the 851 Titrandos: 5.851.0013
- New firmware version for the 852 Titrandos: 5.852.0013

3 Fixed bugs and problems


Configuration program part

In **tiamo 2.5 light**, if in the **Security settings** under **Login policies/Password policies** the **Enforce login with user name** check box was active and then the **[OK]** button was pressed, the program would crash.

4 Compliance

The current software version does not contain any modifications that affect conformity of **tiamo** regarding 21 CFR Part 11 and GAMP.

Herisau, May 19, 2017

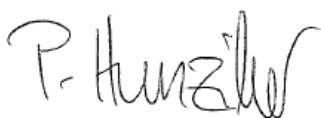

P. Hunziker
Vice President,
Head of Development



Dr. M. Saba
Head of Quality Management