


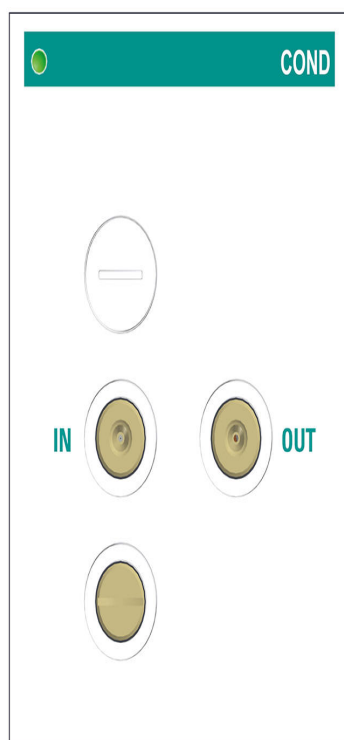

OMNIS Conductivity Detector Unit

2.1056.0810

The OMNIS Conductivity Detector Unit functions as a conductivity detector in the OMNIS IC System.

The conductivity detector continuously measures the conductivity of the pumped liquid. The sample is separated into its ionic components in the separation column. Afterwards, the separated ions generate a change in the measurement signal in the conductivity detector compared to the baseline. This manifests itself in the form of evaluable peaks in the chromatogram.

The conductivity change depends on the type and concentration of the electrolyte and on the sample volume applied. Low baseline noise enables low limits of detection and therefore more precise measurements. The conductivity detector exhibits high thermal stability and thus guarantees reproducible measuring conditions.


Specifications:

- 0–35000 $\mu\text{S}/\text{cm}$ without range switching
- Cell volume: 0.3 μL
- Ring-shaped electrodes made of stainless steel X2CrNiMo17-12-2 (316 L), compatible with MSA
- Maximum operating pressure: 10.0 MPa (100 bar)
- Peltier controlled cell temperature: 15–60 °C in increments of 1 °C
- Temperature stability: 0.001 °C

- Baseline noise: 0.1 nS/cm typical for sequential suppression
- Capillaries: ID 0.18 mm

Below, the accessories are grouped into Scope of delivery and Optional accessories.
Please keep this printout at hand for ordering replacement material.
These lists may be subject to change.

Scope of delivery 2.1056.0810

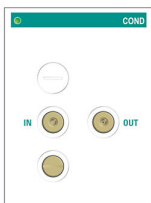Order Number	Description	Quantity
1.1056.0810	OMNIS Conductivity Detector Unit  The OMNIS Conductivity Detector Unit functions as a conductivity detector in the OMNIS IC System. The conductivity detector continuously measures the conductivity of the pumped liquid. The sample is separated into its ionic components in the separation column. Afterwards, the separated ions generate a change in the measurement signal in the conductivity detector compared to the baseline. This manifests itself in the form of evaluable peaks in the chromatogram. The conductivity change depends on the type and concentration of the electrolyte and on the sample volume applied. Low baseline noise enables low limits of detection and therefore more precise measurements. The conductivity detector exhibits high thermal stability and thus guarantees reproducible measuring conditions.	1 PCS
6.01831.190	PEEK capillary COND – MFM/WC, 400/0.18 mm  Pre-cut PEEK capillary with identification label for OMNIS IC	1 PCS

6.01831.200

PEEK capillary COL - COND, 350/0.18 mm

1 PCS

Pre-cut PEEK capillary with identification label for OMNIS IC

