



KFT accessories for IONIC Titrator

6.05609.020

Complete accessory set for volumetric Karl Fischer titration with the IONIC Titrator

Below, the accessories are grouped into Scope of delivery and Optional accessories.
Please keep this printout at hand for ordering replacement material.
These lists may be subject to change.

Scope of delivery 6.05609.020

Order Number	Description	Quantity
6.1903.020	Stirring bar / 16 mm Stirring Bar with magnetic core, PTFE coated, length 16 mm	1 PCS
		
6.1903.030	Stirring bar / 25 mm Stirring bar with magnetic core, PTFE covering, length 25 mm.	1 PCS
		
6.00346.100	Double Pt ring electrode for volumetry Indicator electrode for Karl Fischer titration which, thanks to its two fused platinum rings, is more robust and easier to clean than the conventional platinum pin electrode (6.0338.100).	1 PCS
		

6.01405.010

OMNIS titration vessel lid for volumetric Karl Fischer titration 1 PCS


Titration vessel lid for volumetric Karl Fischer titration with an OMNIS titrator.



6.01406.010

Adsorber tube with lid 1 PCS

Adsorber tube with lid for OMNIS Karl Fischer titrator, volumetric

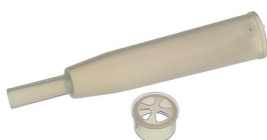

6.01406.220

20 – 90 mL titration vessel 1 PCS

20 – 90 mL titration vessel with screw top for use with 6014050X0 titration vessel lids.



6.01408.010

Set of seals for Karl Fischer titration cell 1 PCS

Complete set of seals, suitable for 6.05609.000 and 6.05609.010.



6.01543.120

Antidiffusion tip

1 PCS

Titration tip with antidiffusion valve, suitable for OMNIS Karl Fischer titrator.



6.02709.010

Stopper

3 PCS

Plastic stoppers, suitable for 2.1001.4220 OMNIS Karl Fischer Titrator.



6.02709.020

Septum stopper

1 PCS

Plastic septum stopper, suitable for 6014050X0 titration vessel lids.



6.02709.030

Screw nipple

2 PCS

Plastic screw nipple, suitable for 6014050X0 titration vessel lids.



6.02801.000

Cotton for adsorber tube

1 PCS

Wad of cotton (4 pieces), suitable for OMNIS 601406010 adsorber tube.



6.1448.010

Septum 12 mm, 5 pieces

1 PCS

Set of 5 items.



6.2816.020

Syringe 10 mL with Luer connection

2 PCS

For various applications in IC and VA



6.2816.030

Needle with Luer connector

2 PCS

Needle for filling in mercury into the Multi-Mode Electrode and for sample addition in Karl Fischer titration.



6.02104.060

Electrode cable (pol.) plug P, 0.75 m

1 PCS



Cable for connecting an analog polarizable electrode (plug-in head G) to an analog measuring module or measuring interface.