



Conductometric Titrator w. Exchange Unit

2.856.3210

Complete package for conductometric titration for the determination of alpha acids in hop products in accordance with EBC 7.4.

Including 856 Conductivity Module, 801 Magnetic Stirrer, 805 Dosimat, 20 mL exchange unit, conductivity measuring cell, and tiamo™ light.

Phạm vi giao hàng 2.856.3210

Qt.	Order no.	Description
1 PCS	1.801.0010	801 Stirrer

Magnetic stirrer without stand for supplementing the Titrino plus, Dosimat plus, Titrandos, Sample Processors, 805 Dosimats and 780/781 pH Meters. With permanently attached cable for MSB (Metrohm Serial Bus).



1 PCS

1.805.0010

805 Dosimat

Dosing device for the Titrando and Sample Processors with read/write hardware for intelligent Exchange Units. With permanently attached cable. Without Exchange Unit.



1 PCS

1.856.0010

856 Conductivity Module

Conductivity measuring module as supplement to an existing Titrando system or "stand-alone" in combination with a 900 Touch Control. With the 856 Conductivity Module, not only conductivity and temperature can be determined, but also TDS and salinity. It supports state-of-the-art conductivity measuring cells, i.e. 5-ring measuring cells.

Thanks to the galvanically isolated measuring input, pH value and conductivity can be measured in the same beaker without interference.

The Conductivity Module is equipped with two USB interfaces for connecting printers, barcode readers or sample changers and four MSB interfaces for stirrers or Dosinos.

Both in conjunction with the 900 Touch Control as stand-alone instrument and also integrated in *tiamo*™ full (from 2.0), it is in compliance with GLP and FDA 21 CFR part 11 requirements.



1 PCS

6.0915.100

5-ring conductivity measuring cell $c = 0.7 \text{ cm}^{-1}$ with Pt1000 (fixed cable)

5-ring conductivity measuring cell with cell constant $c = 0.7 \text{ cm}^{-1}$ (guide value), with integrated Pt1000 temperature sensor and with fixed cable (1.2 m) for connecting to an 856 Conductivity Module.

This sensor is suitable for measurements of medium conductivities ($5 \mu\text{S}/\text{cm}$ to $20 \text{ mS}/\text{cm}$), e.g., in:

- drinking water
- surface water
- wastewater



1 PCS

6.1903.020

Stirring bar / 16 mm

Stirring Bar with magnetic core, PTFE coated, length 16 mm



1 PCS

6.1903.030

Stirring bar / 25 mm

Stirring bar with magnetic core, PTFE covering, length 25 mm.




1 PCS

6.2001.060

Support stand without support rod

To attach an 801, 804, 803 Stirrer to a Titrando, Titrino plus, Dosimat plus or 856 and 867.

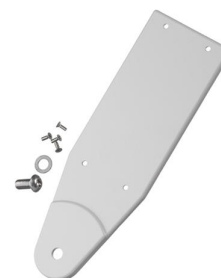

1 PCS

6.2013.010

Clamping ring

For support rods with a diameter of 10 mm.



1 PCS

6.2016.070

Support rod / 400 mm



1 PCS

6.2021.020

Electrode holder

Electrode holder for 4 electrodes and 2 buret tips



1 PCS

6.2065.000

Stacking frame for 846 Dosing Interface, 856 Conductivity Module, 867 pH Module

For fixing the Reagent Organizer on top of the Dosing Interface



1 PCS

6.2151.000

Cable USB A – mini-DIN 8-pin

Controller cable



1 PCS

6.2621.070

5 mm hex key for IC Sample Processors



1 PCS

6.2621.130

Hexagon key 2 mm

2 mm.



1 PCS

6.2739.010

Wrench

For exchange units.



1 PCS

6.3026.220

Exchange Unit 20 mL

Exchange unit with integrated data chip with 20 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, anti-diffusion buret tip and standard amber glass reagent bottle.

